

AD-A191 925

CHARACTERIZATION OF YTTRIUM ALUMINUM GARNET/YAG/USING  
14 MEV NEUTRON ACTIVATION(U) ARMY LAB COMMAND WATERTOWN  
MA MATERIAL TECHNOLOGY LAB F BURNS ET AL. DEC 87

1/1


UNCLASSIFIED

MTL-TR-87-2

F/G 7/2

ML



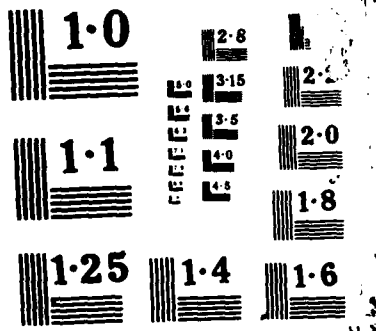

2

MTL TR 87-62

AD

AD-A191 925

DTIC FILE COPY

# CHARACTERIZATION OF YTTRIUM ALUMINUM GARNET /YAG/ USING 14 MeV NEUTRON ACTIVATION

FORREST BURNS, MICHAELINE F. CHEN and JAROSLAV L. CASLAVSKY  
EMERGING MATERIALS DIVISION

December 1987

Approved for public release; distribution unlimited.

DTIC  
ELECTE  
FEB 23 1988  
S H D



US ARMY  
LABORATORY COMMAND  
MATERIALS TECHNOLOGY  
LABORATORY

U.S. ARMY MATERIALS TECHNOLOGY LABORATORY  
Watertown, Massachusetts 02172-0001

88 2 22 174

The findings in this report are not to be construed as an official Department of the Army position, unless so designated by other authorized documents.

Mention of any trade names or manufacturers in this report shall not be construed as advertising nor as an official indorsement or approval of such products or companies by the United States Government.

**DISPOSITION INSTRUCTIONS**

Destroy this report when it is no longer needed.  
Do not return it to the originator.

## UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

REPORT DOCUMENTATION PAGE		READ INSTRUCTIONS BEFORE COMPLETING FORM
1. REPORT NUMBER MIL TR 87-62	2. GOVT ACCESSION NO.	3. RECIPIENT'S CATALOG NUMBER
4. TITLE (and Subtitle) CHARACTERIZATION OF YTTRIUM ALUMINUM GARNET /YAG/ USING 14 MeV NEUTRON ACTIVATION		5. TYPE OF REPORT & PERIOD COVERED
		6. PERFORMING ORG. REPORT NUMBER
7. AUTHOR(s) Forrest Burns, Michaeline F. Chen and Jaroslav L. Caslavsky		8. CONTRACT OR GRANT NUMBER(s)
9. PERFORMING ORGANIZATION NAME AND ADDRESS U.S. Army Materials Technology Laboratory Watertown, Massachusetts 02172-0001 SLCMT-EMC		10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS D/A Project: 1L161102.H42 AMCMS Code: P611102.H42
11. CONTROLLING OFFICE NAME AND ADDRESS U.S. Army Laboratory Command 2800 Powder Mill Road Adelphi, Maryland 20783-1145		12. REPORT DATE December 1987
		13. NUMBER OF PAGES 15
14. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office)		15. SECURITY CLASS. (of this report)  Unclassified
		15a. DECLASSIFICATION/DOWNGRADING SCHEDULE
16. DISTRIBUTION STATEMENT (of this Report)  Approved for public release; distribution unlimited.		
17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report)		
18. SUPPLEMENTARY NOTES  Published in the Journal of Radio-Analytical Chemistry, v. 118, no. 3, 1987, p. 155-167.		
19. KEY WORDS (Continue on reverse side if necessary and identify by block number)  14 MeV neutron, Characterization, Gamma-ray spectroscopy, Yttrium aluminum garnet Oxygen analysis, Fast neutrons		
20. ABSTRACT (Continue on reverse side if necessary and identify by block number)  SEE REVERSE SIDE		

Block No. 20



*alpha*

ABSTRACT

A non-destructive method for the analysis of yttrium, aluminum and oxygen in Yttrium Aluminum Garnet /YAG/ by 14 MeV neutron activation analysis had been developed and evaluated. The fast neutron induces primarily /n,2n/, /n,p/ or /n, $\alpha$ / reactions with yttrium, aluminum, and oxygen to produce isotopes with measurable characteristic gamma-ray spectra. Concentrations of Y, Al and O were determined on the basis of calibrated nuclear decay emission spectra recorded on selected standards. After each analysis the crystal was returned to the furnace for further adjustments of the growth parameters, and any resultant elemental variations were observed in the next analysis. The accuracy of this method was determined by repeated analyses with various "YAG" crystals from a series of separate furnace runs.

*Keywords:*



Accession For	
NTIS GRA&I	<input checked="" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification	
By _____	
Distribution/	
Availability Codes	
Dist	Avail and/or Special
A-1	

CHARACTERIZATION OF YTTRIUM ALUMINUM GARNET /YAG/  
USING 14 MeV NEUTRON ACTIVATION

F.C. Burns, M.F. Chen, J.L. Caslavksy

Materials Characterization Division,  
U.S. Army Materials Technology Laboratory,  
Watertown, Massachusetts, USA

Received 15 April 1987

Accepted 21 April 1987

A non-destructive method for the analysis of yttrium, aluminum and oxygen in Yttrium Aluminum Garnet /YAG/ by 14 MeV neutron activation analysis has been developed and evaluated. The fast neutron induces primarily  $/n,2n/$ ,  $/n,p/$  or  $/n,\alpha/$  reactions with yttrium, aluminum, and oxygen to produce isotopes with measurable characteristic gamma-ray spectra. Concentrations of Y, Al and O were determined on the basis of calibrated nuclear decay emission spectra recorded on selected standards. After each analysis the crystal was returned to the furnace for further adjustments of the growth parameters, and any resultant elemental variations were observed in the next analysis. The accuracy of this method was determined by repeated analyses with various "YAG" crystals from a series of separate furnace runs.

## INTRODUCTION

YAG / $Y_3Al_5O_{12}$ / is widely used in the solid state laser field as the light emitting source. For example, a new holographic interferometer installed at the NASA Ames Research Center is based on a frequency-doubled YAG laser<sup>1</sup>. Its unique physical properties, high peak power, good temporal resolution, and broad spectral coverage are critical to a variety of applications. In order to produce these high quality YAG lasers<sup>3</sup>, it is necessary to synthesize high purity crystals with the correct stoichiometric composition. This is difficult to accomplish due to the extreme environment that is used for the synthesis and subsequent single crystal growth. In this process, yttrium oxide and aluminum oxide are mixed and then melted at temperatures above 2000 °C /Ref. 2/. It is important to find the best method to mix the oxides prior to heating to the final temperature so that large compositional variation in the crystals will not occur. The resulting crystal must be characterized in order to determine the correct chemical composition of the melt. X-ray powder diffraction<sup>4</sup> can be helpful in the detection of various phases present.

Accurate analysis for oxygen determination using 14 MeV neutrons has been a standard technique for many years<sup>5</sup>. Aluminum and yttrium each form isotopes which decay with the emission of gamma rays and can be detected and resolved with a Ge/Li/ detector. Thus the concentration of all three elements can be determined simultaneously on the same crystal.

It is possible to use standard wet chemical techniques to analyze YAG quantitatively for aluminum and yttrium, but not oxygen<sup>6</sup>. Also, such analyses are destructive.



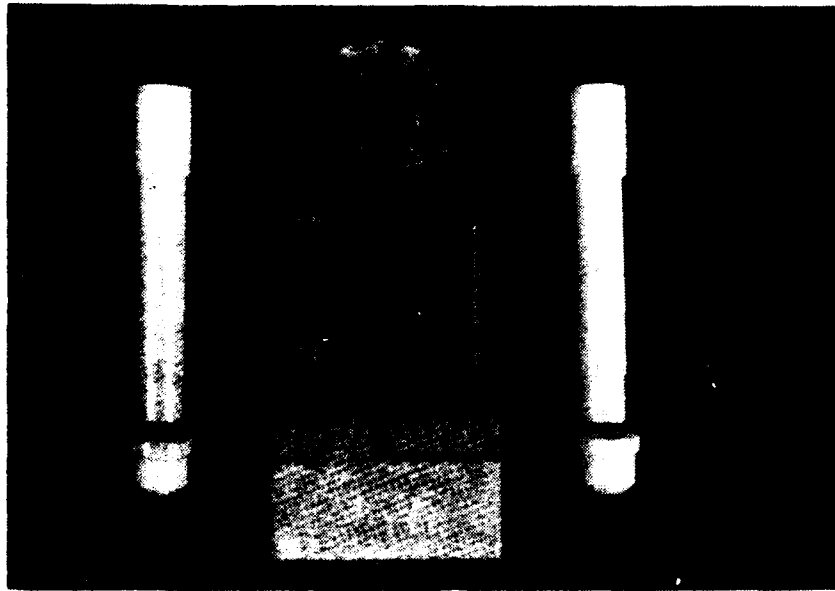


Fig. 1.

Using the 14 MeV neutron activation method, all three elements can be determined non-destructively while allowing the same analyzed crystal to be returned to the furnace for further adjustments of the growth parameters.

It is known that yttrium aluminum garnet /YAG/ can accept impurity ions into its lattice<sup>9,10</sup>. By judiciously selecting and incorporating optically active impurity ions, the properties of solid state materials can be tailored to specific applications.

#### EXPERIMENTAL

##### Procedure for yttrium and aluminum analysis

Samples were core drilled from the YAG boule using a 1/4" diamond core drill in Fig. 1. Sample and standard were vibrated and tamped into individual 1.250" x

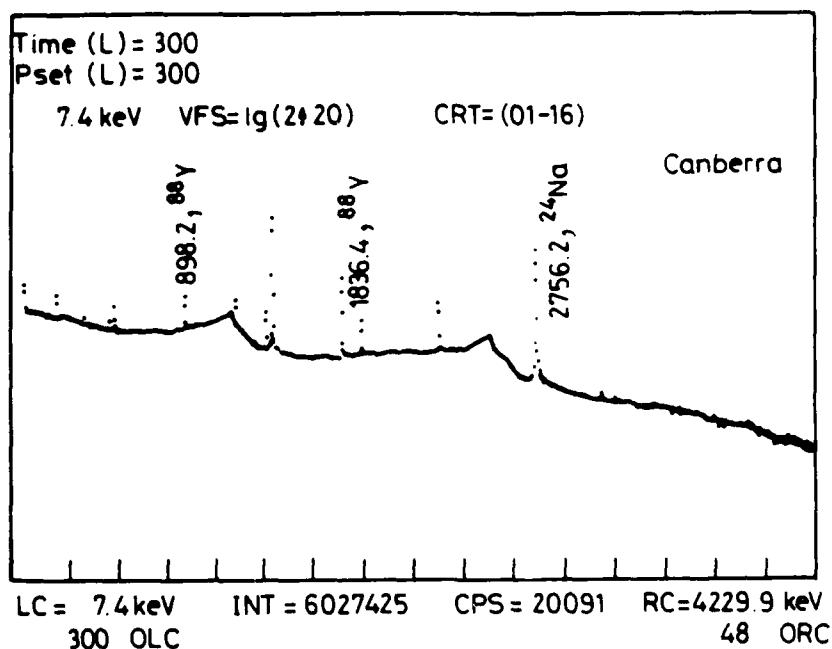


Fig. 2.

0.304" diameter polystyrene capsules. The concentrations of yttrium and aluminum were determined by counting the <sup>88</sup>Y and <sup>24</sup>Na activities produced by a 10 min irradiation using the 14 MeV neutron generator and special irradiation equipment. A Canberra Series-80 multichannel analyzer in conjunction with a high resolution lithium drifted germanium detector was used as the counting system. A typical gamma-ray spectrum of <sup>88</sup>Y, and <sup>24</sup>Na is shown in Fig. 2. By comparing the counts of the sample vs. the standard the quantitative results for yttrium and aluminum were obtained.

#### Procedures for oxygen analyses

The same samples used for yttrium and aluminum analyses were used for the oxygen determination. The oxygen standard used was a Lucite TM cylinder of the same dimensions as the sample. Sample and standard were placed in

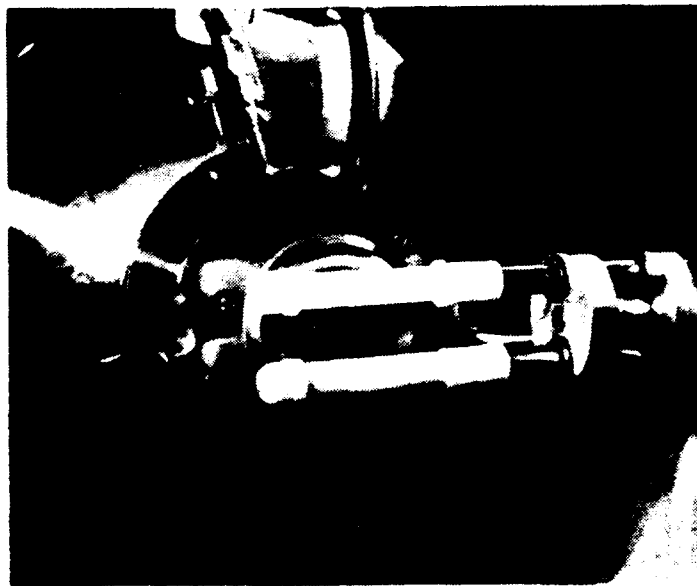


Fig. 3.

separate polyethylene holders for irradiation. In order to obtain a uniform neutron flux both sample and standard underwent three motions; rotation, revolution and traverse<sup>5</sup> as shown in Fig. 3. After irradiation for 20 sec the capsules were ejected from the neutron generator to the counting station /Fig. 4/ through a pair of flexible polyethylene tubes /3/4" diameter/ by means of compressed air. The samples were removed from the polyethylene holders by a specially designed stripper and dropped between two NaI/Tl/ detectors where a photocell was activated to start the counting sequence<sup>9</sup>. Samples and standard were sequentially dropped into the counting position. The complete irradiation and counting process took only 2 min. Such rapid and consistent transport enhances the accuracy of the oxygen analyses.

A 6.1 MeV gamma-ray is produced when  $^{16}\text{O}$  reacts with a neutron to form 7.1 sec  $^{16}\text{N}$ . As can be seen

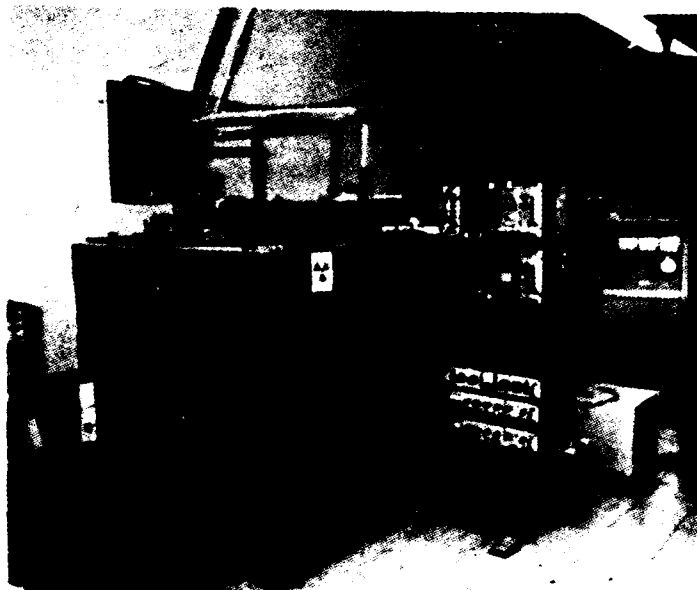


Fig. 4.

from Table 1, the short-lived  $^{16}\text{N}$  can be analyzed as many times as necessary without any interference from yttrium or aluminum.

The counting system consisted of two 5"x5" sodium iodide thallium crystals. A vertical adjustment system ensures that the sample is centered between the crystals to maximize the  $4\pi$  geometry. The electronics of the counting equipment are shown in block diagram of Fig. 5. This system allows an energy window to be set up using a single channel analyzer. The counts in the energy window are used for the analytical data. By using a calibrated pulser any drift in the energy peak can be detected and adjusted.

From the radioactive decay equation it is seen that the rate of disintegration depends on the number of radioactive nuclei present at any time. Integration of the relationship.

TABLE I  
Nuclear reactions

Element	Reaction	Isotopic abundance, %	Cross section, m.barn	Half-life	Energy, MeV
0	$^{16}_0\text{n,p}/^{16}_7\text{N}$	99.76	42	7.13 s	6.13
Al	$^{27}_{13}\text{Al}/\text{n,p}/^{27}_{12}\text{Mg}$	100	65	9.46 m	0.844
	$^{27}_{13}\text{Al}/\text{n},\gamma/^{28}_{13}\text{Al}$	100	230	2.24 m	1.779
	$^{27}_{13}\text{Al}/\text{n},\alpha/^{24}_{11}\text{Na}$	100	118	14.96 h	1.369
Y	$^{89}_{39}\text{Y}/\text{n},2\text{n}/^{88}_{39}\text{Y}$	100	600	106.6 d	1.836
Y	$^{89}_{39}\text{Y}/\text{n,p}/^{89}_{38}\text{Sr}$	100		50.6 d	none
					2.754

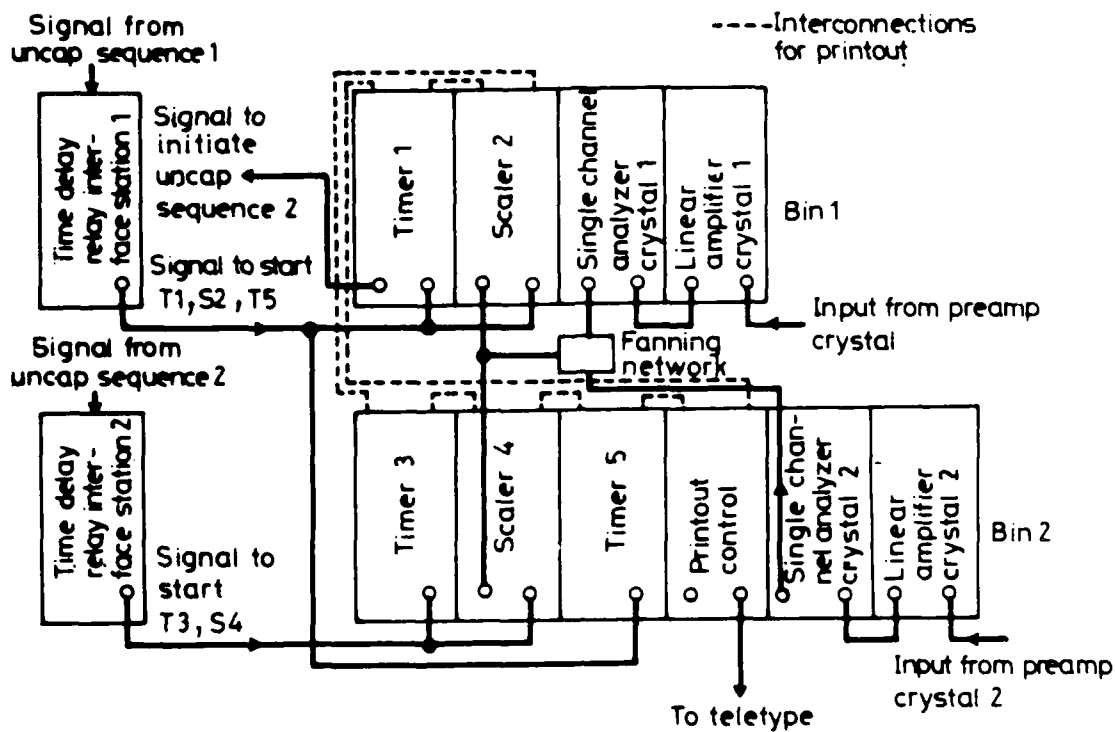


Fig. 5.

$$\frac{dN}{dt} = -\lambda N$$

gives  $N = N_0 e^{-\lambda t}$

where  $N$  = number of radioactive atoms at any time,

$N_0$  = number of radioactive atoms at  $t=0$ ,

$\lambda$  = radioactive desintegration constant =

$$0.693/t_{1/2}$$

$N_0$  is calculated for both the unknown and the standard.

The counts in each counting interval is a definite integral.

$$N_0 = \frac{N}{e^{-\lambda t}}$$

$$\int_{ta}^{tb} e^{-\lambda t} dt = \frac{e^{-\lambda tb}}{-\lambda} - \frac{e^{-\lambda ta}}{-\lambda}$$

$$N_{O/unk/} = \frac{N/unk/}{\frac{e^{-\lambda t/unk/}}{-\lambda} - \frac{e^{-\lambda t/std/}}{-\lambda}}$$

$$N_{O/std/} = \frac{N/std/}{\frac{e^{-\lambda t/std/}}{-\lambda} - \frac{e^{-\lambda t/std/}}{-\lambda}}$$

The ratio of  $N_{O(unk)}/N_{O(std)}$  is the ratio of oxygen in the standard and the unknown, and

$$\text{Oxygen/g/in unk} = \frac{N_{O unk}}{N_{O std}} \times \text{Oxygen/g/in std.}$$

This method of calculation does not depend on the exact reproducibility of time, only the interval of time measurement. The sensitivity using this system is of the order of 1200-1300 counts per milligram of oxygen. The accuracy of this method is considerably enhanced by the removal of the irradiated capsule before counting.

#### Standard

High purity  $Al_2O_3$ , 99.992%, from Adolf Mueller Co., and  $Y_2O_3$ , 99.999% pure, in powder form from Moly Corp., were vibrated and tamped into a weighed polystyrene capsule. A 0.304" diameter x 1.250" long piece of Lucite was used for the oxygen standard. Once it was established that the YAG crystal was stoichiometric it was used as the standard. This simplified the analysis considerably.

### Apparatus

A Kaman A-711 neutron generator was employed in this work. High energy neutrons were produced by the  $^3\text{H}/\text{d}, \text{n}/^4\text{He}$  reaction using a 7 Ci tritium on a titanium target. The neutron flux was  $5 \times 10^9 \text{ n cm}^{-2} \text{ sec}^{-1}$  at the sample position with deuteron beam intensity of 4 mA.

### RESULTS AND DISCUSSION

The quantitative results for yttrium and aluminum were based on the measurement of the activity of the 1.84 MeV gamma peak of  $^{88}\text{Y}$  and the 1.369 MeV gamma peak for aluminum for both sample and standard. A simple comparison between sample count and standard count yielded the quantitative result. The values of these determinations are given in Table 2. As can be seen the relative standard deviation falls well within the quantitative region.

Sources of error resulting from variation of the neutron flux<sup>9</sup> from the generator, counting loss, and decay correction due to multichannel dead time and pulse pile up losses have been taken into consideration and eliminated.

### CONCLUSION

This neutron activation technique is a powerful and useful tool and also unique for the stoichiometric analysis of YAG because it is nondestructive and allows macro-concentrations of three elements to be determined. There is no other analytical method that can accomplish this. In our particular experiments we had several



TABLE 2  
Elemental analysis of YAG /wt %/

Element	Calculated theoretical	*No. of crystal	Content% mean + Std. dev. found
Yttrium	44.93	1	44.62+0.41
		2	44.84+0.39
		3	44.73+0.44
Aluminum	22.73	1	22.61+0.24
		2	22.72+0.21
		3	22.63+0.25
Oxygen**	32.34	1	32.12+0.16
		2	32.22+0.12
		3	32.26+0.21

\*At least 10,000 counts were accumulated for each isotope.

\*\*Each sample analyzed 5 times.

grams of sample. However, the technique is not concentration dependent as are most wet chemical methods.

Recently strong luminescence with a wavelength at 390 nm was discovered in YAG undoped single crystals<sup>10</sup>. Oxygen deficiency is suspected to be a cause of this fluorescence. By using the neutron activation analysis technique in determination of oxygen concentration in YAG can prove invaluable.

A value of  $Y_3Al_5O_{12}$  was obtained from each of the three YAG samples indicating proper premixing before the furnace treatment. It was also found that a pure YAG crystal can be used as a comparative standard in the analysis of newly prepared YAG crystals.

#### REFERENCES

1. J.E. Craig, G. Lee, SPIE 26 Annual International Technical Symposium, Aug. 1981.
2. D. Viechnicki, F. Schmid, J. of Crystal Growth, 26 /1984/ 162.
3. D. Viechnicki, J. Caslavksy, Solid State Formation of Nd:  $Y_3Al_5O_{12}$ , Army Materials Technology Laboratory, MTL, Tech. Report, 78-7.
4. J.S. Abell, J.R. Harris, B. Cockayne, B. Lent, J. Mater. Sci., 9 /1974/ 527.
5. T.M. Rossell, L.D. Hulett, K.J. Northcut, J. Radioanal. Chem., 75 /1982/ 259.
6. W.S. Lyon, Jr., Ed., Guide to Activation Analysis, D. Van Nostrand Co. Inc., Princeton, 1964.
7. W.J. Croft, American Mineralogist, 50, 1634.

BURNS et al.: CHARACTERIZATION OF YTTRIUM ALUMINUM GARNET

8. H.F. Priest, F. Burns, G.L. Priest, Nuclear Instr. and Methods, 50 /1967/ 141.
9. H. Priest, F. Burns, C. Forrest, G. Priest, An Irradiation Transfer and Counting System for Neutron Activation Analysis of Short-lived Components and Inhomogeneous Samples Army Materials Technology Laboratory. MTL TR 70-21. 1970.
10. J. Caslavsky, U.S. Army Materials Technology Laboratory, Watertown Ma., /private communication/.

## DISTRIBUTION LIST

No. of Copies	To	No. of Copies	To
1	Office of the Under Secretary of Defense for Research and Engineering, The Pentagon, Washington, DC 20301	1	Director, U.S. Army Research & Technology Labs, Ames Research Center, Moffet Field, CA 94035
1	ATTN: Mr. J. Persh	1	ATTN: DAVDL-D, Dr. R. Carlson
1	Dr. L. Young	1	DAVOL-AL-D, Dr. I. C. Statler, MS215-1, Aeromechanics Laboratory
1	Mr. K. R. Foster		
2	Commander, U.S. Army Laboratory Command, 2800 Powder Mill Road, Adelphi, MD 20783-1145	1	Commander, U.S. Army Missile Command, Redstone Scientific Information Center, Redstone Arsenal, AL 35898-5241
1	ATTN: SLCIS-IM-TL	1	ATTN: AMSMI-RD-CS-R/ILL Open Lit
1	SLCIS-TD	1	AMSMI-RL, Dr. J. J. Richardson
1	SLCIS-TD-A	1	AMSMI-R, Dr. W. C. McCorkle
1	SLCIS-PA		
1	SLCIS-TP		
2	Commander, Defense Technical Information Center, Cameron Station, Building 5, 5010 Duke Street, Alexandria, VA 22304-6145	1	Commander, U.S. Army Aviation Systems Command, P.O. Box 209 St. Louis, MO 63120
2	ATTN: DTIC-FDAC	1	ATTN: AMSAV-NS, Mr. M. L. Bauccio
1	National Technical Information Service, 5285 Port Royal Road, Springfield, VA 22161	1	Technical Library
1	Director, Defense Advanced Research Projects Agency, 1400 Wilson Boulevard, Arlington, VA 22209	1	Commander, U.S. Army Natick Research, Development, and Engineering Center, Natick, MA 01760
1	ATTN: Dr. P. Parrish	1	ATTN: Technical Library
1	Dr. B. Wilcox	1	Dr. J. A. Sousa
1	Dr. K. Hardmann-Rhyne	1	Dr. R. J. Byrne
		1	Dr. R. Lewis
1	Battelle Columbus Laboratories, Metals and Ceramics Information Center, 505 King Avenue, Columbus, OH 43201	1	Commander, U.S. Army Satellite Communications Agency, Fort Monmouth, NJ 07703
1	ATTN: Mr. W. Duckworth	1	ATTN: Technical Document Center
1	Dr. D. Niesz		
1	Department of the Army, Office of the Assistant Secretary of the Army (ROA), Washington, DC 20310	1	Commander, U.S. Army Science and Technology Center Far East Office, APO San Francisco, CA 96328
1	ATTN: Dr. J. G. Prather, Dep for Sci & Tech	1	ATTN: Terry L. McAfee
1	Dr. J. R. Sculley, SARD		
1	Deputy Chief of Staff, Research, Development, and Acquisition, Headquarters, Department of the Army, Washington, DC 20310	1	Commander, U.S. Army Communications and Electronics Command, Fort Monmouth, NJ 07703
1	ATTN: DAMA-ZE, Mr. C. M. Church	1	ATTN: AMSEL-TDD, Mr. T. A. Pfeiffer, Technical Dir.
1	Commander, U.S. Army Research and Development Office, Chief Research and Development, Washington, DC 20315	1	Director, Electronic Technology and Devices Lab, Fort Monmouth, NJ 07703
1	ATTN: Physical and Engineering Sciences Division	1	ATTN: DELET-D, Dr. C. G. Thornton
1	Commander, Army Research Office, P.O. Box 12211, Research Triangle Park, NC 27709-2211	1	Commander, U.S. Army Tank-Automotive Command, Warren, MI 48090
1	ATTN: Information Processing Office	1	ATTN: Dr. W. Bryzik
1	Dr. G. Mayer	1	D. Rose
1	Dr. J. Hurt	1	AMSTA-RKA
1	Dr. A. Crowson	1	AMSTA-UL, Technical Library
1	Dr. R. Reeber	1	AMSTA-R
1	Dr. R. Shaw	1	AMSTA-NS, Dr. H. H. Dobbs
1	Dr. R. E. Weigle		
		1	Commander, U.S. Army Armament, Munitions and Chemical Command, Dover, NJ 07801
1	Commander, U.S. Army Materiel Command, 5001 Eisenhower Avenue, Alexandria, VA 22333	1	ATTN: Mr. J. Lannon
1	ATTN: AMCLD, Dr. L. Hagan	1	Mr. H. E. Peibly, Jr., PLASTEC, Director
1	AMCDE, Mr. D. L. Griffin	1	Technical Library
1	AMCQA-EQ, Mr. H. L. Light	1	Dr. T. Davidson
1	AMCQA, Mr. S. J. Lorber	1	Dr. B. Ebihara
		1	AMSMC-LC(D), Dr. J. T. Frasier
1	Commander, U.S. Army Electronics Research and Development Command, Fort Monmouth, NJ 07703	1	Commander, U.S. Army Armament, Munitions and Chemical Command, Rock Island, IL 61299
1	ATTN: AMDET-ES, Dr. A. Tauber	1	ATTN: Technical Library
1	Director, Electronics Warfare Laboratory, Fort Monmouth, NJ 07703	1	Commander, U.S. Army Armament, Munitions and Chemical Command, Aberdeen Proving Ground, MD 21010
1	ATTN: AMDEW-D, Mr. M. Adler	1	ATTN: AMSMC-CLN-ST, Mr. S. Shukis
1	Commander, U.S. Army Materiel Systems Analysis Activity, Aberdeen Proving Ground, MD 21005	1	Commander, Aberdeen Proving Ground, MD 21005
1	ATTN: AMXSY-MP, H. Cohen	1	ATTN: AMDAR-CLB-PS, Mr. J. Vervier
1	Commander, U.S. Army Night Vision Electro-Optics Laboratory, Fort Belvoir, VA 22060	1	U.S. Army Corps of Engineers, Construction Engineering Research Lab, P.O. Box 4005, Champaign, IL 61820
1	ATTN: DELNV-S, Mr. P. Travesky	1	ATTN: Dr. Robert Quattrone
1	DELNV-L-D, Dr. R. Buser		
1	DELNV-D, Dr. L. Cameron	1	Commander, U.S. Army Belvoir RD&E Center, Fort Belvoir, VA 22060-5606
		1	ATTN: STRBE-FS, Mr. W. McGovern, Fuel & Wtr Sup Div
		1	AMDME-V, Mr. E. York
		1	AMDME-HS, Dr. K. H. Steinbach
		1	AMDME-ZT, Mr. T. W. Lovelace, Tech Dir
		1	Mr. M. Lepera
1	Commander, Harry Diamond Laboratories, 2800 Powder Mill Road, Adelphi, MD 20783	1	Director, U.S. Army Ballistic Research Laboratory, Aberdeen Proving Ground, MD 21005
1	ATTN: Technical Information Office	1	ATTN: AMDAR-BLT, Dr. A. M. Dietrich
1	AMSLC-RAE	1	AMDAR-BLF, Dr. A. Miller
		1	AMDAR-BLI, Mr. L. Watermeier
		1	AMSMC-BL(A), Dr. R. J. Eichelberger

No. of Copies	To	No. Copies	To
1	Commander, Rock Island Arsenal, Rock Island, IL 61299 ATTN: SARRI-EN	1	Commander, Air Force Armament Center, Eglin Air Force Base, FL 32542 ATTN: Technical Library
1	Director, U.S. Army Industrial Base Engineering Activity, Rock Island, IL 61299 ATTN: AMX18-MT, Mr. G. B. Ney	1	National Aeronautics and Space Administration, Lewis Research Center, 21000 Brookpark Road, Cleveland, OH 44135 ATTN: J. Accurio, USAMRDL
1	Chemical Research and Development Center, Aberdeen Proving Ground, MD 21010 ATTN: AMSMC-CLD(A), Dr. B. Richardson	1	Dr. H. B. Probst, MS 49-1
1	Commander, U.S. Army Test and Evaluation Command, Aberdeen Proving Ground, MD 21005 ATTN: AMSTE-ME	1	Dr. S. Dutta
1	AMSTE-TD, Mr. H. J. Peters	1	National Aeronautics and Space Administration, Washington, DC 20546 ATTN: AFSS-AD, Office of Scientific and Technical Info.
1	Commander, U.S. Army Foreign Science and Technology Center, 220 7th Street, N.E., Charlottesville, VA 22901 ATTN: Military Tech	1	National Aeronautics and Space Administration, Langley Research Center, Hampton, VA 23665 ATTN: Mr. J. Buckley, MS 387
1	Mr. J. Crider	1	Dr. J. Heyman, MS 231
1	Ms. P. Durrer	1	Mr. R. L. Long, MS 266
1	Mr. P. Greenbaum	1	Commander, White Sands Missile Range, Electronic Warfare Laboratory, OMEW, ERADCOM, White Sands, NM 88002 ATTN: Mr. Thomas Reader, AMSEL-WLM-ME
1	Chief, Benet Weapons Laboratory, Watervliet, NY 12189 ATTN: AMDAR-LCB-TL	1	Department of Energy, Division of Transportation, 20 Massachusetts Avenue, N.W., Washington, DC 20545 ATTN: Mr. G. Thur
1	Dr. G. D'Andrea	1	Dr. R. J. Gottschall, ER-131, GTN
1	AMDAR-LCB, Dr. F. Sautter	1	Department of Transportation, 400 Seventh Street, S.W., Washington, DC 20590 ATTN: Mr. M. Lauriente
1	Director, Eustis Directorate, U.S. Army Mobility Research and Development Laboratory, Fort Eustis, VA 23604 ATTN: SAVDL-E-MOS (AMCCOM)	1	Mechanical Properties Data Center, Belfour Stulen Inc., 13917 W. Bay Shore Drive, Traverse City, MI 49684
1	Commander, U.S. Army Engineer Waterways Experiment Station, Vicksburg, MS 39180 ATTN: Research Center Library	1	National Bureau of Standards, Washington, DC 20234 ATTN: E. S. Etz, Bldg. 222, Rm A-121
1	Project Manager, Munitions Production Base, Modernization and Expansion, Dover, NJ 07801 ATTN: AMCPM-PBM-P	1	D. L. Hunston, Bldg. 224, Rm A-209
1	Technical Director, Human Engineering Laboratories, Aberdeen Proving Ground, MD 21005 ATTN: Technical Reports Office	1	Dr. D. H. Reneker, Dep. Dir., Ctr for Mat'l's Sci.
1	AMXHE-D, Dr. J. D. Weisz	1	Dr. Lyle Schwartz
1	Chief of Naval Research Arlington, VA 22217 ATTN: Code 471	1	Dr. Stephen Hsu
1	Dr. A. Diness	1	Dr. Allan Draggoo
1	Dr. R. Pohanka	1	U.S. Bureau of Mines, Mineral Resources Technology, 2401 E. Street, N.W., Washington, DC 20241 ATTN: Mr. M. A. Schwartz
1	Naval Research Laboratory, Washington, DC 20375 ATTN: Code 5830	1	National Bureau of Standards, Gaithersburgh, MD 20760 ATTN: Dr. S. Wiederhorn
1	Headquarters, Naval Air Systems Command, Washington, DC 20360 ATTN: Code 5203	1	Dr. J. B. Wachtman
1	Headquarters, Naval Sea Systems Command, 1941 Jefferson Davis Highway, Arlington, VA 22376 ATTN: Code 035	1	Dr. N. Tighe
1	Headquarters, Naval Electronics Systems Command, Washington, DC 20360 ATTN: Code 504	1	National Research Council, National Materials Advisory Board, 2101 Constitution Avenue, Washington, DC 20418 ATTN: Dr. K. Zwilsky
1	Commander, Naval Ordnance Station, Louisville, KY 40214 ATTN: Code 85	1	D. Groves
1	Commander, Naval Material Industrial Resources Office, Building 537-2, Philadelphia Naval Base, Philadelphia, PA 19112 ATTN: Technical Director	1	R. M. Spriggs
1	Commander, Naval Weapons Center, China Lake CA 93555 ATTN: Mr. F. Markarian	1	J. Lane
1	Commander, U.S. Army Wright Aeronautical Labs, Wright-Patterson Air Force Base, OH 45433 ATTN: Dr. M. Tallan	1	National Science Foundation, Materials Division, 1800 G Street, N.W., Washington, DC 20006 ATTN: Dr. L. Toth
1	Dr. H. Graham	1	Dr. J. Hurt
1	Dr. R. Ruh	1	AiResearch Manufacturing Company, AiResearch Casting Company, 2525 West 190th Street, Torrance, CA 90505 ATTN: Mr. K. Styhr
1	Aero Propulsion Labs, Mr. R. Marsh	1	AVCO Corporation, Applied Technology Division, Lowell Industrial Park, Lowell, MA 01887 ATTN: Dr. T. Vasilos
1	Dr. H. M. Burte	1	Case Western Reserve University, Department of Metallurgy, Cleveland, OH 60605 ATTN: Prof. A. H. Meuer
1	AFWL/MLLP, Mr. D. Forney	1	Defence Research Establishment Pacific, FMO, Victoria, B.C., VOS 1B0, Canada ATTN: R. D. Barer
1	AFML/MLLM, Mr. H. L. Gegel	1	European Research Office, 223 Old Marylebone Road, London, NW1 - 5th, England ATTN: Dr. I. Ahmad
1	AFSC/MLLM, Dr. A. Katz		

No. of Copies	To	No. of Copies	To
1	Ford Motor Company, Turbine Research Department, 20000 Rotunda Drive, Dearborn, MI 48121	1	Rensselaer Polytechnic Institute, Department of Materials Engineering, Troy, NY 12181
1	ATTN: Mr. A. F. McLean	1	ATTN: R. J. Diefendorf
1	Mr. J. A. Mangels		
1	Ford Motor Company, P.O. Box 2053, Dearborn, MI 48121	1	Oak Ridge National Laboratory, P.O. Box X
1	ATTN: Dr. D. Compton, Vice President, Research	1	Oak Ridge, TN 37830
		1	ATTN: P. F. Becher
		1	V. J. Tennery
		1	R. Johnson
1	General Electric Company, Research and Development Center, Box 8, Schenectady, NY 12345	1	Sandia Laboratories, Albuquerque, NM 87185
1	ATTN: Dr. R. J. Charles	1	ATTN: Dr. F. Gerstle, Div 5814
1	Dr. C. D. Greskovich		
1	Dr. S. Prochazka		
1	Georgia Institute of Technology, EES, Atlanta, GA 30332	1	The John Hopkins University, Department of Civil Engineering/
1	ATTN: Mr. J. D. Walton	1	Materials Science and Engineering, Baltimore, MD 28218
		1	ATTN: Dr. R. E. Green, Jr.
1	GTE Sylvania, Waltham Research Center, 40 Sylvania Road, Waltham, MA 02154	1	Director, Office of Science and Technology Policy, Old Executive Office Building, Washington, DC 20223
1	ATTN: Dr. W. H. Rhodes		
1	Martin Marietta Laboratories, 1450 South Rolling Road, Baltimore, MD 21227	1	Subcommittee on Science, 2319 Rayburn House Office Building, Washington, DC 20515
1	ATTN: Dr. J. Venables	1	ATTN: Mr. P. C. Maxwell
1	Massachusetts Institute of Technology, Department of Metallurgy and Materials Science, Cambridge, MA 02139	1	Aerospace Corporation, Materials Science Laboratory, 2350 East El Segundo Boulevard, El Segundo, CA 90245
1	ATTN: Prof. R. L. Coble	1	ATTN: Dr. L. R. McCreight
1	Prof. H. K. Bowen		
1	Prof. W. D. Kingery	1	IBM Corporation, Thomas B. Watson Research Center, Yorktown Heights, NY 10598
1	Prof. J. Vander Sande	1	ATTN: Dr. G. Onoda
1	Materials Research Laboratories, P.O. Box 50, Ascot Vale, VIC 3032, Australia	1	Corning Glass Works, Research and Development Division, Corning, NY 14830
1	ATTN: Dr. C. W. Weaver	1	ATTN: Dr. W. R. Prindle
1	Midwest Research Institute, 425 Volker Boulevard, Kansas City, MO 64110	1	3M Company, New Products Department, 219-01-01, 3M Center, St. Paul, MN 55144
1	ATTN: Mr. G. W. Gross, Head, Physics Station	1	ATTN: R. E. Richards
1	Pennsylvania State University, Materials Research Laboratory, Materials Science Department, University Park, PA 16802	1	Technology Strategies, Inc., 10722 Shingle Oak Ct., Burke, VA 22015
1	ATTN: Prof. R. Roy	1	ATTN: Dr. E. C. Van Reuth
1	Prof. R. E. Newnham		
1	Prof. R. E. Tressler	1	Rutgers University, Center for Ceramics, Rm A274, P.O. Box 909, Piscataway, NJ 08854
1	Dr. C. Pantano	1	ATTN: Prof. J. B. Wachtman, Jr., Director
1	Mr. C. O. Ruud		
1	State University of New York at Albany, Department of Physics, Albany, NY 12222	1	Syracuse University, 304 Administration Building, Syracuse, NY 13210
1	ATTN: Prof. W. A. Lanford	1	ATTN: Dr. V. Weiss
1	State University of New York at Stony Brook, Department of Materials Science, Long Island, NY 11790	1	Lehigh University, Materials Research Center #32, Bethlehem, PA 18015
1	ATTN: Prof. F. F. Y. Wang	1	ATTN: Dr. D. M. Smyth
1	Stanford Research International, 333 Ravenswood Avenue, Menlo Park, CA 94025d	1	Alfred University, New York State College of Ceramics, Alfred, NY 14802
1	ATTN: Dr. P. Jorgensen	1	ATTN: Dr. R. L. Snyder
1	Dr. D. Rowcliffe		
1	United Technologies Research Center, East Hartford, CT 06108	1	University of California, Center for Advanced Materials, D58, Hildebrand Hall, Berkeley, CA 94720
1	ATTN: Dr. J. Brennan	1	ATTN: Prof. G. Somorjai
1	Dr. K. Prewo		
1	University of California, Lawrence Livermore Laboratory, P.O. Box 808, Livermore, CA 94550	1	Boeing Aerospace Company, 11029 Southeast 291, Auburn, MA 98002
1	ATTN: Mr. R. Landingham	1	ATTN: W. E. Strobel
1	Dr. C. F. Cline		
1	Dr. J. Birch Holt	1	University of California, Materials Science and Mineral Engineering, Heart Mining Building, Rm 284, Berkeley, CA 94720
1	University of Florida, Department of Materials Science and Engineering, Gainesville, FL 32611	1	ATTN: Prof. G. Thomas
1	ATTN: Dr. L. Hench		
1	University of Washington, Ceramic Engineering Division, FB-10, Seattle, WA 98195	2	Director, U.S. Army Materials Technology Laboratory, Watertown, MA 02172-0001
1	ATTN: Prof. R. Bradt	3	ATTN: SLCMT-IML
			Authors
1	Westinghouse Electric Corporation, Research Laboratories, Pittsburgh, PA 15235		
1	ATTN: Dr. R. J. Bratton		
1	Battelle Pacific Northwest Lab, NDT Section, Richland, WA 99353		
1	ATTN: Mr. A. Birks, Associate Manager		

U.S. Army Materials Technology Laboratory  
Watertown, Massachusetts 02172-0001  
CHARACTERIZATION OF YTTRIUM ALUMINUM  
GARNET /YAG/ USING 14 MeV NEUTRON  
ACTIVATION - Forrest Burns, Michaeline F. Chen  
and Jaroslav L. Caslavsky

Technical Report MIL TR 87-62, December 1987, 15 pp -  
Illustrations, D/A Project IL161102.H42

AD  
UNCLASSIFIED  
UNLIMITED DISTRIBUTION  
Key Words  
14 MeV neutron  
Gamma-ray spectroscopy  
Oxygen analysis

A non-destructive method for the analysis of yttrium, aluminum and oxygen in Yttrium Aluminum Garnet (YAG) by 14 MeV neutron activation analysis has been developed and evaluated. The fast neutron induces primarily  $n,2n$ ,  $n,p$ , or  $n,\alpha$  reactions with yttrium, aluminum, and oxygen to produce isotopes with measurable characteristic gamma-ray spectra. Concentrations of Y, Al and O were determined on the basis of calibrated nuclear decay emission spectra recorded on selected standards. After each analysis the crystal was returned to the furnace for further adjustments of the growth parameters, and any resultant elemental variations were observed in the next analysis. The accuracy of this method was determined by repeated analyses with various "YAG" crystals from a series of separate furnace runs.

U.S. Army Materials Technology Laboratory  
Watertown, Massachusetts 02172-0001  
CHARACTERIZATION OF YTTRIUM ALUMINUM  
GARNET /YAG/ USING 14 MeV NEUTRON  
ACTIVATION - Forrest Burns, Michaeline F. Chen  
and Jaroslav L. Caslavsky

Technical Report MIL TR 87-62, December 1987, 15 pp -  
Illustrations, D/A Project IL161102.H42

AD  
UNCLASSIFIED  
UNLIMITED DISTRIBUTION  
Key Words  
14 MeV neutron  
Gamma-ray spectroscopy  
Oxygen analysis

A non-destructive method for the analysis of yttrium, aluminum and oxygen in Yttrium Aluminum Garnet (YAG) by 14 MeV neutron activation analysis has been developed and evaluated. The fast neutron induces primarily  $n,2n$ ,  $n,p$ , or  $n,\alpha$  reactions with yttrium, aluminum, and oxygen to produce isotopes with measurable characteristic gamma-ray spectra. Concentrations of Y, Al and O were determined on the basis of calibrated nuclear decay emission spectra recorded on selected standards. After each analysis the crystal was returned to the furnace for further adjustments of the growth parameters, and any resultant elemental variations were observed in the next analysis. The accuracy of this method was determined by repeated analyses with various "YAG" crystals from a series of separate furnace runs.

U.S. Army Materials Technology Laboratory  
Watertown, Massachusetts 02172-0001  
CHARACTERIZATION OF YTTRIUM ALUMINUM  
GARNET /YAG/ USING 14 MeV NEUTRON  
ACTIVATION - Forrest Burns, Michaeline F. Chen  
and Jaroslav L. Caslavsky

Technical Report MIL TR 87-62, December 1987, 15 pp -  
Illustrations, D/A Project IL161102.H42

AD  
UNCLASSIFIED  
UNLIMITED DISTRIBUTION  
Key Words  
14 MeV neutron  
Gamma-ray spectroscopy  
Oxygen analysis

A non-destructive method for the analysis of yttrium, aluminum and oxygen in Yttrium Aluminum Garnet (YAG) by 14 MeV neutron activation analysis has been developed and evaluated. The fast neutron induces primarily  $n,2n$ ,  $n,p$ , or  $n,\alpha$  reactions with yttrium, aluminum, and oxygen to produce isotopes with measurable characteristic gamma-ray spectra. Concentrations of Y, Al and O were determined on the basis of calibrated nuclear decay emission spectra recorded on selected standards. After each analysis the crystal was returned to the furnace for further adjustments of the growth parameters, and any resultant elemental variations were observed in the next analysis. The accuracy of this method was determined by repeated analyses with various "YAG" crystals from a series of separate furnace runs.

U.S. Army Materials Technology Laboratory  
Watertown, Massachusetts 02172-0001  
CHARACTERIZATION OF YTTRIUM ALUMINUM  
GARNET /YAG/ USING 14 MeV NEUTRON  
ACTIVATION - Forrest Burns, Michaeline F. Chen  
and Jaroslav L. Caslavsky

Technical Report MIL TR 87-62, December 1987, 15 pp -  
Illustrations, D/A Project IL161102.H42

AD  
UNCLASSIFIED  
UNLIMITED DISTRIBUTION  
Key Words  
14 MeV neutron  
Gamma-ray spectroscopy  
Oxygen analysis

A non-destructive method for the analysis of yttrium, aluminum and oxygen in Yttrium Aluminum Garnet (YAG) by 14 MeV neutron activation analysis has been developed and evaluated. The fast neutron induces primarily  $n,2n$ ,  $n,p$ , or  $n,\alpha$  reactions with yttrium, aluminum, and oxygen to produce isotopes with measurable characteristic gamma-ray spectra. Concentrations of Y, Al and O were determined on the basis of calibrated nuclear decay emission spectra recorded on selected standards. After each analysis the crystal was returned to the furnace for further adjustments of the growth parameters, and any resultant elemental variations were observed in the next analysis. The accuracy of this method was determined by repeated analyses with various "YAG" crystals from a series of separate furnace runs.

END  
DATE  
FILMED  
5-88  
DTIC