

UNCLASSIFIED

<b>AD NUMBER</b>
ADB204477
<b>NEW LIMITATION CHANGE</b>
<b>TO</b> Approved for public release, distribution unlimited
<b>FROM</b> Distribution authorized to DoD only; Specific Authority; Oct 95. Other requests shall be referred to Commander, U.S. Army Medical Research and Materiel Command, Attn: MCMR-RMI-S, Ft. Detrick, MD 21702-5012.
<b>AUTHORITY</b>
USAMRMC ltr dtd 21 Jan 2000

THIS PAGE IS UNCLASSIFIED

AD \_\_\_\_\_

CONTRACT NUMBER: DAMD17-95-C-5041


TITLE: Non-Invasive Deep Tissue Oxygenation Sensor

PRINCIPAL INVESTIGATOR: Stephen D. Walker

CONTRACTING ORGANIZATION: Vixel Corporation  
Broomfield, Colorado 80021

REPORT DATE: October 1995

TYPE OF REPORT: Final, Phase I

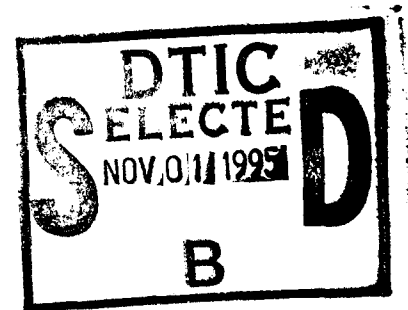

PREPARED FOR: U.S. Army Medical Research and Materiel Command  
Fort Detrick, Maryland 21702-5012

DISTRIBUTION STATEMENT: Distribution authorized to DoD Components only, Specific Authority. Other requests shall be referred to the Commander, U.S. Army Medical Research and Materiel Command, ATTN: MCMR-RMI-S, Fort Detrick, MD 21702-5012 (oct95)

The views, opinions and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy or decision unless so designated by other documentation.

19951027 106

DTIC QUALITY INSPECTED 8

# REPORT DOCUMENTATION PAGE

Form Approved  
OMB No. 0704-0188

Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503.

1. AGENCY USE ONLY (Leave blank)	2. REPORT DATE October 1995	3. REPORT TYPE AND DATES COVERED Final, Phase I 15 Mar 95 - 14 Sep 95	
4. TITLE AND SUBTITLE Non-Invasive Deep Tissue Oxygenation Sensor		5. FUNDING NUMBERS DAMD17-95-C-5041	
6. AUTHOR(S) Stephen D. Walker			
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) Vixel Corporation Broomfield, Colorado 80021		8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING / MONITORING AGENCY NAME(S) AND ADDRESS(ES) U.S. Army Medical Research and Materiel Command Fort Detrick, Maryland 21702-5012		10. SPONSORING / MONITORING AGENCY REPORT NUMBER	
11. SUPPLEMENTARY NOTES			
12a. DISTRIBUTION / AVAILABILITY STATEMENT Distribution authorized to DoD components only, Specific Authority. Other requests shall be referred to the Commander, U.S. Army Medical Research and Materiel Command, ATTN: MCMR-RMI-S, Fort Detrick, MD 21702-5012 <b>Oct 95</b>		12b. DISTRIBUTION CODE	
13. ABSTRACT (Maximum 200 words)  This Small Business Innovation Research Phase I project proved the feasibility of a compact and low cost instrument to measure tissue oxygenation based on near-infrared (NIR) spectroscopy. NIR light in the .7-1.1 μm range penetrates human tissue to a depth of several centimeters, and is attenuated during transmission by oxyhemoglobin, deoxyhemoglobin, melanin, water and cytochrome aa3, each with its own characteristic absorption spectrum. NIR absorption spectroscopy offers a number of advantages in comparison with other oxygen measurement modalities: it is non-invasive, uses non-ionizing radiation, offers high spatial and temporal resolution, and supplies new types of metabolic and functional information for low cost, rapid, bedside diagnostic procedures. The Phase I project determined the most practical means to implement a hand-held NIR spectroscopic oxygenation monitor, specified the components needed, defined the measurement method and assessed the practicality of the instrument for Phase II development.			
14. SUBJECT TERMS  Tissue oxygenation, Near Infrared Spectroscopy Oxyhemoglobin, Deoxyhemoglobin, Oxygen, Cytochrome aa3		15. NUMBER OF PAGES 12	
		16. PRICE CODE	
17. SECURITY CLASSIFICATION OF REPORT Unclassified	18. SECURITY CLASSIFICATION OF THIS PAGE Unclassified	19. SECURITY CLASSIFICATION OF ABSTRACT Unclassified	20. LIMITATION OF ABSTRACT limited

## GENERAL INSTRUCTIONS FOR COMPLETING SF 298

The Report Documentation Page (RDP) is used in announcing and cataloging reports. It is important that this information be consistent with the rest of the report, particularly the cover and title page. Instructions for filling in each block of the form follow. It is important to *stay within the lines* to meet *optical scanning requirements*.

**Block 1. Agency Use Only (Leave blank).**

**Block 2. Report Date.** Full publication date including day, month, and year, if available (e.g. 1 Jan 88). Must cite at least the year.

**Block 3. Type of Report and Dates Covered.** State whether report is interim, final, etc. If applicable, enter inclusive report dates (e.g. 10 Jun 87 - 30 Jun 88).

**Block 4. Title and Subtitle.** A title is taken from the part of the report that provides the most meaningful and complete information. When a report is prepared in more than one volume, repeat the primary title, add volume number, and include subtitle for the specific volume. On classified documents enter the title classification in parentheses.

**Block 5. Funding Numbers.** To include contract and grant numbers; may include program element number(s), project number(s), task number(s), and work unit number(s). Use the following labels:

C - Contract	PR - Project
G - Grant	TA - Task
PE - Program Element	WU - Work Unit Accession No.

**Block 6. Author(s).** Name(s) of person(s) responsible for writing the report, performing the research, or credited with the content of the report. If editor or compiler, this should follow the name(s).

**Block 7. Performing Organization Name(s) and Address(es).** Self-explanatory.

**Block 8. Performing Organization Report Number.** Enter the unique alphanumeric report number(s) assigned by the organization performing the report.

**Block 9. Sponsoring/Monitoring Agency Name(s) and Address(es).** Self-explanatory.

**Block 10. Sponsoring/Monitoring Agency Report Number. (If known)**

**Block 11. Supplementary Notes.** Enter information not included elsewhere such as: Prepared in cooperation with...; Trans. of...; To be published in.... When a report is revised, include a statement whether the new report supersedes or supplements the older report.

**Block 12a. Distribution/Availability Statement.** Denotes public availability or limitations. Cite any availability to the public. Enter additional limitations or special markings in all capitals (e.g. NOFORN, REL, ITAR).

**DOD** - See DoDD 5230.24, "Distribution Statements on Technical Documents."

**DOE** - See authorities.

**NASA** - See Handbook NHB 2200.2.

**NTIS** - Leave blank.

**Block 12b. Distribution Code.**

**DOD** - Leave blank.

**DOE** - Enter DOE distribution categories from the Standard Distribution for Unclassified Scientific and Technical Reports.

**NASA** - Leave blank.

**NTIS** - Leave blank.

**Block 13. Abstract.** Include a brief (*Maximum 200 words*) factual summary of the most significant information contained in the report.

**Block 14. Subject Terms.** Keywords or phrases identifying major subjects in the report.

**Block 15. Number of Pages.** Enter the total number of pages.

**Block 16. Price Code.** Enter appropriate price code (*NTIS only*).

**Blocks 17. - 19. Security Classifications.** Self-explanatory. Enter U.S. Security Classification in accordance with U.S. Security Regulations (i.e., UNCLASSIFIED). If form contains classified information, stamp classification on the top and bottom of the page.

**Block 20. Limitation of Abstract.** This block must be completed to assign a limitation to the abstract. Enter either UL (unlimited) or SAR (same as report). An entry in this block is necessary if the abstract is to be limited. If blank, the abstract is assumed to be unlimited.

FOREWORD

Opinions, interpretations, conclusions and recommendations are those of the author and are not necessarily endorsed by the US Army.

Where copyrighted material is quoted, permission has been obtained to use such material.

Where material from documents designated for limited distribution is quoted, permission has been obtained to use the material.

*sw* Citations of commercial organizations and trade names in this report do not constitute an official Department of Army endorsement or approval of the products or services of these organizations.

In conducting research using animals, the investigator(s) adhered to the "Guide for the Care and Use of Laboratory Animals," prepared by the Committee on Care and Use of Laboratory Animals of the Institute of Laboratory Resources, National Research Council (NIH Publication No. 86-23, Revised 1985).

For the protection of human subjects, the investigator(s) adhered to policies of applicable Federal Law 45 CFR 46.

In conducting research utilizing recombinant DNA technology, the investigator(s) adhered to current guidelines promulgated by the National Institutes of Health.

In the conduct of research utilizing recombinant DNA, the investigator(s) adhered to the NIH Guidelines for Research Involving Recombinant DNA Molecules.

In the conduct of research involving hazardous organisms, the investigator(s) adhered to the CDC-NIH Guide for Biosafety in Microbiological and Biomedical Laboratories.

<b>Accession For</b>	
NTIS GRA&I	<input type="checkbox"/>
DTIC TAB	<input checked="" type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification _____	
By _____	
Distribution/ _____	
<b>Availability Codes</b>	
<b>Dist</b>	<b>Avail and/or Special</b>
E-4	

*swallen* 9/27/95  
PI - Signature Date

## Table of Contents


INTRODUCTION	1
BODY	2
CONCLUSIONS	6
REFERENCES	7
APPENDIX	8

# Non Invasive Deep Tissue Oxygenation Sensor

## Phase I Final Report

### INTRODUCTION

The aim of this Phase I project was to assess the practicality of using laser diodes specifically for tissue oxygenation near infrared spectroscopy. Near Infrared (NIR) light penetrates human tissue to a depth of several centimeters. Light is attenuated during transmission by dominant absorbers, each of which absorbs at a specific wavelength in the NIR region. These dominant absorbers are involved with cellular metabolism and tissue oxygenation. The amount of light absorbed is proportional to the concentration of the absorber. Figure 1 illustrates the basic principle.

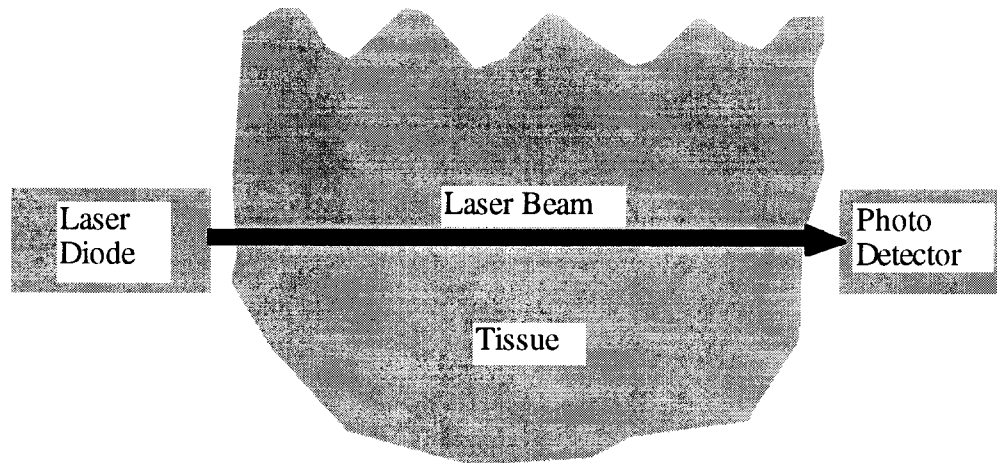

**Figure 1.** NIR laser light is directed through tissue onto a photodetector. For a specific wavelength, the amount of light absorbed is proportional to the amount of absorber in the optical path.

The project developed a multiple discrete wavelength approach to perform NIR spectroscopy. First, we determined the absorption spectra of absorbers that are involved in tissue oxygenation. Potential absorbers were oxyhemoglobin, deoxyhemoglobin, metoxyhemoglobin, carboxyhemoglobin, cytochrome aa3, melanin, and water. By using two discrete lasers with emission wavelengths corresponding to absorptive bands of two different absorbers, the amounts of the two different absorbers can be measured. Amounts of four different absorbers can be determined with four different emission wavelengths.

Conventional edge-emitting laser diode are prevalent in commercial products such as compact disk players, desktop laser printers, and handheld barcode scanners. Vertical-Cavity Surface-Emitting Lasers (VCSELs) are an emerging form of laser diode which have intrinsic characteristics far superior to edge-emitters in terms of manufacturing, compact low-cost packaging, and spectral stability. While edge emitters are several hundred microns long and emit in the plane of the chip, VCSELs are circular, about 10  $\mu\text{m}$  in diameter, and emit perpendicular to the plane of the chip, providing very easy access to the beam. VCSELs are manufactured by standard GaAs processing techniques and the emitted beams are circular, low-divergence and aberration free. The fabricated VCSEL is a nearly-planar structure which can be tested while still on the wafer by completely automated means, thereby greatly reducing the cost of the laser. Required operating currents are only 3-5 milliamps, compared to 30-50 milliamps in a compact-disk laser. Thus the power required for the VCSEL-based sensor will be much lower than an edge-emitter based one, which is an important consideration for hand-held, battery-operated instruments utilizing multiple lasers.

# Non Invasive Deep Tissue Oxygenation Sensor

## Phase I Final Report

### BODY

#### Phase I Technical Objectives

1. Analyze spectra of primary absorbers.
2. Determine the optimum set of wavelengths for monitoring oxygenation.
3. Define the measurement method for measuring oxygenation.
4. Assess the overall practicality of the proposed instrument.
5. Report results and develop a Phase II plan to develop a prototype to be commercialized in a Phase III program.


#### Analyze spectra of primary absorbers

1. Oxyhemoglobin. Oxyhemoglobin is oxygenated hemoglobin and equivalent to arterial blood. Knowledge of oxyhemoglobin concentration is essential for the measurement of tissue oxygenation because it is the primary oxygen source. Low concentrations of oxyhemoglobin are caused by insufficient oxygenation of the blood in the pulmonary system. This can be the result of insufficient cardiac output or improper ventilation. Localized high concentrations are an indicator of hematomas or necrotic tissue because cellular metabolism is not taking place. Oxyhemoglobin has a broad absorption band from 600-810 nm.
2. Deoxyhemoglobin. Deoxyhemoglobin is the reduced hemoglobin molecule and equivalent to venous blood. Low concentrations of deoxyhemoglobin indicate cellular metabolism is not taking place or the hemoglobin molecule is immediately combining with interstitial oxygen. Deoxyhemoglobin has a broad absorption band from 800-1000 nm.
3. Cytochrome aa3. Cytochrome aa3 is a mitochondrial enzyme whose presence indicates cellular metabolism is taking place. It has a moderate absorption band from 830-870 nm.
4. Melanin. Melanin is a skin pigment and essentially an interfering absorber. Different skin colors absorb NIR light to different degrees. It has a broad absorption band from 600-1000 nm.
5. Water. Water is everywhere and has several very high absorption bands. These bands are specified in the HITRAN 92 database.

Potential interfering NIR absorbers found in tissue were also identified and are shown in Figure 2.




**Non Invasive Deep Tissue Oxygenation Sensor  
Phase I Final Report**

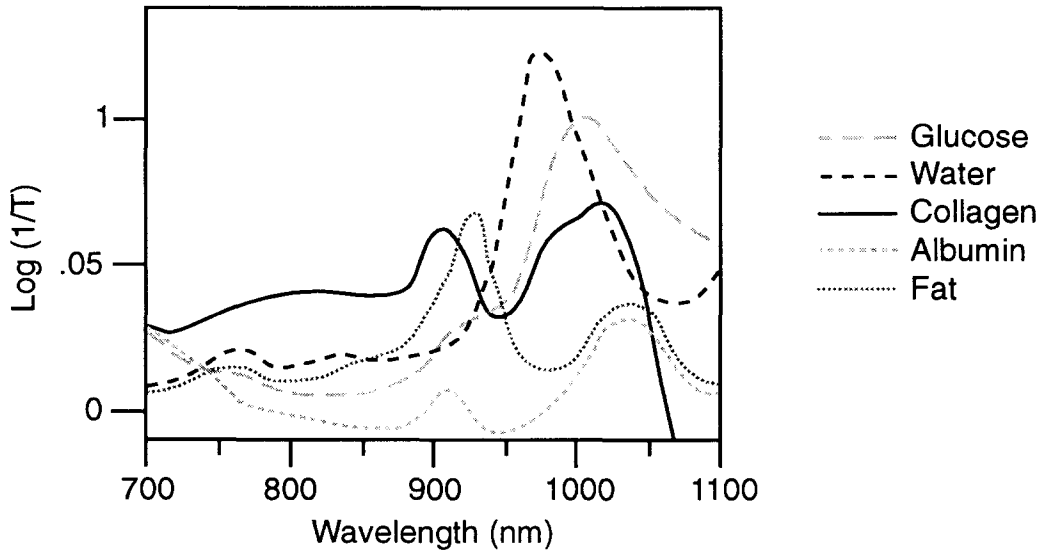

**Figure 2.** Near IR absorption spectra of glucose and interfering absorbers.

**Determine the optimum set of wavelengths for monitoring oxygenation**

Neuro Shell 2, a commercial neural network program, was used to develop a software model of the proposed spectrometer. The objective of the simulation was to identify five laser diode wavelengths which when used as tissue oxygenation spectrometer light sources could quantify the targeted absorbers. Fabricating laser diodes at a specified wavelength is very expensive and the advantage of this simulator is obvious. This simulation ran for 22 hours on a 25 MHz 486 computer and selected 5 wavelengths shown in Figure 3.

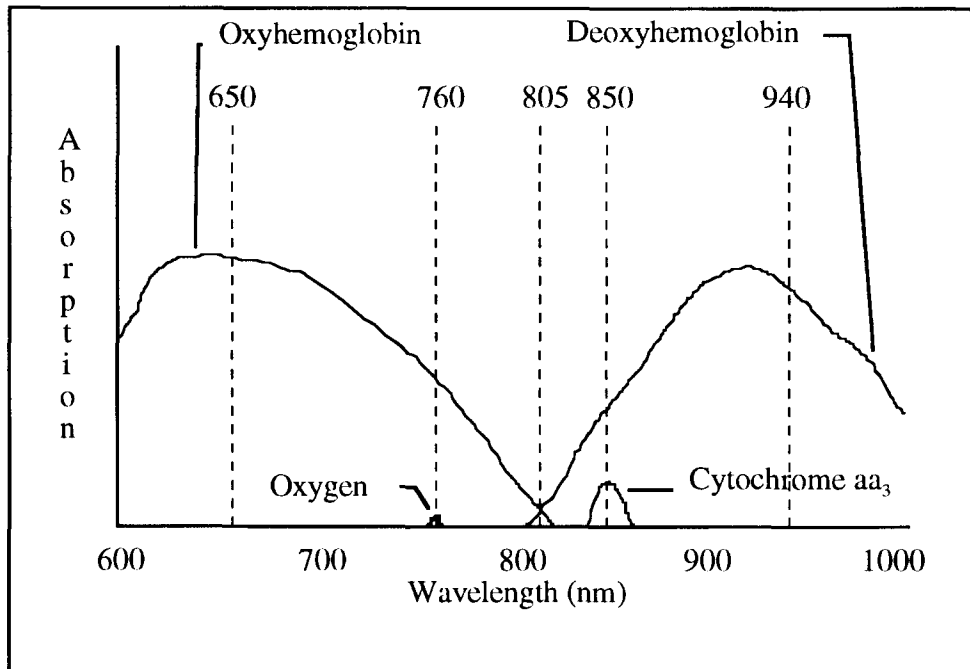

**Figure 3.** NIR absorption spectra of oxyhemoglobin, deoxyhemoglobin, oxygen and cytochrome aa3. Laser diode wavelengths selected by simulation are shown as dotted lines.

# Non Invasive Deep Tissue Oxygenation Sensor

## Phase I Final Report

The simulator was first trained to recognize the primary absorbers and oxygen by giving it their individual spectra as inputs. It was then trained to recognize combinations of 2 spectra. Finally, it was allowed to choose only five wavelengths in the spectra as inputs. The simulator varied these five wavelengths until it found the highest recognition rate. Melanin and water spectra were not simulated because they are not directly associated with cellular metabolism.

### Define the measurement method for measuring oxygenation

A prototype spectrometer was constructed and tested on the oxygen wavelength to validate the WMS method. The block diagram is shown in Figure 4. A thermoelectric cooler is used to adjust the VCSEL temperature for an emitted wavelength of 760.3 nm. Wavelength modulation of 0.2 nm is achieved by applying a sawtooth drive current to the VCSEL. The laser beam is directed onto a beamsplitter where 50% of the energy is reflected through a reference cell containing pure nitrogen and onto a detector whose output is the inverting input to a differential amplifier. The 50% of the energy which passed through the beamsplitter is reflected off a mirror through the sample chamber containing pure oxygen and onto a detector whose output is the non-inverting input to the differential amplifier. The output of the amplifier is the oxygen spectra.

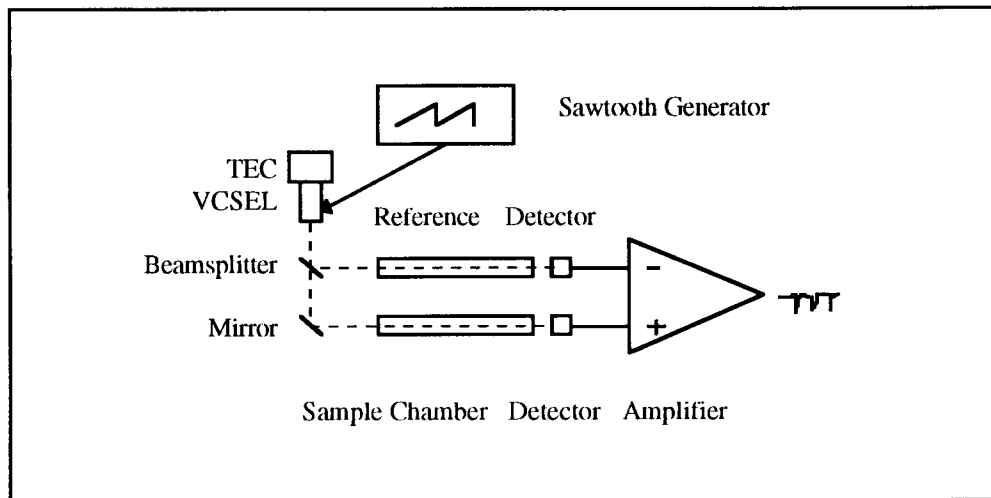



Figure 4. Block diagram of the prototype Wavelength Modulated Spectrometer.

Figure A1, in the appendix, shows oscilloscope waveforms from the prototype spectrometer. The resolution shown compares favorably with FTIR spectrometers and laboratory Laser Spectrometers at a fraction of the cost and size.

The oxygen absorption spectrum in the 761 nm region is shown in Figure 5. Individual oxygen resonance linewidths are on the order of 1 GHz. VCSEL's have 65 MHz linewidth [Olbright, 1992], which is much smaller than an oxygen resonance. The VCSEL's narrow beam is swept from outside to inside and back outside the oxygen resonance by varying the emission wavelength by varying the VCSEL temperature with a sawtooth waveform. The differential amplitude of the sidebands is zero outside the resonance and significant inside the resonance. How significant depends on the slope of the resonance.

# Non Invasive Deep Tissue Oxygenation Sensor

## Phase I Final Report

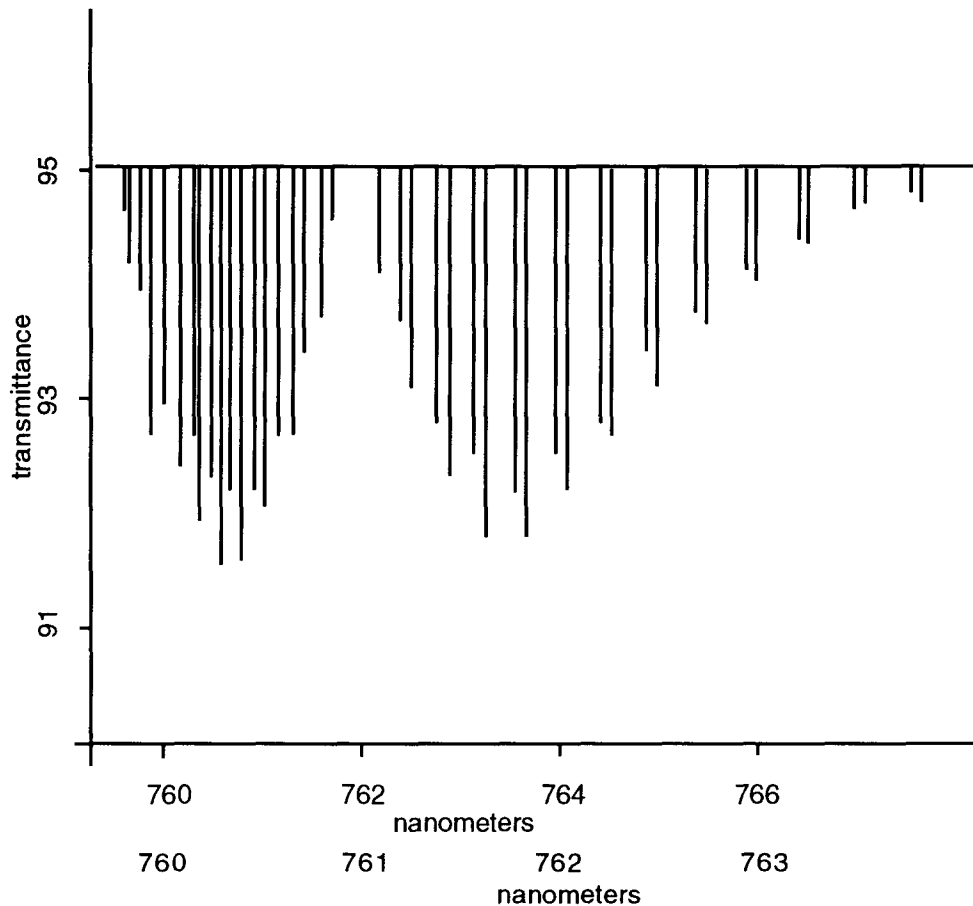


**Figure 5.** Oxygen absorption spectrum. The laser diode has an emission wavelength within this region and is tuned to an absorption line by varying the VCSEL temperature.

The percentage of light absorbed can be very small when sensing gases dissolved in interstitial fluid. This necessitates the use of differential and self-referencing techniques. A particularly successful technique which is well-suited for laser diodes is wavelength modulation in which a laser diode is modulated at over 100 MHz by modulations in its current supply produced by a crystal oscillator. The modulation produces sidebands in the spectral emission as shown in Figure 6. The sidebands are separated from the optical carrier frequency by integral multiples of the modulation frequency. Under weak modulation, the laser spectrum can be approximated by the carrier frequency,  $\omega_0$ , and single upper and lower sidebands at  $\omega_0 - \omega_m$  and  $\omega_0 + \omega_m$  [see Figure 6(a)]. When light with this spectrum is detected in a square-law photodiode, each of the sidebands mixes with the carrier to generate a signal at the modulation frequency. These two signals are equal in amplitude and 180° out of phase, thus they exactly cancel each other out. If one of the sidebands is absorbed prior to detection [see Figure 6(b)], the amplitudes are no longer equal [see Figure 6(c)] and a signal at the modulation frequency  $\omega_m$  appears in the detector current. The amplitude of this signal is directly proportional to the amount of light absorbed. By filtering the detector output at the modulation frequency  $\omega_m$ , the absorption is obtained. Tunability of the laser diode wavelength tunes all three emission wavelengths together. Tuning such that one of the sidebands tunes throughout the absorption line is useful for two main reasons: 1) one insures that the maximum absorption is measured; and 2) the shape of the absorption line can be traced. For the oxygen sensor, this provides an absolute measure of the amount of oxygen in the path of the

# Non Invasive Deep Tissue Oxygenation Sensor

## Phase I Final Report

laser beam. Because the carrier and the sidebands traverse exactly the same optical path, the technique is self-referencing. Since the sideband signals cancel each other out in the absence of differential absorption of one sideband with respect to the other, the technique is inherently differential. Excellent sensitivity of  $10^{-4}$  is therefore obtained.

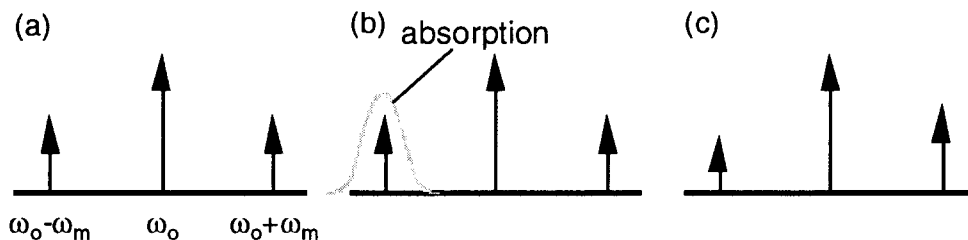

**Figure 6.** Frequencies of light emitted by a modulated laser diode under WM operation. (a) Carrier frequency and equal-amplitude sidebands emitted. (b) Presence of absorption at the lower sideband frequency. (c) Resultant inequality in the sideband amplitudes.

### Assess the overall practicality of the proposed instrument

The goal of the Phase I project was to establish the feasibility of a low-cost, handheld instrument for noninvasive measurement of tissue oxygenation. Spectral analysis of oxyhemoglobin, deoxyhemoglobin, cytochrome aa3, melanin and water showed oxyhemoglobin, deoxyhemoglobin and cytochrome aa3 have useful spectral features. Melanin and water are essentially constants and are not directly involved in cellular metabolism. Oxygen was added to the list of targeted absorbers because it has useful spectral features and are directly involved in tissue oxygenation. Simulation of the proposed instrument indicated 5 VCSEL wavelengths which could resolve levels of oxyhemoglobin, deoxyhemoglobin, cytochrome aa3 and oxygen. The feasibility of using Wavelength Modulated Spectroscopy for the measurement method was established with a prototype spectrometer for oxygen.

### CONCLUSIONS

The proposed instrument for tissue oxygenation is practical. The targeted absorbers have been shown to have absorbances within the range of VCSEL technology. VCSELs are low cost semiconductors with the necessary power, sensitivity and resolution to penetrate a dense scattering media like tissue. The measurement method chosen, WMS, is easily implemented with off-the-shelf components and does not involve any complex, high frequency circuitry.

Three additional instruments based on this technology were proven to be practical during the Phase I project. Bioreactors have a need for low-cost, rapid noninvasive measurement of oxygen and glucose sensors. A noninvasive instrument for oxygen blood gas measurement based on the bioreactor oxygen sensor is practical as is a noninvasive blood sugar meter based on the bioreactor glucose sensor.

# Non Invasive Deep Tissue Oxygenation Sensor


## Phase I Final Report

### REFERENCES

- Chance, B. (1993) Current State of Methodology on Hemoglobin Oximetry in Tissues. ISOTT (Mainz, Germany) Adv. Exp. Biol. and Med. in press
- Chance, B., Kang, K., He, L., Weng, J., Sevick, E. (1993) Highly Sensitive Object Location in Tissue Models with Linear In-Phase and Anti-Phase Multi-element Optical Arrays in One and Two Dimensions. Proc. Natl. Acad. Sci. USA 90:3423-3427.
- Chance, B., Zhuang, Z., UnAh Chu, . Alter, C., L. Lipton (1993) Cognition Activated Low Frequency Modulation of Light Absorption in Human Brain. Proc. Natl. Acad. Sci. USA. 90:3770-3774
- Chance, B., (1992) Noninvasive Approaches to Oxygen Delivery and Cell Bioenergetics in Functioning Muscle. Clin. J. Sport Med. 2:132-138.
- Chance, B. (1991) The Optical Method In Annu. Rev. Biophys. Biophys. Chem. 20:1-18.
- Chance, B., Maris, M., Sorge, J. and Zhang, M.Z. (1990) A Phase Modulation System for Dual Wave length Difference Spectroscopy of Hemoglobin Deoxygenation in Tissues. In Time-Resolved Laser Spectroscopy in Biochemistry II, (J.R. Lakowicz, ed.) Proc. SPIE 1204:481-491.
- Chance, B., Zhang, M.Z., Sorge, J., Leigh, J.S., Kurth, D. and Nioka, S.(1990) A 100 MH Dual Wavelength Phase Modulation Spectrometer. Biophysical J. 57:189.
- Chance, B. (1989) ed. Photon Migration in Tissues, Plenum Press, New York.
- Haselgrove, J., Wang, N.G. and Chance, B. (1992) Investigation of the Non-Linear Aspects of Imaging through a Highly Scattering Medium. Med. Physics. 19:17-23.
- HITRAN 92 Database, Phillips Laboratory, Geophysics Directorate, 1992.
- Shah, A.R., Kurth, C.S., Gwiazdowski, S.G., Chance, B. and Delivoria-Papadopoulos, M. (1992) Fluctuations in cerebral oxygenation and blood volume during endotracheal suctioning in premature infants. J. Pediatrics. 120:769-774.
- Kurth, C.D., Steven, J.M., Nicolson, S.C., Chance, B. and Delivoria-Papadopoulos, M. (1992) Kinetics of Cerebral Deoxygenation during Deep Hypothermic Circulatory Arrest in Neonates. Anesthesiology 77:656-661.
- Mancini, D.M., Bolinger, L., Li, H., Kendrick, K., Chance, B. and Wilson, J.R. (1994) Validation of Near-Infrared Spectroscopy in Man" J. Appl. Physiol in press
- Nioka, S., Chance, B., Smith, D.S., Mayevsky, A., Reilly, M.P., Alter, C. and Asakura, T. (1990) Cerebral Energy Metabolism and Oxygen State during Hypoxia in Neonate and Adult Dog. Ped. Res. 28:54-62.
- Patterson, M.S., Moulton, J.D., Wilson, B.C. and Chance, B. (1990) Applications of Time-Resolved Light Scattering Measurements to Photodynamic Therapy Dosimetry. SPIE (J. Appl. Optics)1203:62-75.
- Patterson, M.S., Chance, B. and Wilson, B.C. (1989) Time Resolved Reflectance and Transmittance for the Noninvasive Measurement of Tissue Optical Properties. J. Appl. Optics 28:2331-36.
- Sevick, E.M., Chance, B., Leigh, J., Nioka, S. and Chance, B. (1991) Quantitation of Time and Frequency Resolved Optical Spectra for the Determination of Tissue Oxygenation. Anal. Biochem. 195:330-351.
- Wilson, B.C., Sevick, E., Patterson, M.S. and Chance, B. (1992) Time-Dependent Optical Spectroscopy and Imaging for Biomedical Applications. Proc. IEEE, 80:918-930.
- Wray S, Cope M, Delphy DT, et al, Characterization of infrared absorption spectra of cytochrome aa3 and hemoglobin for the noninvasive monitoring of cerebral oxygenation, Biochem Biophys Acta, 933:184-192, 1988.

# Non Invasive Deep Tissue Oxygenation Sensor Phase I Final Report

## APPENDIX

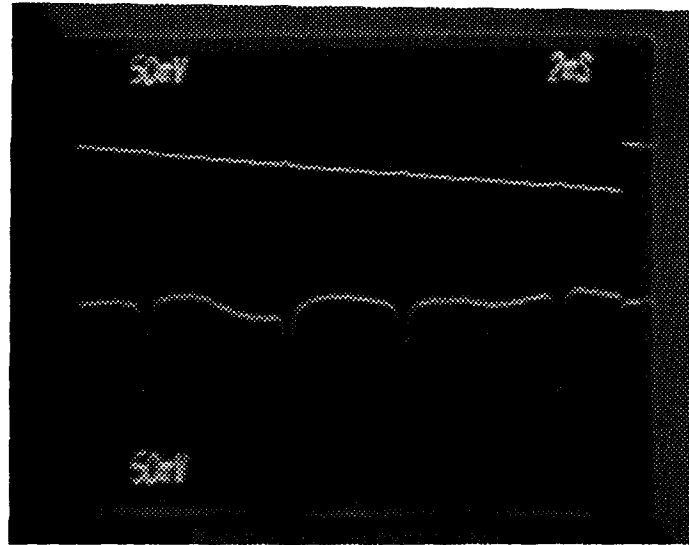

**Figure A1.** Photograph of VCSEL drive current (top) and oxygen spectra (bottom).

Received 2/8/00



DEPARTMENT OF THE ARMY

US ARMY MEDICAL RESEARCH AND MATERIEL COMMAND  
504 SCOTT STREET  
FORT DETRICK, MARYLAND 21702-5012

REPLY TO  
ATTENTION OF:

MCMR-RMI-S (70-1y)

21 Jan 00

MEMORANDUM FOR Administrator, Defense Technical Information  
Center, ATTN: DTIC-OCA, 8725 John J. Kingman  
Road, Fort Belvoir, VA 22060-6218


SUBJECT: Request Change in Distribution Statement

1. The U.S. Army Medical Research and Materiel Command has reexamined the need for the limitation assigned to technical reports written for the attached Awards. Request the limited distribution statements for Accession Document Numbers listed be changed to "Approved for public release; distribution unlimited." These reports should be released to the National Technical Information Service.

2. Point of contact for this request is Ms. Virginia Miller at DSN 343-7327 or by email at virginia.miller@det.amedd.army.mil.

FOR THE COMMANDER:

Encl  
as

  
PHYLIS M. RINEHART  
Deputy Chief of Staff for  
Information Management