

AD 608339

AMRL-TR-64-84


| | | | | |
|------------|---|---|---------|---|
| COPY | 2 | 0 | 3 | 1 |
| HARD COPY | | | | |
| MICROFICHE | | | | |
| | | | \$ 1.00 | |
| | | | \$ 0.50 | |

14p

RESPIRATORY GAS ANALYSIS WITH A TIME-OF-FLIGHT MASS SPECTROMETER

J. B. GRAVES, SENIOR MASTER SERGEANT, USAF
P. R. B. CALDWELL, CAPTAIN, USAF, MC

SEPTEMBER 1964

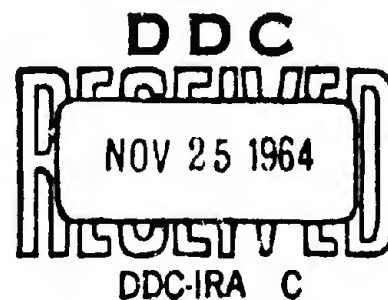

BIOMEDICAL LABORATORY
AEROSPACE MEDICAL RESEARCH LABORATORIES
AEROSPACE MEDICAL DIVISION
AIR FORCE SYSTEMS COMMAND
WRIGHT-PATTERSON AIR FORCE BASE, OHIO

ARCHIVE COPY

NOTICES

When US Government drawings, specifications, or other data are used for any purpose other than a definitely related government procurement operation, the government thereby incurs no responsibility nor any obligation whatsoever; and the fact that the government may have formulated, furnished, or in any way supplied the said drawings, specifications, or other data is not to be regarded by implication or otherwise, as in any manner licensing the holder or any other person or corporation, or conveying any rights or permission to manufacture, use, or sell any patented invention that may in any way be related thereto.

Qualified requesters may obtain copies from the Defense Documentation Center (DDC), Cameron Station, Alexandria, Virginia 22314. Orders will be expedited if placed through the librarian or other person designated to request documents from DDC (formerly ASTIA).

Do not return this copy. Retain or destroy.

Stock quantities available at Office of Technical Services, Department of Commerce, Washington, D. C. 20230.

Change of Address

Organizations receiving reports via the 6570th Aerospace Medical Research Laboratories automatic mailing lists should submit the addressograph plate stamp on the report envelope or refer to the code number when corresponding about change of address.

RESPIRATORY GAS ANALYSIS WITH A TIME-OF-FLIGHT MASS SPECTROMETER

*J. B. GRAVES, SENIOR MASTER SERGEANT, USAF
P. R. B. CALDWELL, CAPTAIN, USAF, MC*

FOREWORD

This study was made in the Altitude Protection Branch of the Aerospace Medical Research Laboratories. The work was done in support of Project No. 7164, "Biomedical Criteria for Aerospace Flight," and Task No. 716404, "Physiological Criteria for Altered Environments." The evaluation reported was studied during a three year period from June 1961 to June 1964.

This technical report has been reviewed and is approved.

WAYNE H. McCANDLESS
Technical Director
Biomedical Laboratory

ABSTRACT

A modified remote sampling system has been devised for the Bendix Model 17-210 time-of-flight mass spectrometer. This has been adapted for use in measuring respiratory gases breath by breath. The stability of the mass spectrometer output shows a variation of $\pm 5\%$. The sensitivity of the system for oxygen is in the range of 0.2% absolute concentration. The accuracy for oxygen is within $\pm 5\%$.

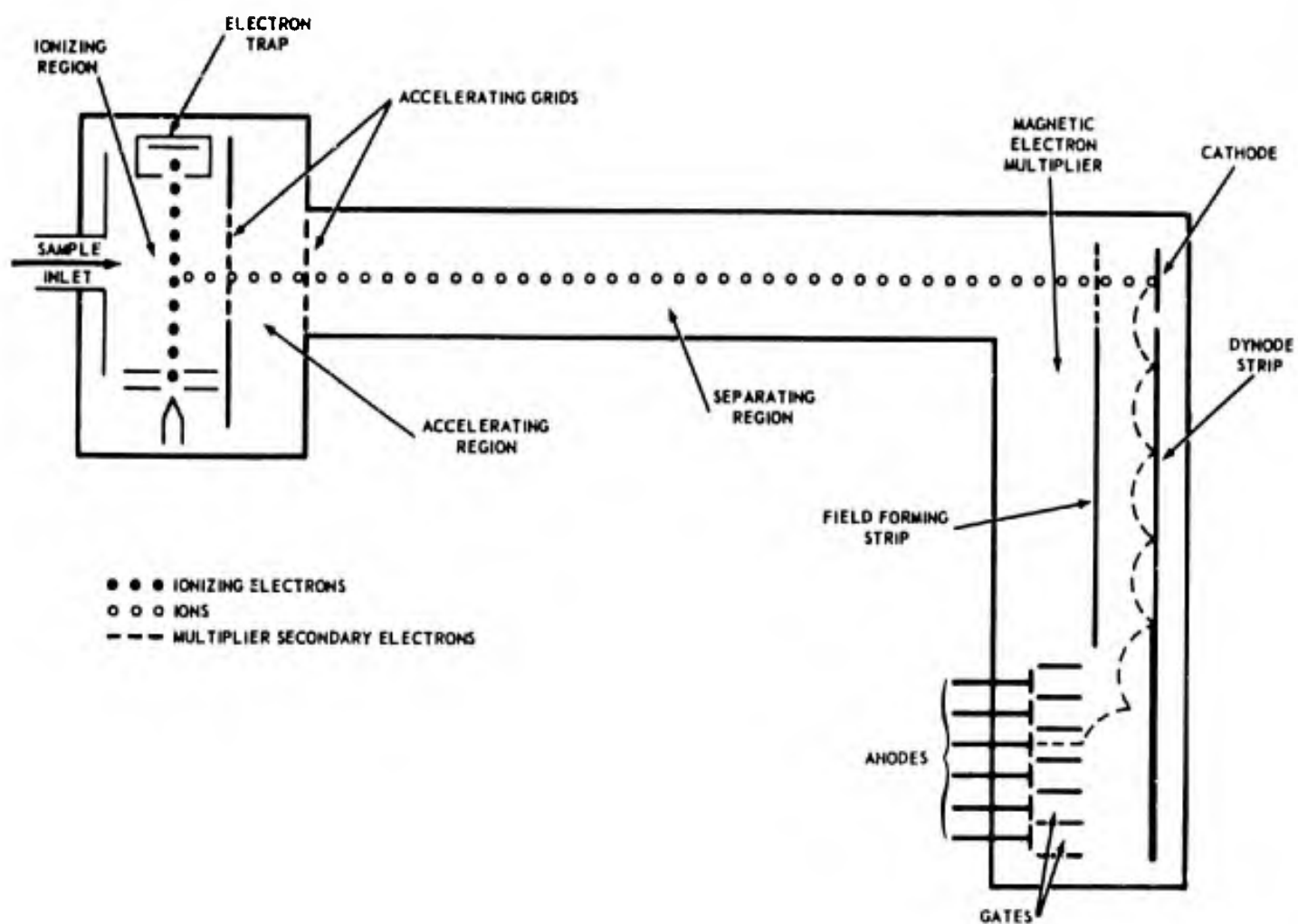


Figure 1. Schematic of Bendix Model 17-210 Time-Of-Flight Mass Spectrometer.

RESPIRATORY GAS ANALYSIS WITH A TIME-OF-FLIGHT MASS SPECTROMETER

SECTION I

INTRODUCTION

This report describes our experience using a Bendix model 17-210 time-of-flight mass spectrometer built by the Research Laboratories Division of the Bendix Corporation. As with any highly complex instrument, the reliability of the results depends on the way it is used. This study, then, was not intended to evaluate the instrument, but was performed solely to apply to measurements of respiratory gases.

SECTION II

DESCRIPTION OF MASS SPECTROMETER

The essential features of this model 17-210 Bendix time-of-flight mass spectrometer are shown in figure 1. The gas to be analyzed is carried from its source by a remote sample system to the sample inlet. The gas enters the low pressure ionizing region through a capillary molecular leak of one to two one thousandth inch diameter. The gas is ionized by a pulsed electron beam emitted from a tungsten filament. The positive ions are ejected from the source region by a pulse applied to the accelerating grids. All ions leave the source region with the same energy and achieve velocities proportional to their masses. Groups of ions separate according to mass during passage through the drift tube. At the end of the drift tube the ion pulses impinge on the cathode emitting a proportionate number of secondary electrons. These are multiplied during their travel down the dynode strip and gated by a pulse to the gate electrode, which is timed with the arrival of the electron multiplier pulse corresponding to the desired mass peak.¹

¹G. J. O'Halloran and W. R. Wyess. The Model 17 Series Bendix Time-Of-Flight Mass Spectrometer. Presented at Eleventh Annual Conference on Mass Spectrometry and Allied Topics. May 1963. San Francisco, California.

SECTION III

MODIFICATION OF SAMPLE SYSTEM

Certain modifications were necessary to adapt this unit to continuous breath by breath measurements of respiratory gases. The sample head was modified as shown in figure 2. This decreased the dead space volume by one half to permit a more rapid response time. The remote sample line consisted of a length up to 20 feet of polyethylene tubing fitted with modified hypodermic needles as shown in figure 3. To construct the sample line, a 25-gage 0.25-inch needle was filed to a blunt end. A stylet was placed in the needle to prevent blocking during its insertion into one end of the desired length of polyethylene tubing. This joint was sealed with household cement, and then adapted to the modified sample head. The connection was lubricated with a silicone gel to make it air tight. The other end of the sample line was similarly fitted with a needle that was then inserted into a mouthpiece. So long as care was taken to avoid plugging the line with silicone gel at the sample head fitting, difficulties with obstruction of the capillary molecular leak were very infrequent.

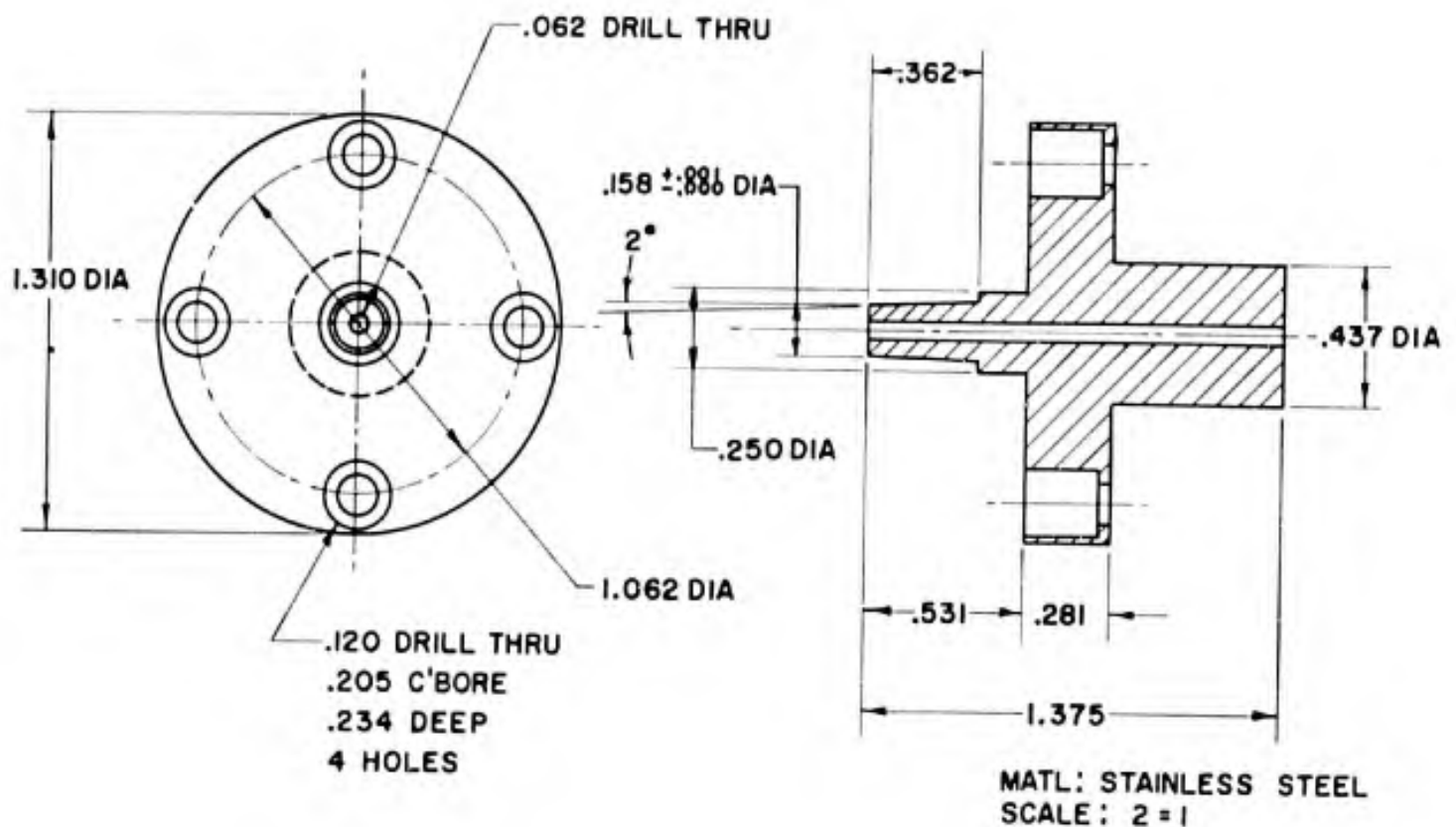


Figure 2. Modified Sample System Head.

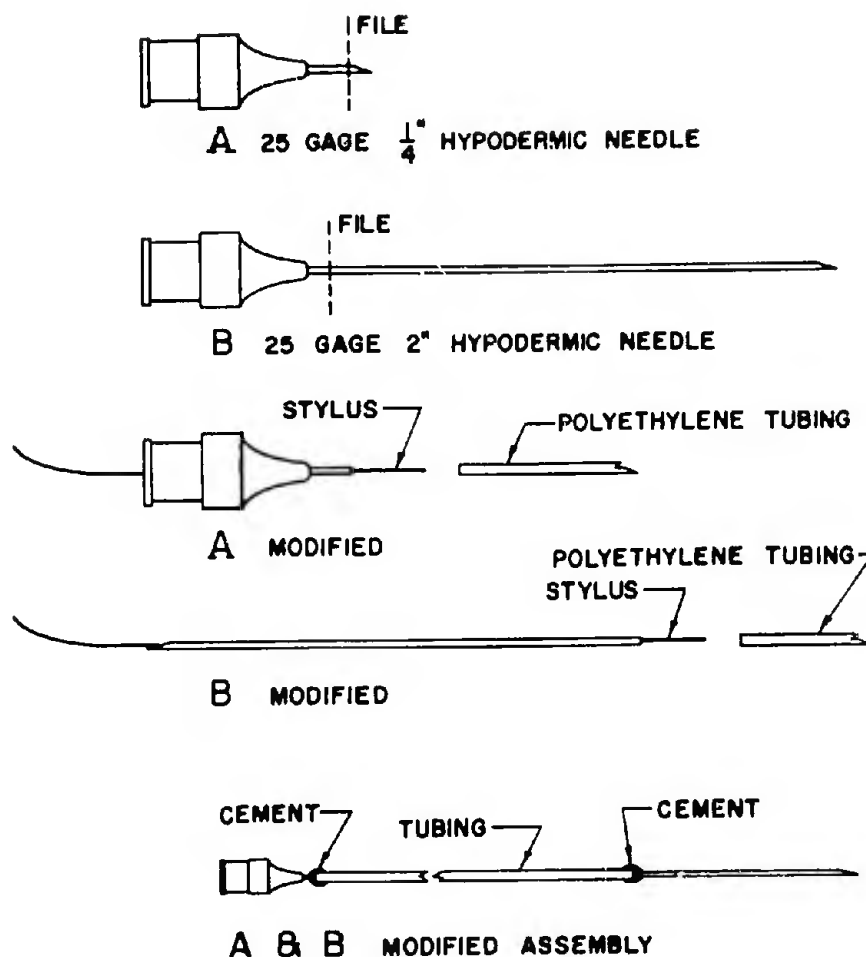


Figure 3. Remote Sample Line Construction.

For most purposes a sample line length of 5-10 feet with a 0.015-inch ID polyethylene tubing maintained at a pressure between 0.8 and 1.0 mm Hg gave the desired result. The output voltage for any given gas in this system is, of course, linearly related to the sample system pressure. See figure 4. This data was obtained by sampling from a closed tube connected to the reduction valve on a tank of pure nitrogen. By altering the reduction valve aperture, the pressure in the sample line could be regulated and maintained while voltage output on mass 28 was recorded.

The recording system used was a D-C amplifier and a model 5-124 Consolidated Electronics Corporation oscillograph. A 90% response time for the entire system of 1 second was obtained with a mass spectrometer time constant setting of 0.5 second. This gave a stable recording and was an appropriate response time for breath-by-breath measurements at a frequency of 10 breaths per minute. The 90% response time could be shortened to 0.3 second by adjustment of the time constant setting. The 90% response time was determined by measuring the time for 90% of a total response after abruptly removing the sample line from a pure oxygen source to room air.

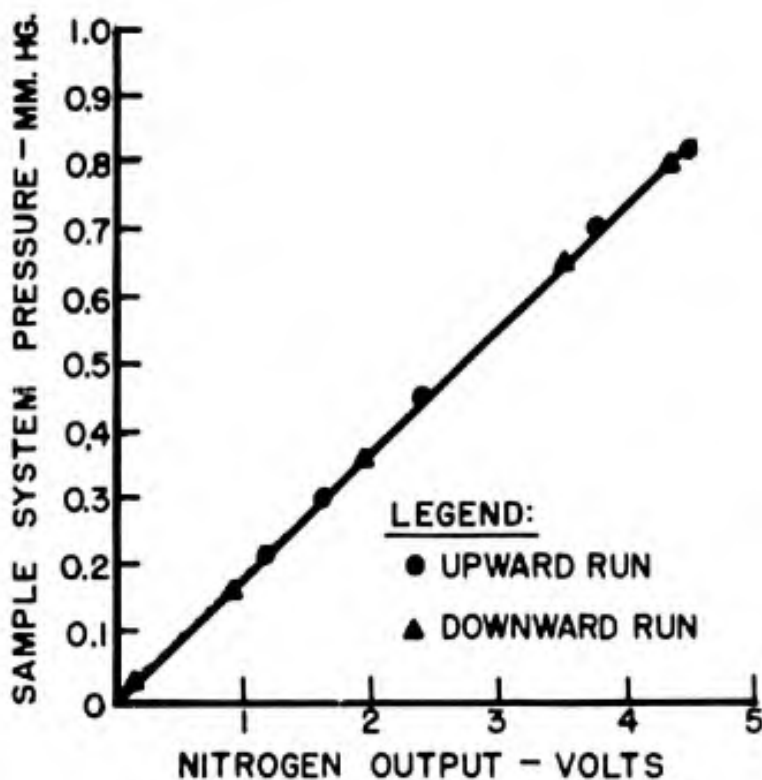


Figure 4.
Relation of Nitrogen Output
to Sample System Pressure.

SECTION IV

OVERALL STABILITY OF SYSTEM

The overall system stability was determined by continuous sampling from a 100-liter Tissot spirometer containing pure nitrogen. The voltage output on mass 28 was measured periodically and found to vary by about $\pm 5\%$. A representative run is shown in figure 5. The gate was checked prior to each determination so that gate drift should not be a factor in this instability. These results illustrate the need for frequent calibration in the use of this machine.

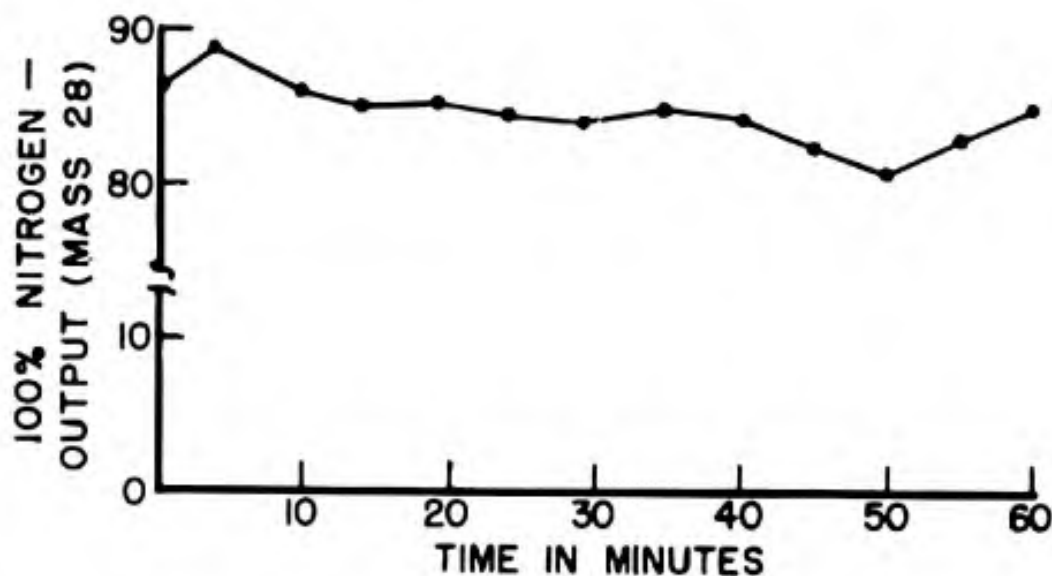


Figure 5. Periodic Determinations of a Single Sample of Pure Nitrogen.

SECTION V

ACCURACY OF OXYGEN MEASUREMENT

The accuracy of gas measurement was tested by comparing the mass spectrometer analysis of various oxygen in nitrogen mixtures to that of the micro-Scholander gas analyzer. The mixtures were prepared in a 100-liter Tissot spirometer and mixed by fan for at least 20 minutes. A sample was collected in a greased syringe, then the mass spectrometer analysis was made at the collecting site, and finally a second syringe sample was obtained. The two syringe samples were analyzed in duplicate on the micro-Scholander. The standard deviation for these analyses was $\pm 0.1\%$.

The mass spectrometer was calibrated with gaseous aviator's breathing oxygen (99.5%) and gaseous water pumped nitrogen (99.5%). Only mass 32 was measured because previous analysis indicated that the splitting of oxygen into masses 32 and 16 was in constant 5:1 ratio over the range of 1-100% concentrations. See figure 6. To avoid any possible error due to gate drift, the spectrum was scanned at a scan rate of 7, time constant setting of 0.02 second and recorded on the oscillograph. The output on mass 32 was then measured on the paper deflection. Background levels of mass 32 were determined at the multiplier voltage and any range settings used during the measurement and subtracted from the measured values of mass 32. The concentrations were then calculated from the measured outputs from the relationship,

$$\frac{99.5}{K_1 - B} = \frac{X}{K_2 - B}$$

where

X is the unknown concentration

B is the background output of mass 32 under operating conditions

K₁ is the output on aviator's breathing oxygen

K₂ is the output of the unknown gas.

The results of this comparison are shown in figure 7 and indicate an accuracy within about $\pm 5\%$ over the range of concentrations measured.

Below concentrations of oxygen of 0.2% the mass spectrometer analysis usually indicated zero concentration. The sensitivity of the mass spectrometer for oxygen, then, appears to be in this range of concentration.

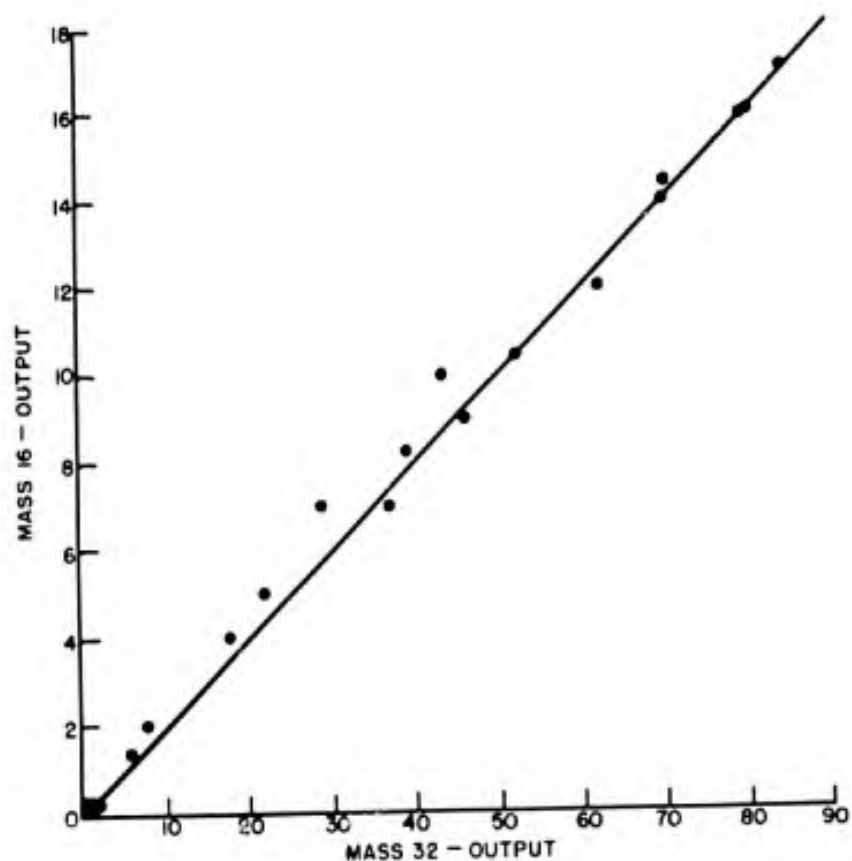


Figure 6. Ratio of Outputs on Masses 32 and 16 for Oxygen Mixtures ranging from 1-100% Concentration.

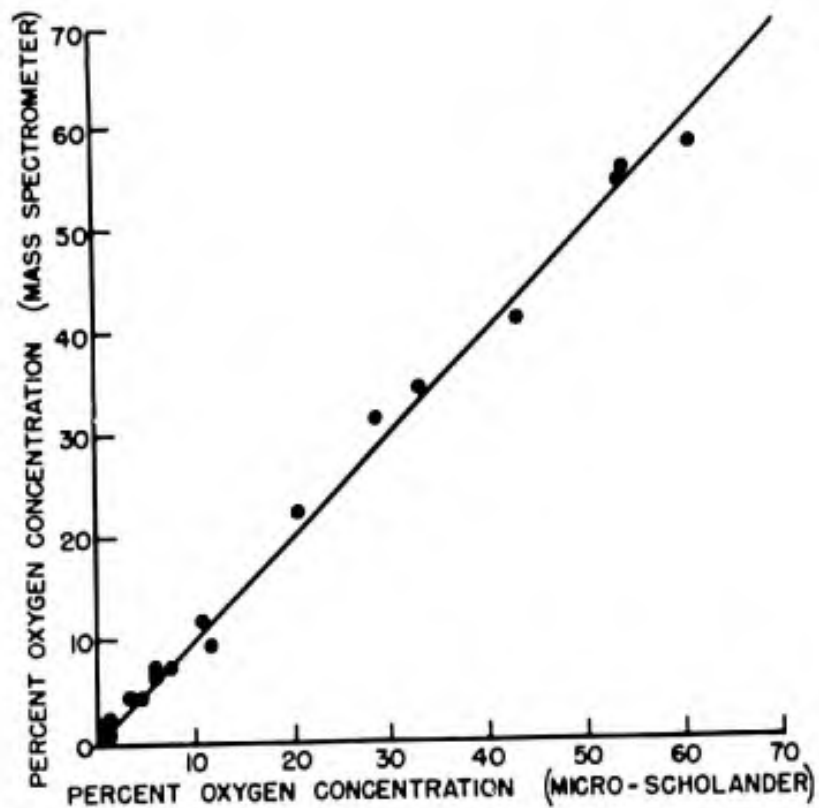


Figure 7. Comparison of Measurements of Oxygen Concentration in Various Mixtures by the Mass Spectrometer and Micro-Scholander Techniques.

| DOCUMENT CONTROL DATA - R&D | | |
|--|--|---|
| <small>(Security classification of title, body of abstract and indexing annotation must be entered when the overall report is classified)</small> | | |
| 1 ORIGINATING ACTIVITY (Corporate author) Aerospace Medical Research Laboratories Wright-Patterson Air Force Base, Ohio | | 2a REPORT SECURITY CLASSIFICATION UNCLASSIFIED |
| | | 2b GROUP N/A |
| 3 REPORT TITLE RESPIRATORY GAS ANALYSIS WITH A TIME-OF-FLIGHT MASS SPECTROMETER | | |
| 4 DESCRIPTIVE NOTES (Type of report and inclusive dates) Final report, June 1961 - June 1964 | | |
| 5 AUTHOR(S) (Last name, first name, initial) Graves, J. B., Senior Master Sergeant, USAF Caldwell, P. R. B., Captain, USAF, MC | | |
| 6 REPORT DATE September 1964 | 7a TOTAL NO. OF PAGES 12 | 7b NO. OF REFS 1 |
| 8a CONTRACT OR GRANT NO. | 9a ORIGINATOR'S REPORT NUMBER(S) AMPL-TR-64-84 | |
| b PROJECT NO 7164 c Task No. 716404 d | 9b OTHER REPORT NO(S) (Any other numbers that may be assigned this report) | |
| 10 AVAILABILITY/LIMITATION NOTICES Qualified requesters may obtain copies of this report from DDC. Available, for sale to the public, from the Office of Technical Services, U.S. Department of Commerce, Washington, D. C. 20230. | | |
| 11 SUPPLEMENTARY NOTES | 12. SPONSORING MILITARY ACTIVITY Aerospace Medical Research Laboratories Wright-Patterson Air Force Base, Ohio | |
| 13 ABSTRACT A modified remote sampling system has been devised for the Bendix Model 17-210 time-of-flight mass spectrometer. This has been adapted for use in measuring respiratory gases breath by breath. The stability of the mass spectrometer output shows a variation of $\pm 5\%$. The sensitivity of the system for oxygen is in the range of 0.2% absolute concentration. The accuracy for oxygen is within $\pm 5\%$. | | |

| 14 KEY WORDS | LINK A | | LINK B | | LINK C | |
|---|--------|----|--------|----|--------|----|
| | ROLE | WT | ROLE | WT | ROLE | WT |
| Spectrum Analyzers, Spectroscopy Spectroscopy, Mass Spectroscopy Respiration, physiology Samplers, instrumentation Medical research | | | | | | |

INSTRUCTIONS

1. ORIGINATING ACTIVITY: Enter the name and address of the contractor, subcontractor, grantee, Department of Defense activity or other organization (*corporate author*) issuing the report.

2a. REPORT SECURITY CLASSIFICATION: Enter the overall security classification of the report. Indicate whether "Restricted Data" is included. Marking is to be in accordance with appropriate security regulations.

2b. GROUP: Automatic downgrading is specified in DoD Directive 5200.10 and Armed Forces Industrial Manual. Enter the group number. Also, when applicable, show that optional markings have been used for Group 3 and Group 4 as authorized.

3. REPORT TITLE: Enter the complete report title in all capital letters. Titles in all cases should be unclassified. If a meaningful title cannot be selected without classification, show title classification in all capitals in parenthesis immediately following the title.

4. DESCRIPTIVE NOTES: If appropriate, enter the type of report, e.g., interim, progress, summary, annual, or final. Give the inclusive dates when a specific reporting period is covered.

5. AUTHOR(S): Enter the name(s) of author(s) as shown on or in the report. Enter last name, first name, middle initial. If military, show rank and branch of service. The name of the principal author is an absolute minimum requirement.

6. REPORT DATE: Enter the date of the report as day, month, year, or month, year. If more than one date appears on the report, use date of publication.

7a. TOTAL NUMBER OF PAGES: The total page count should follow normal pagination procedures, i.e., enter the number of pages containing information.

7b. NUMBER OF REFERENCES: Enter the total number of references cited in the report.

8a. CONTRACT OR GRANT NUMBER: If appropriate, enter the applicable number of the contract or grant under which the report was written.

8b, 8c, & 8d. PROJECT NUMBER: Enter the appropriate military department identification, such as project number, subproject number, system numbers, task number, etc.

9a. ORIGINATOR'S REPORT NUMBER(S): Enter the official report number by which the document will be identified and controlled by the originating activity. This number must be unique to this report.

9b. OTHER REPORT NUMBER(S): If the report has been assigned any other report numbers (*either by the originator or by the sponsor*), also enter this number(s).

10. AVAILABILITY/LIMITATION NOTICES: Enter any limitations on further dissemination of the report, other than those

imposed by security classification, using standard statements such as:

- (1) "Qualified requesters may obtain copies of this report from DDC."
- (2) "Foreign announcement and dissemination of this report by DDC is not authorized."
- (3) "U. S. Government agencies may obtain copies of this report directly from DDC. Other qualified DDC users shall request through _____."
- (4) "J. S. military agencies may obtain copies of this report directly from DDC. Other qualified users shall request through _____."
- (5) "All distribution of this report is controlled. Qualified DDC users shall request through _____."

If the report has been furnished to the Office of Technical Services, Department of Commerce, for sale to the public, indicate this fact and enter the price, if known.

11. SUPPLEMENTARY NOTES: Use for additional explanatory notes.

12. SPONSORING MILITARY ACTIVITY: Enter the name of the departmental project office or laboratory sponsoring (*paying for*) the research and development. Include address.

13. ABSTRACT: Enter an abstract giving a brief and factual summary of the document indicative of the report, even though it may also appear elsewhere in the body of the technical report. If additional space is required, a continuation sheet shall be attached.

It is highly desirable that the abstract of classified reports be unclassified. Each paragraph of the abstract shall end with an indication of the military security classification of the information in the paragraph, represented as (TS), (S), (C), or (U).

There is no limitation on the length of the abstract. However, the suggested length is from 150 to 225 words.

14. KEY WORDS: Key words are technically meaningful terms or short phrases that characterize a report and may be used as index entries for cataloging the report. Key words must be selected so that no security classification is required. Identifiers, such as equipment model designation, trade name, military project code name, geographic location, may be used as key words but will be followed by an indication of technical context. The assignment of links, rules, and weights is optional.

BLANK PAGE