

SECURITY CLASSIFICATION OF THIS PAGE

REPORT DOCUMENTATION PAGE

Form Approved
OMB No. 0704-0188

1a. REPORT SECURITY CLASSIFICATION Unclassified		1b. RESTRICTIVE MARKINGS None	
2a. SECURITY CLASSIFICATION AUTHORITY DTIC ELECTE FEB 07 1990		3. DISTRIBUTION/AVAILABILITY OF REPORT Approved for public release; distribution unlimited.	
2b. DECLASSIFICATION/DOWNGRADING SCHEDULE AD-A217 683		5. MONITORING ORGANIZATION REPORT NUMBER(S) AFOSR-TR-90-0115	
6a. NAME OF PERFORMING ORGANIZATION Michigan State University College of Human Medicine		6b. OFFICE SYMBOL (If applicable)	7a. NAME OF MONITORING ORGANIZATION Air Force Office of Scientific Research/NL
6c. ADDRESS (City, State, and ZIP Code) Dept. of Pediatrics/Human Development East Lansing, Michigan 48824		7b. ADDRESS (City, State, and ZIP Code) Building 410 Bolling AFB, D.C. 20332-6448	
8a. NAME OF FUNDING/SPONSORING ORGANIZATION AFOSR	8b. OFFICE SYMBOL (If applicable) NL	9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER AFOSR-89-0114	
8c. ADDRESS (City, State, and ZIP Code) Building 410 Bolling AFB, D.C. 20332-6448		10. SOURCE OF FUNDING NUMBERS	
		PROGRAM ELEMENT NO. 61104D	PROJECT NO. 3842
		TASK NO. A4	WORK UNIT ACCESSION NO.
11. TITLE (Include Security Classification) DURIP - Upgrade of the Meridian ACAS-470 for Toxicological Research			
12. PERSONAL AUTHOR(S) James E. Trosko, Ph.D.			
13a. TYPE OF REPORT Final	13b. TIME COVERED FROM 12/01/88 TO 1/30/89	14. DATE OF REPORT (Year, Month, Day) 01/18/90	15. PAGE COUNT 7
16. SUPPLEMENTARY NOTATION Meridian Upgrade; Application notes & Preprints; ACAS Image data and Abstract.			
17. COSATI CODES		18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number)	
FIELD	GROUP	SUB-GROUP	
		Meridian ACAS-570; Laser image analyzer; Quantitative fluorescence detection; Mechanisms of chemical toxicities	
19. ABSTRACT (Continue on reverse if necessary and identify by block number) This proposal was to have our current Meridian ACAS 570 laser cytometer upgraded with a 5 watt coherent laser, Olympus microscope with a CCD camera and video monitor, a mass data storage computer and dual detector system. This upgrade allowed us new dimensions of research potentials for studying mechanisms of chemical-induced toxicities in living cells, specifically enabling us to detect if chemicals caused toxicities by altering free Ca ⁺⁺ level, causing pH changes or induction of free radicals in single cells. During this year, we have pioneered the application of this new instrument in the demonstration that (a) heptachlor, a toxic pesticide/neurotoxicant, caused increased of free Ca ⁺⁺ ; (b) several cytotoxic chemicals could generate free radicals; and (c) oncogene products could be detected using fluorescent antibodies; and (d) UV-induced DNA damage could be quantified using fluorescent antibodies. Original contains color plates. All DTIC reproductions will be in black and white.			
20. DISTRIBUTION/AVAILABILITY OF ABSTRACT ☒ UNCLASSIFIED/UNLIMITED ☐ SAME AS RPT ☐ DTIC USERS		21. ABSTRACT SECURITY CLASSIFICATION	
22a. NAME OF RESPONSIBLE INDIVIDUAL T. Jan Carvery, Lt. Col., USAF		22b. TELEPHONE (Include Area Code) 202-747-5021	22c. OFFICE SYMBOL NL

DD Form 1473, JUN 86

Previous editions are obsolete.

SECURITY CLASSIFICATION OF THIS PAGE

20030206008

88 JAN 1990

TERMINAL REPORT APO-RK. 90-0115

DURIP-Upgrade of the Meridian
ACAS-470 for Toxicological Research

12/01/88 - 11/30/89

James E. Trosko, Ph.D.
Dept. of Pediatrics/Human Development
Michigan State University
East Lansing, Michigan 48824

REPRODUCTION QUALITY NOTICE

This document is the best quality available. The copy furnished to DTIC contained pages that may have the following quality problems:

- **Pages smaller or larger than normal.**
- **Pages with background color or light colored printing.**
- **Pages with small type or poor printing; and or**
- **Pages with continuous tone material or color photographs.**

Due to various output media available these conditions may or may not cause poor legibility in the microfiche or hardcopy output you receive.

☐ **If this block is checked, the copy furnished to DTIC contained pages with color printing, that when reproduced in Black and White, may change detail of the original copy.**

I. Terminal Report of Progress

The purpose of this proposal was to upgrade our Meridian ACAS 470 laser cytometer, a novel and powerful interactive laser/computer assisted image analyzer. The device was designed to detect, using a microlaser beam, tracer fluorescent molecules in a cell and generate a false-color digitized image which could be quantitatively analyzed [see Appendix A]. This grant enabled us to bring this instrument up to the "state of the art", by allowing us to have a laser which could generate a powerful UV laser beam, as well as the visible laser beam we had. This gave us the potential to study molecules and mechanisms which could only be detected with UV-excitable dyes, as well as store and analyze more complex data, as well as manipulate the experiments with more precision.

The three original aims, i.e., (a) to detect intracellular free Ca^{++} [an indicator of both adaptive and maladaptive cellular physiology]; (b) to detect changes in intracellular pH with a pH sensitive dye (also a measure of both adaptive/maladaptive changes in cellular physiology caused by chemicals; and (c) to detect intracellular free radical production in single cells with a fluorescent dye caused by many toxic chemicals, were achieved

II. Results of Upgrade

The Meridian company modified our instrument within months after the grant was awarded. Because of this upgrade, we were able, in collaboration with the Meridian company, to develop several new applications with this new instrument leading to several manuscripts which are now in press or in preparation, as well as the development of new "Application Notes" [Appendix B].

III. Benefits to Research, Department and University

To re-iterate, since my USAFOSR grant [AFOSR-89-0325] is to study the mechanisms of chemical toxicants which act by inhibiting intercellular communication, the upgrade of this unique interactive laser image analyzer has made it possible to (a) do quantitative studies, on a single cell basis, which were either impossible or extremely difficult to do and inaccurate prior to its development; (b) increase our productivity by allowing us to generate data and to analyze data, simultaneously, without tying up the machine for one or the other; and (c) develop new potential applications for toxicological research [i.e., our recent use to detect specific oncogene products in normal or cancer cells].

As a result of the upgrade, both the number and diversity of studies using the Meridian ACAS-570 have increased. For example, the instrument is signed up for use almost every day, including many evenings and weekends. Free calcium studies are being done by scientists (Drs. Kalimi and Wilson) to study the effects of



Dist	Avail and/ Special
10-1	

hormone imbalance on calcium homeostasis in testicular tissue; the effect of anti-parasitic drugs is being studied in schistosomes by Dr. Bennett (Pharmacology and Toxicology Dept., MSU); oncogene detection with fluorescent antibodies in cancer cells by E. de Feijter (a foreign scholar); stem-cell lineages in human breast tissue (Dr. S. Nakatsuka); DNA damage with fluorescent antibodies to UV-induced pyrimidine dimers by Dr. C.C. Chang (Dept. Pediatrics/Human Development); the reaction rates of enzymatic functions in bio-engineered matrices (Dr. M. Worden, Chemical Engineering Dept., MSU); a new way to measure heterologous intercellular communication (Dr. G. Kalimi, Dept. Pediatrics/Human Development) [see Appendix C].

In addition, benefit to the University accrues in many ways: (a) Collaboration with other research groups studying toxicology at different levels; for example, Dr. K. Honn of the Chemistry Department at Wayne State University has used our upgrade to initiate studies on the detection of free radicals after cells were exposed to eiconosoids; Dr. James Klaunig of the Medical College of Ohio in Toledo is collaborating with us to study chemical inhibition of intercellular communication in liver slices; and Dr. Phil Liu of Case/Western-Reserve University (Cleveland) is collaborating with us to detect DNA polymerases in single cells; (b) Because we were the first group using the Meridian ACAS and its upgrade, we are now communicating with others across the country and sharing new information at applications; (c) Because of the major technological breakthroughs represented by this instrument, I also have been asked to give workshops and seminars [i.e., At the latest Cell Biology meeting in Houston, Texas, I conducted a workshop].

Finally, this upgrade has also helped the University directly by attracting other grants utilizing this instrumentation [NSF grant to Dr. Worden in Chemical Engineering; NIEHS grant to Dr. L. Fischer, Center for Environmental Toxicology, MSU]. Several papers, abstracts and multiple speaking invitations, including one to Germany [International Gap Junction Meeting]; to England [International Society of Toxicology Meeting]; and recently to France [Workshop on Alternatives to Animal Toxicity Testing].

IV. Papers, Meetings and Seminars

1. J.E. Trosko and C.C. Chang, "Stem cell theory of carcinogenesis", Toxicology Letters, in press.
2. J.E. Trosko, "Towards understanding carcinogenic hazards: A crisis in paradigms", J. Amer. College of Toxicol., 8:1121-1132, 1989.
3. J.E. Trosko, C.C. Chang and B.V. Madhukar, "Cell-cell communication: Relationship of stem cells to the carcinogenic process", In: Mouse Liver Carcinogenesis;

Mechanisms and Species Comparisons, T.J. Slaga, Ed., Alan R. Liss, Inc., New York, in press.

4. M. El-Fouly, J.E. Trosko, C.C. Chang and S.T. Warren, "Potential role of the human Ha-ras oncogene in the inhibition of gap junctional intercellular communication", Molecular Carcinogenesis, 2:131-135, 1989.
5. J.E. Trosko, C.C. Chang and B.V. Madhukar, "In vitro analysis of modulators of intercellular communication: Implications for biologically based risk assessment models for chemical exposure", Toxicology In Vitro, in press.
6. S.Y. Oh, B.V. Madhukar, C.C. Chang, J.E. Trosko and E. Beyer, "Characterization of gap junctional communication deficient mutants from a HG-PRT deficient rat liver epithelial cell line", provisionally accepted to
7. A.W. de Feijter, J.S. Ray, C.M. Weghorst, J.E. Klaunig, J.I. Goodman, C.C. Chang and J.E. Trosko, "Correlation between the methylation status of the v-Ha-ras oncogene, gap junctional communication and tumorigenicity in rat liver epithelial cells", Molecular Carcinogenesis, accepted for publication.
8. B.V. Madhukar, J.E. Trosko, and C.C. Chang, "Chemical, oncogene and growth factor modulation of gap junctional communication in carcinogenesis". In: Cell Interactions and Gap Junctions. Vol. 1, N. Sperelakis and W.C. Cole, eds., CRC Press, Boca Raton, FL, pp. 143-157, 1989.

Abstracts

1. B.V. Madhukar, S.Y. Oh, E. de Feijter and J.E. Trosko, "Inhibition of intercellular communication by toxic xenobiotic chemicals in vitro in a human epithelial cell system", Soc. of Toxicology, Atlanta, GA, Feb. 27 - March 3, 1989. The Toxicologist 9:4, 1989.
2. S.G. Lilly, E. de Feijter, C.M. Weghorst, C.C. Chang, B.V. Madhukar, M. El-Fouly, J.E. Trosko, and J.E. Klaunig, "Reduced gap junctional intercellular communication and tumorigenicity of rat liver epithelial F344-WB cells". Soc. of Toxicol., Atlanta, GA, Feb. 27-March 3, 1989. The Toxicologist 9:123, 1989.
3. J.E. Trosko, "Mechanisms of chemical and oncogene modulation of gap junction and communication during carcinogenesis", Workshop on the Molecular and Cell Biology of Gap Junctions, Irsee, Germany, July 18-23,

1989.

4. J.E. Trosko, C.C. Chang, and B.V. Madhukar, "In vitro analysis of modulators of intercellular communication: Implication to mechanisms of tumor promotion and to predictions of potential tumor promoters". 2nd Intern. Conf. on Practical In Vitro Toxicology, Nottingham, U.K., July 23-27, 1989.
5. J.E. Trosko, J.E. Klaunig, B.V. Madhukar, C.C. Chang, E. de Feijter, and G. Kalimi, "Chemical and oncogene modulation of intercellular communication during carcinogenesis". Symposium on Molecular Cell Biology of Liver Growth and Function, Lake Placid, NY, Aug. 13-16, 1989.
6. J.E. Trosko, J.E. Klaunig, M. Yeager, A. Koestner, M. El-Fouly, J. Buboltz, B. Cool, and C.C. Chang, "Correlation between the reduction of gap junctional communication and tumorigenesis in rat glial and liver epithelial cells containing the expressed neu oncogene". Society of Toxicology Meeting, Miami Beach, FL, Feb. 12-16, 1990.
7. A.W. de Feijter, M.W. Lieberman, and J.E. Trosko, "Up-regulation of a metallothionein-ras T24 fusion gene in rat liver epithelial cells is correlated with the down regulation of gap junction function". Society of Toxicology Meeting, Miami Beach, FL, Feb. 12-16, 1990.
8. R.J. Ruch, J.E. Trosko, B.V. Madhukar, P. Somani, and J.E. Klaunig, "Restoration of gap junction intercellular communication by lovastatin in WB-H-rasII in rat liver epithelial cells". Society of Toxicology Meeting, Miami Beach, FL, Feb. 12-16, 1990.
9. B.V. Madhukar, H. Hong Xu, B. Lockwood, and J.E. Trosko, "Inhibition of intercellular communication by environmental chemicals in rat Leydig cells in vitro". Society of Toxicology Meeting, Miami Beach, FL, Feb. 12-16, 1990.
10. W.J. Paradee, B.V. Madhukar, and J.E. Trosko, "Mezerein inhibition of intercellular communication and activation of protein kinase C in human kidney epithelial cells". Society of Toxicology Meeting, Miami Beach, FL, Feb. 12-16, 1990.
11. G. Kalimi, J.E. Trosko, L.L. Hampton, S.S. Thorgerirsson, and A.C. Huggett, "Gap junctional intercellular communication in raf and raf/myc transformed rat liver epithelial cell lines". American Association of Cancer

Research Meeting, Washington, D.C., May 23-26, 1990.

Meetings Attended and Seminars Given

1. Invited seminar speaker, "Chemical modification of intercellular communication: Implications for nongenotoxic toxicity". Program in Toxicology seminar series, Rutgers University, New Jersey, March 8, 1989.
2. Invited seminar speaker, "Chemical promoters, oncogenes, and growth factor modulation of intercellular communication: A unified theory". Loma Linda, California, March 21, 1989.
3. Invited symposium speaker, "Inhibition of intercellular communication during radiation and chemical carcinogenesis. Radiation Research meeting, Seattle, Washington, March 27, 1989.
4. Invited workshop speaker, "Promotion as a factor in risk assessment". Carcinogen Risk Assessment, Society for Risk Analysis, Wash. D.C., April 4, 1989.
5. Invited symposium speaker, "Oncogenesis and abnormal intercellular communication and its implication to the cause and prevention of carcinogenesis", and "Clinical implications of intercellular communication on wound healing". Quebec Association of General Surgeons, Quebec City, Canada, May 4-5, 1989.
6. Invited seminar speaker, "Oncogenes, tumor promoters and growth factor inhibition of intercellular communication: Integrative theory of cancer". University of Hamburg, Germany, July 17, 1989.
7. Invited symposium speaker, "Mechanisms of chemical and oncogene modulation of gap junctional communication during carcinogenesis". Molecular Cell Biology of Gap Junctions Symposium, Irsee, Germany, July 18-23, 1989.
8. Invited symposium speaker, "In vitro analysis of modulators of intercellular communication: Implications to mechanisms of tumor promotion and to predictions of potential tumor promoters". Practical In Vitro Toxicology Symposium, University of Nottingham, England, July 23-27, 1989.
9. Invited symposium speaker, "Chemical and oncogene modulation of intercellular communication during carcinogenesis". Symposium on Molecular, Cellular Biology of Liver Growth and Function, Lake Placid, NY, Aug. 13-16, 1989.

10. Invited seminar speaker, "Oncogenes, tumor promoters, and growth factors: An integrated theory of carcinogenesis". Eppley Cancer Center, Univ. of Nebraska, Omaha, Sept. 28, 1989.
11. Invited seminar speaker, "Role of chemical inhibition of gap junctional communication in pharmacological and toxicological effects of drugs". Bristol Myers, Syracuse, NY, Oct. 9, 1989.
12. Invited symposium speaker, "Modulated intercellular communication: Consequence of extracellular communication". Eiconosoids and Bioactive Lipids in Cancer and Radiation Injury, Detroit, MI, Oct. 11-14, 1989.
13. Invited seminar speaker, "Modulation of gap junctional communication: Pharmacological and toxicological implications". Park-Davis Warner Lambert Pharmaceuticals, Ann Arbor, MI, Oct. 17, 1989.
14. Invited seminar speaker, "A crisis in paradigms: The roles of genotoxic versus epigenetic mechanisms in toxicology". Environmental and Industrial Health Department, Univ. of Michigan, Ann Arbor, MI, Dec. 1, 1989.

APPENDIX A

Information on New Meridian ACAS 570 Upgrade

APPENDIX B

Application Notes on new techniques that the Upgrade of ACAS can perform:

1. Calcium measurements.
2. Free Radicals.
3. Reprint "In situ (6-4) photoproduct determination by laser cytometry and autoradiography".
4. Reprint "Characterization of an in vitro human kidney epithelial system to study gap junctional intercellular communication".

APPENDIX C

New technique to measure heterologous and homologous cell-cell communication using the ACAS:

1. Actual ACAS-Image data.
2. Abstract.

ACAS 570

The Interactive Laser Cytometer

*Has Your Number
For...*

- *Fluorescence Intensity/Cell*
- *Concentration*
- *Mutant Selection*
- *Transport*
- *Diffusion*
- *Ionic Flux*
- *Migration*

And More...

MERIDIAN™

ACAS 570...

The Instrument for Quantitative Adherent Cell Analysis and Sorting

The ACAS 570 Interactive Laser Cytometer is uniquely designed for researchers who use quantitative fluorescence techniques to study adherent cells.

Conceived and developed by cell biologists for cell biologists, the ACAS — **Adherent Cell Analysis and Sorting** — provides fluorescence measurements of structure, function, and response at the single cell level. Applications range from cell and molecular biology to toxicology, immunology, and neurobiology.

Designed to insure you get the **numbers** you need, the ACAS 570 is the only instrument that provides capabilities to perform all of the following:

- Measure fluorescence distribution within cells
- Quantitate cellular responses using ratio-metric methods for Ca^{2+} , pH and other ions
- Measure macromolecular diffusion, permeation, and interaction
- Isolate and clone cells without removing them from culture
- Sort adherent cells based on physical and biochemical measurements
- Perform quantitative confocal imaging

Investigate these and other innovative applications that make the ACAS 570 — **The Cytometer for All Reasons.**

REASON — INTERACTIVE LASER CYTOMETRY


The ACAS 570 makes cytometry **interactive**. Only the ACAS lets you

- Isolate individual cells without removing them from their growth surface
- Create cellular fluorescence gradients by photobleaching or photoactivation
- Measure macromolecular diffusion or flux
- Monitor intracellular communications
- Perform microsurgery on cells


With the ACAS, you can use the laser for more than illumination. Now you can use the laser as a "photon gun," to photobleach or photoactivate for real-time measurements of cellular dynamics. Turn up the power and eliminate unwanted contaminating cells, or isolate and retain rare cells.




Two-Dimensional, pseudo-color fluorescent image of WISH cell stained with the membrane potential dye $\text{DiOC}_2(3)$. (Courtesy of Dr. Bruce Jensen, Zynaxis Cell Science, Inc.)

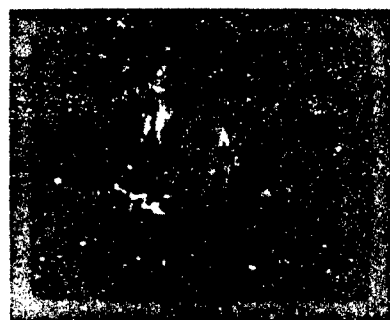

Intensity dimension plot of the membrane potential dye $\text{DiOC}_2(3)$ in WISH Cell. (Courtesy of Dr. Bruce Jensen, Zynaxis Cell Science, Inc.)

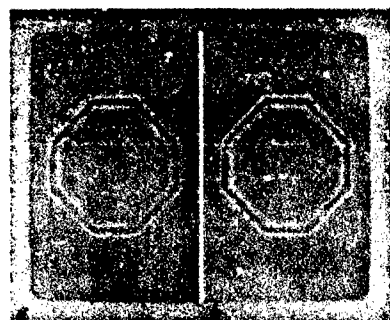

a. CHO cell isolated by the "Cookie Cutter"™ technique. (Courtesy of Dr. Margaret Wade, Meridian Instruments, Inc.)
b. Cell out-growth after two days

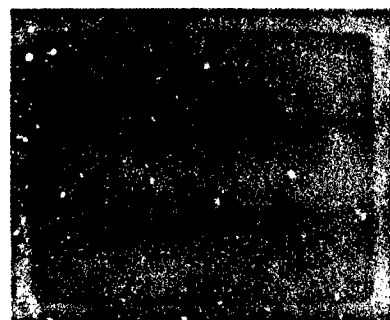

Fluorescence photobleaching of carboxyfluorescein diacetate labeled human teratocarcinoma cells as a measurement of cell-cell communication. Polygon enclosed cells are bleached and monitored for fluorescence recovery. (Courtesy of Dr. Margaret Wade, Meridian Instruments, Inc.)

Sorting Cells the ACAS Way

The ACAS is the only instrument that lets you select anchorage-dependent cells directly from heterogeneous cultures for isolation, propagation, and analysis.

The unique "Cookie-Cutter"™ technique automatically isolates single cells growing in film-lined dish based on fluorescence measurements or morphology.

A high intensity laser beam encircles the selected cell, fusing octagon-shaped "cookies" to the culture dish. When the excess film is pulled from the dish, the selected cell remains attached viable for continued propagation.

Use this powerful **interactive** technique for sorting:

- Mutant cells
- Transfected cells
- Hybridomas
- Rare cells for cloning

You can also use the ACAS laser as a "photon scalpel" to interact with cells grown on specially treated, heat-absorptive glass cover slips or culture dishes.

Command the laser to produce localized heat of target areas to:

- Select homogeneous cell populations
- Ablate single cells
- Eliminate large areas of cells


FRAP — Dynamic Parameters of Living Cells.

Fluorescence Redistribution After Photobleaching (FRAP) techniques involve photobleaching fluorescent dye molecules in discrete sub-cellular areas with a high intensity laser pulse creating non-fluorescent patches. The rate at which the patches are filled by unbleached molecules from surrounding areas is measured by periodic, low intensity laser scanning.

Use this unique, non-destructive tool to measure a wide range of cellular dynamics, such as:

- Receptor mobility
- Membrane transport
- Macromolecular assembly
- Cell-cell communication

™ "Cookie Cutter" is a trademark of Meridian Instruments, Inc.

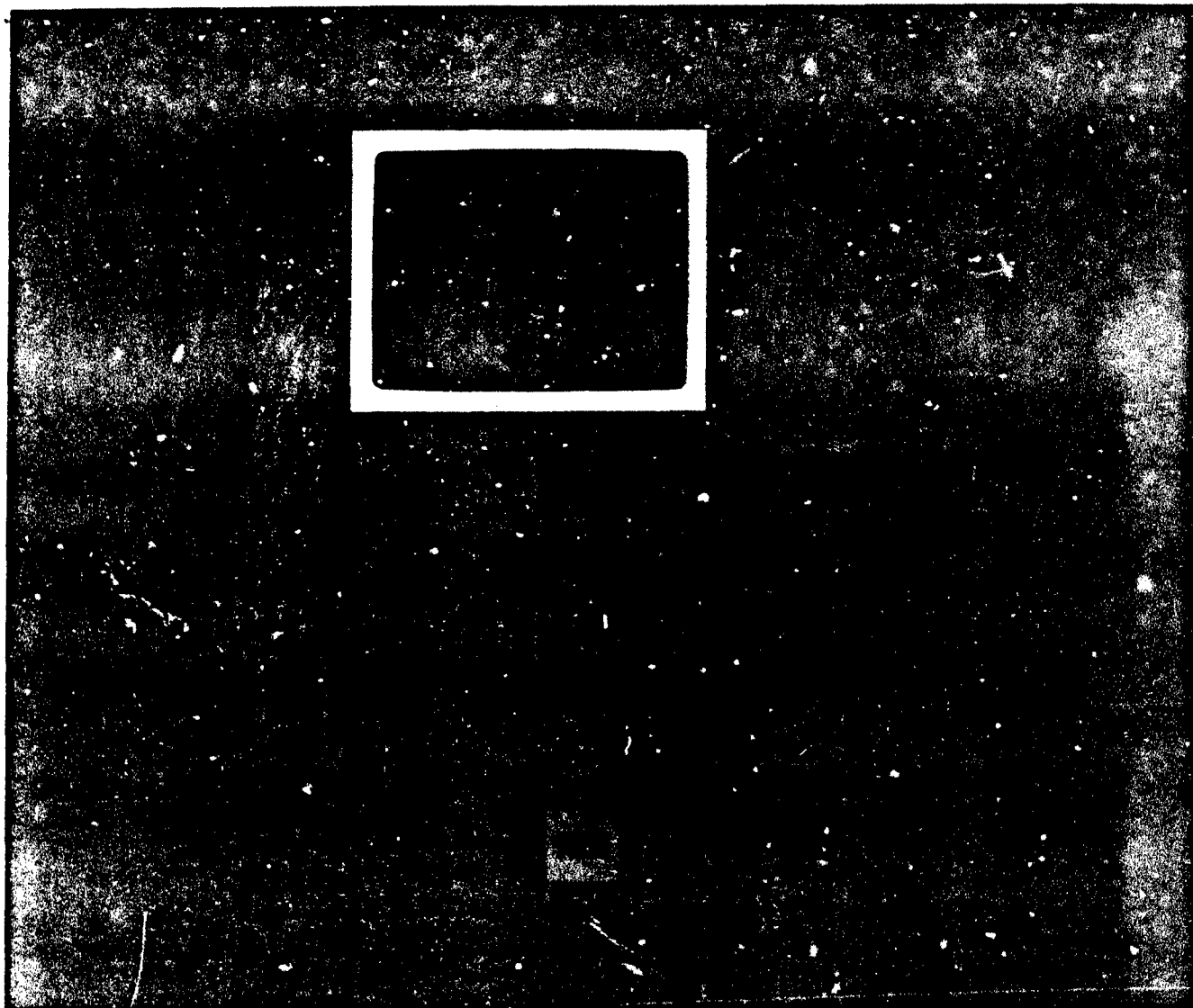

From Signal to Image—The Basics

The schematic above shows you the basic optical and signal paths. The beam is directed along the center of the optical path and focused to a small spot by the microscope objective onto a fluorescently labeled cell.

Fluorescence emissions are collected by the objective, and directed to photomultiplier tubes by a series of mirrors and lenses. Appropriate wavelengths are selected by the choice of filters.

The following sequence describes how the ACAS acquires and displays fluorescence data:

1. Cells of interest are positioned above the objective by the X-Y scanning stage. A touch of the finger starts the stage moving to begin the scan.
2. As the stage moves, the focused laser beam is briefly pulsed by the XOM. Data is collected at intervals as small as 0.25 micron.
3. The photometric detection system captures the fluorescence emissions during the laser pulse. The resulting electronic signals are digitized and stored by the computer along with the stage coordinates.
4. The computer displays the signal resulting from a single laser pulse as a pixel with a color or gray value that represents the intensity of the emission at that point. Single or two dimensional scans are obtained by repeating this process while the state continues to move until the selected number of data points are acquired. (For example, see WISH cell on facing page.)
5. Fluorescence intensities are converted to numbers (i.e., concentration, rate, of diffusion) depending on the software package selected.

REASON • QUANTITATIVE FLUORESCENCE MEASUREMENTS

Specifically designed for low level fluorescence detection and quantitation in living cells, the ACAS 570 gives you numbers you can count on.

Only with the ACAS can you:

- Measure physical and biochemical changes within cells
- Automatically acquire data and calculate histograms of fluorescence variations within cell populations
- Directly obtain macromolecular rates of diffusion
- Measure Ca^{2+} , pH and other ions using ratiometric methods
- Obtain rapid (microseconds) temporal measurements

Sensitive, Quantitative Measurements with Minimal Photobleaching


The ACAS 570 has harnessed the unique properties of coherent, monochromatic laser illumination to provide researchers with a powerful tool for fluorescence quantitation.

Fluorescent probes provide a high degree of sensitivity and specificity, but are susceptible to photobleaching, especially at the low concentrations required for the study of living cells.


The ACAS is designed to resolve this problem through controlled laser intensity, duration, and spot size. The laser beam illuminates **only** the area to be measured, and **only** during the time of measurement.

A computer controlled acousto-optic modulator (AOM) provides rapid changes in the laser spot intensity and duration. Synchronized with the XY scanning stage, the AOM sequentially exposes small areas of a cell to brief (microseconds) laser pulses.


During a scan, therefore, each point in the field is illuminated for only a fraction of the time as required by full-field illumination systems, and with better sensitivity. The results: **photobleaching under control.**

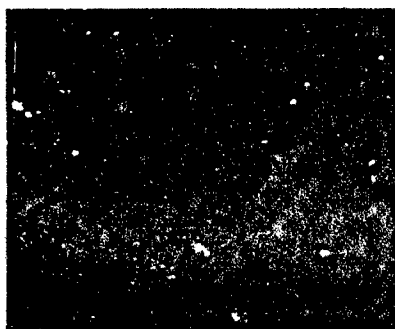


Serum-stimulated increase in intracellular calcium levels in human teratocarcinoma cells labeled with Indo-1. (Courtesy of Dr. Margaret Wade, Meridian Instruments, Inc.)

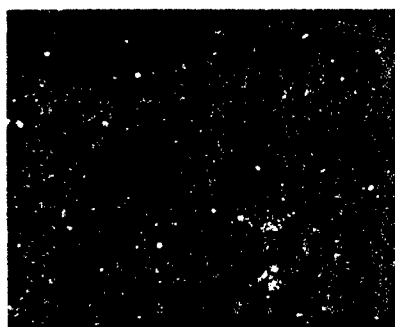

Response of 3T3 cells labeled with the fluorescent pH indicator SNARF-1 to nigericin. (Courtesy of Asimina Jura, Meridian Instruments, Inc.)

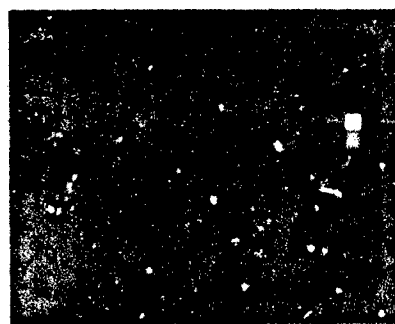

Simultaneous dual label image of chick myoblasts stained with propidium iodide and fluorescent anti-cytosine antibody. (Courtesy of Dr. Lettess Quinn, Univ. of Washington School of Medicine.)

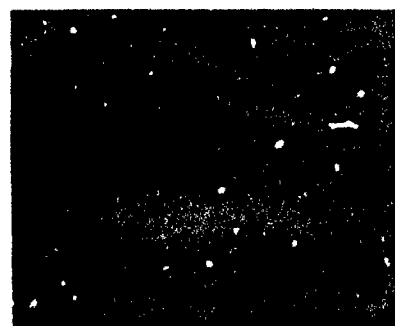

Detection of 3000 copies of ABL mRNA by in-situ hybridization in the leukemia cell line K562. (Courtesy of Dr. Joel Fluorine MD, Anderson Hospital.)

REASON • THE ACAS REVEALS THAN ME THE EYE

The ACAS uses a focused laser spot & advanced photometric detection, to sensitive, accurate fluorescence mea-

With the ACAS you can:

- Minimize interference with cell functions by using lower dye concentrations
- Detect fluorescence even when small number of sites are labeled
- Perform multiple measurement cell with minimal photobleaching
- Simultaneously measure two or with the dual detector option

Photomultiplier tubes provide the signal gain, linear response, and large dynamic range not found in video-based systems. Under most experimental conditions, the ACAS 570 provides you with reliable data from cells with minimal photobleaching, without the need for signal averaging.

With the ACAS 570, the Only is Your Imagination.

From cloning cells to rapid calcium measurements, the ACAS Interactive Cytometer gives you more information about your cells than pretty pictures.

Designed by cell researchers, the ACAS integrates laser illumination, photo detection, and comprehensive, interactive software, for the broadest range of applications in cell and molecular biology available in a single instrument.

Options such as quantitative confocal video processing, and dual photomultiplier offer you the flexibility and versatility to accommodate present and future needs.

The ACAS also comes with the Meridian commitment to total customer support, assuring maximum performance and productivity from your instrument. Dedicated staff of cell biologists, engineers, and computer scientists are immediately available with a phone call to answer your questions and instrumentation questions.

These are a few of the many reasons the ACAS 570 is **The Cytometer For 1000 Reasons.**

To discover **YOURS**, call us at 1-800-451-7200. In Michigan call 517-349-7200.

MERIDIAN™
Instruments, Inc.

The Cytometer for All Reasons

Tunable Argon-ion Laser

Photometric Detection

Acousto-Optic Modulator (AOM)

Dual Monitors

Inverted Fluorescence
Microscope

80386 Computer System

Keyboard and Mouse Interface

Motorized X-Y Stage

your

10-10-68

1. NAME
 2. ADDRESS
 3. CITY
 4. STATE
 5. ZIP
 6. PHONE
 7. TELETYPE
 8. FAX
 9. EMAIL
 10. DATE
 11. TIME
 12. LOCATION
 13. REMARKS
 14. INITIALS
 15. SIGNATURE
 16. DATE
 17. TIME
 18. LOCATION
 19. REMARKS
 20. INITIALS
 21. SIGNATURE
 22. DATE
 23. TIME
 24. LOCATION
 25. REMARKS
 26. INITIALS
 27. SIGNATURE
 28. DATE
 29. TIME
 30. LOCATION
 31. REMARKS
 32. INITIALS
 33. SIGNATURE
 34. DATE
 35. TIME
 36. LOCATION
 37. REMARKS
 38. INITIALS
 39. SIGNATURE
 40. DATE
 41. TIME
 42. LOCATION
 43. REMARKS
 44. INITIALS
 45. SIGNATURE
 46. DATE
 47. TIME
 48. LOCATION
 49. REMARKS
 50. INITIALS
 51. SIGNATURE
 52. DATE
 53. TIME
 54. LOCATION
 55. REMARKS
 56. INITIALS
 57. SIGNATURE
 58. DATE
 59. TIME
 60. LOCATION
 61. REMARKS
 62. INITIALS
 63. SIGNATURE
 64. DATE
 65. TIME
 66. LOCATION
 67. REMARKS
 68. INITIALS
 69. SIGNATURE
 70. DATE
 71. TIME
 72. LOCATION
 73. REMARKS
 74. INITIALS
 75. SIGNATURE
 76. DATE
 77. TIME
 78. LOCATION
 79. REMARKS
 80. INITIALS
 81. SIGNATURE
 82. DATE
 83. TIME
 84. LOCATION
 85. REMARKS
 86. INITIALS
 87. SIGNATURE
 88. DATE
 89. TIME
 90. LOCATION
 91. REMARKS
 92. INITIALS
 93. SIGNATURE
 94. DATE
 95. TIME
 96. LOCATION
 97. REMARKS
 98. INITIALS
 99. SIGNATURE
 100. DATE
 101. TIME
 102. LOCATION
 103. REMARKS
 104. INITIALS
 105. SIGNATURE
 106. DATE
 107. TIME
 108. LOCATION
 109. REMARKS
 110. INITIALS
 111. SIGNATURE
 112. DATE
 113. TIME
 114. LOCATION
 115. REMARKS
 116. INITIALS
 117. SIGNATURE
 118. DATE
 119. TIME
 120. LOCATION
 121. REMARKS
 122. INITIALS
 123. SIGNATURE
 124. DATE
 125. TIME
 126. LOCATION
 127. REMARKS
 128. INITIALS
 129. SIGNATURE
 130. DATE
 131. TIME
 132. LOCATION
 133. REMARKS
 134. INITIALS
 135. SIGNATURE
 136. DATE
 137. TIME
 138. LOCATION
 139. REMARKS
 140. INITIALS
 141. SIGNATURE
 142. DATE
 143. TIME
 144. LOCATION
 145. REMARKS
 146. INITIALS
 147. SIGNATURE
 148. DATE
 149. TIME
 150. LOCATION
 151. REMARKS
 152. INITIALS
 153. SIGNATURE
 154. DATE
 155. TIME
 156. LOCATION
 157. REMARKS
 158. INITIALS
 159. SIGNATURE
 160. DATE
 161. TIME
 162. LOCATION
 163. REMARKS
 164. INITIALS
 165. SIGNATURE
 166. DATE
 167. TIME
 168. LOCATION
 169. REMARKS
 170. INITIALS
 171. SIGNATURE
 172. DATE
 173. TIME
 174. LOCATION
 175. REMARKS
 176. INITIALS
 177. SIGNATURE
 178. DATE
 179. TIME
 180. LOCATION
 181. REMARKS
 182. INITIALS
 183. SIGNATURE
 184. DATE
 185. TIME
 186. LOCATION
 187. REMARKS
 188. INITIALS
 189. SIGNATURE
 190. DATE
 191. TIME
 192. LOCATION
 193. REMARKS
 194. INITIALS
 195. SIGNATURE
 196. DATE
 197. TIME
 198. LOCATION
 199. REMARKS
 200. INITIALS
 201. SIGNATURE
 202. DATE
 203. TIME
 204. LOCATION
 205. REMARKS
 206. INITIALS
 207. SIGNATURE
 208. DATE
 209. TIME
 210. LOCATION
 211. REMARKS
 212. INITIALS
 213. SIGNATURE
 214. DATE
 215. TIME
 216. LOCATION
 217. REMARKS
 218. INITIALS
 219. SIGNATURE
 220. DATE
 221. TIME
 222. LOCATION
 223. REMARKS
 224. INITIALS
 225. SIGNATURE
 226. DATE
 227. TIME
 228. LOCATION
 229. REMARKS
 230. INITIALS
 231. SIGNATURE
 232. DATE
 233. TIME
 234. LOCATION
 235. REMARKS
 236. INITIALS
 237. SIGNATURE
 238. DATE
 239. TIME
 240. LOCATION
 241. REMARKS
 242. INITIALS
 243. SIGNATURE
 244. DATE
 245. TIME
 246. LOCATION
 247. REMARKS
 248. INITIALS
 249. SIGNATURE
 250. DATE

...the ...
...the ...
...the ...
...the ...
...the ...

100-375

50) Negative

1. The first step is to identify the problem or question that needs to be answered. This involves understanding the context and the specific requirements of the task.

1. CONFIDENTIAL (If the information is not confidential, delete this marking.)

4. 5/27/68 - 1st processing work done

* Over 1000 pieces and graphics packing.

2000/05/25 07:51

What is a registered trademark of The McGraw-Hill Companies.
All specifications are subject to change without notice.
The McGraw-Hill Companies and its subsidiaries are not responsible for any errors or omissions in this publication.

MERIDIAN™

Instruments, Inc.

1360 Science Parkway • Okemos, Michigan 48864

027 349-7239 • 800-467-8339

REF ID: A66544

Calcium Measurements in Single Cells with the ACAS Interactive Laser Cytometer

Margaret H. Wade, Ph.D. and Susan A. McQuiston, MT(ASCP), Meridian Instruments, Inc.

Intracellular free calcium in single cells can be quantitated using the ACAS Interactive Laser Cytometer. Small changes in free calcium can have dramatic effects on many cellular processes including cell motility, muscle contraction, gap junction function, fertilization and cell division (1-5).

Free intracellular calcium levels are typically between 50 and 150 nM. These levels are approximately 20,000 times lower than the total calcium concentration estimated to be bound in mammalian cells or in the external milieu. Most of the calcium within a cell is sequestered within nuclei, mitochondria, and endoplasmic and sarcoplasmic reticuli. Multiple mechanisms in the cell maintain this enormous concentration gradient. Other mechanisms are responsible for releasing free calcium ions into the cytosol in response to a variety of environmental stimuli (5,6).

Calcium ions are therefore regarded as "one of the most important intracellular messengers in biology" (7,8). Since the changes in free calcium levels which trigger cell activity are minute and often very rapid, quantitating these changes has been technically very difficult. Recently, however, Roger Tsien and his associates (9) have developed a fluorescent calcium indicator, Indo-1, which is particularly well suited for quantitating rapid calcium flux in single cells.

The fluorescent probe, Indo-1, is used to monitor calcium in its bound and free states by exciting the dye with the 351-363 nm lines of the argon laser, and simultaneously detecting two emissions at 405 nm and 485 nm (9). A ratio of these emissions is computed automatically by the integrated software and allows quantitation of free calcium within the cell. The ratio method eliminates most of the problems associated with dye leakage and concentration of fluorescent probe within the cell, and the simultaneous detection provides for extremely accurate quantitation of temporal and spatial calcium changes. The laser excitation source provides a small beam (about 1 micron in diameter) which both allows for extremely low levels of fluorescence to be detected with minimum laser power, and greatly reduces the amount of photobleaching that might occur with repeated scanning. The X-Y scanning stage allows for .25

micron data acquisition to provide maximum spatial analysis.

Several types of experimental data can be collected using the ACAS. Single point data, integrated cross sectional line scans as well as two dimensional spatial scans can be measured within a cell.

Calcium is quantitated directly in cells by comparing the ratio of the emissions to that generated in a standard curve. The standard curve is produced by mixing the Indo-1 free acid with varying amounts of added calcium in a physiological buffer containing EGTA. A typical standard curve is shown in Figure 1 where the Detector 2/Detector 1 ratio is plotted against free calcium concentration.

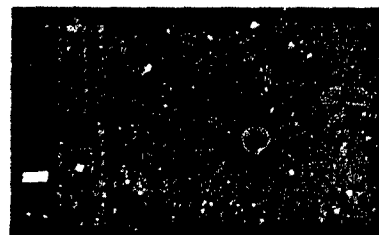


Figure 1. Standard curve of the ratio of the two emissions generated by adding calcium to EGTA containing solution of Indo-1 free acid.

Examples of line and image scan data are presented using serum starved human teratocarcinoma cells (HT) labelled with 1 μ M Indo-1 AM, the cell permeant form of the dye. The addition of fetal calf serum to these cells causes a sharp 2-3 fold increase in calcium concentration which then returns to a baseline level. The particular component(s) in serum responsible for the calcium spike have not yet been identified.

Figures 2 and 3 are typical examples of line scan data. Figure 2 represents one of the consecutive scans across the cell. Figure 3 is a ratio plot of the integrated value under each scan as a function of time. The vertical mark denotes the time when serum was added. The temporal shift has occurred very rapidly and the calcium levels return to baseline.



Figure 2. Line scan data showing scans across the cell with fluorescence at both 405 and 485 displayed.

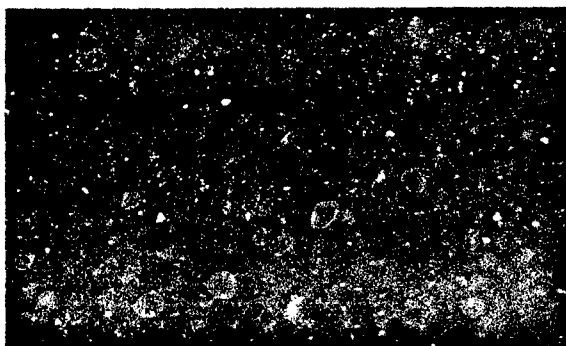


Figure 3. Ratio calculated from the integrated values determined from 400 line scans.

Figure 4-6 show an example of the 2-dimensional image scans that can be generated in a timed mode. The relative fluorescence intensity at both emissions in the first of twenty scans is plotted using the color scale denoted in the center of Figure 4. In Figure 5, serum was added, and by using the standard curve, the ratio of Detector 2 over Detector 1 is plotted as calcium values (calcium map). This map can be produced for each of the 20 images generated. In addition, a plot of the calcium concentration versus time in the designated cells can be produced as in Figure 6.

Use of Indo-1 provides high sensitivity of intracellular free calcium concentration. Controls should be used to insure that the fluorescence detected is calcium sensitive (9-12), including manganese quenching of the Indo-1 free acid, assessment of any Indo-1 ester fluorescence and polarity adjustments to the standard curve (13).

Quantitation of transient calcium values can easily be obtained in single cells using the ACAS Interactive Laser Cytometer. Spatial and temporal resolution of calcium can be analyzed using multiple scanning modes. Production of a standard curve allows for direct quantitation of calcium within cells.

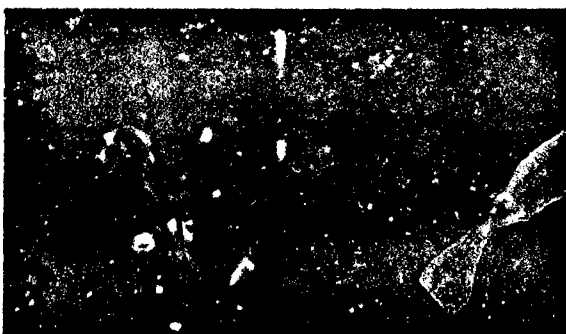


Figure 4 Two dimensional plot of data showing both emissions for a number of single cells.

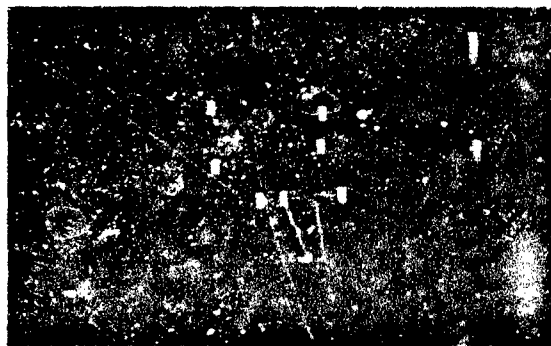


Figure 5. A plot of the ratio of detector 2/detector 1 after addition of serum.


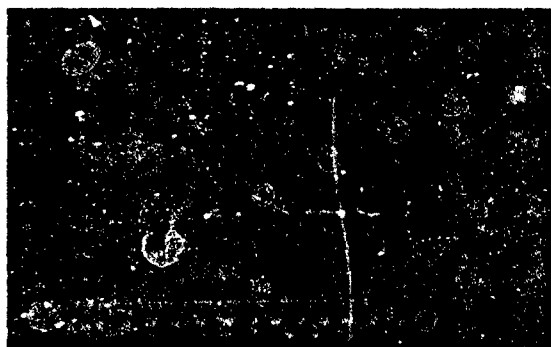


Figure 6. Time plot of the ratio values expressed as calcium concentration from the cells in Figure 5.

REFERENCES


1. Blinks, J. R., W. G. Wier, P. Hess and F. G. Prendergast (1982) *Prog. Biophys. Molec. Biol.* 40:1-114.
2. Poenie, M., J. Alderton, R. Steinhardt and R. Tsien (1986) *Science* 22 886-889.
3. Tsien, R. and M. Poenie (1986) *Trends in Biochem. Sci.* 11:450-455.
4. Keith, C. H., R. Ratan, R. F. Maxfield, A. Bajer and M. L. Shelanski (1985) *Nature* 316 848-850.
5. A. K. Campbell (1983) In, *Intracellular Calcium*, John Wiley and Sons, New York.
6. Carafoli, E. and J. T. Penniston (1985) *Sci. Amer.* 253:70-78.
7. Berridge, M. J. and R. F. Irvine (1984) *Nature* 312:315-321.
8. Nishizuka, Y. (1984) *Nature* 308 693-698.
9. Grynkiewicz, G., M. Poenie and R. Y. Tsien (1985) *J. of Biol. Chem.* 260:3440-3450.
10. Chused, T. M., A. Wilson, D. Greenblatt, Y. Ishida, L. J. Edison, R. Y. Tsien and R. D. Finkleman (1987) *Cytometry* 8:396-404.
11. Luckoff, A. (1986) *Cell Calcium* 7:233-248.
12. Scanlon, M., D. A. Williams and F. S. Fay, (1987) *J. Biol. Chem.* 262 6303-6312.
13. Popov, E. G., I. Y. Gavrilov, E. Y. Posin and A. Z. Gabbasov (1988) *Arch. of Biochem. and Biophys.* 262:91-96.

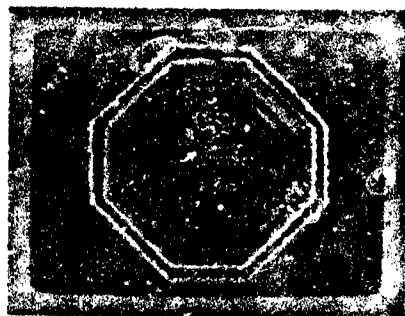


Outgrowth of CHO cell isolated by the "Cookie Cutter" technique. (Courtesy of Dr. Margaret Wade, Meridian Instruments, Inc.)



Fluorescence photobleaching of carboxyfluorescein diacetate labeled human teratocarcinoma cells as a measurement of cell-cell communication. Polygon enclosed cells are bleached and monitored for fluorescence recovery. (Courtesy of Dr. Margaret Wade, Meridian Instruments, Inc.)

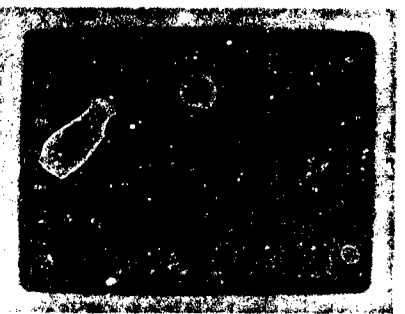

FRAP analysis of human fibroblasts stained with NBD phosphatidylcholine. Graph shows recovery of fluorescence following spot bleaching of the disc. (Courtesy of Dr. Margaret Wade, Meridian Instruments, Inc.)

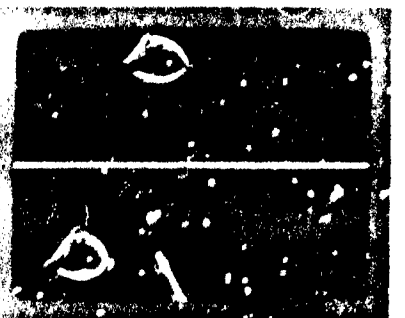

Laser transection of PC-12 neurites. (Courtesy of Dr. Steven Heidermann, Michigan State Univ. Dept. of Physiology.)

ACAS 570...

The Cytometer for All Reasons

Reason #1 INTERACTIVE LASER CYTOMETRY

The ACAS 570 Interactive Laser Cytometer lets you go beyond observation to Interaction with your cells in culture.

You command the power of the laser for new research capabilities including interactive techniques such as rare event selection, microsurgery, and fluorescence photobleaching.

CELL SELECTION IN CULTURE

The ACAS is the only instrument that lets you select cells directly from adherent cell cultures for isolation, propagation, and analysis.

• Sorting Cells the ACAS Way

Sort adherent cells without removing them from culture, using the unique, patented "Cookie-Cutter"™ process. This protocol automatically isolates a single cell based on analytical measurements, fluorescence, or morphology from a population growing on a film-lined dish.

A high intensity laser beam encircles the selected cells, leaving octagon-shaped "cookies" to the culture dish. When the film is peeled off the dish, the selected cells remain attached and continue to propagate.

This powerful **Interactive** technique is ideal when working with small quantities of cells and rare cells such as mutants, transfected cells and hybridomas, even at frequencies of one in a million.

And with "cell-list" memory, your selected cell locations are stored for repeated scans of the same cell or population of cells over time.

• Laser Ablation

Elate cell sub-populations using automatic laser irradiation to destroy either fluorescent or non-fluorescent cells, leaving desired cells

undisturbed for continued growth as a homogeneous culture.

Cells are grown on specially treated glass cover slips which absorb high intensity laser light to produce localized heating and killing of only the target cells.

Cell Selection applications include:

- Mutant selection based on physical and biochemical characteristics including membrane protein lateral mobility, gap junction communication, and transfection
- Selection and cloning of transformed cell lines such as smooth muscle, Chinese hamster ovary, and human teratocarcinoma stem cells
- Isolation of subpopulations of human melanocytes, T cells, and parenchymal hepatocytes

CELLULAR MICROSURGERY

With the ACAS 570 you can use the laser as a "photon scalpel" to perform cell surgery on individual cells, create transient pores in cell membranes, and ablate sub-cellular organelles such as mitochondria and lysosomes. Laser surgery can be coupled with quantitative fluorescence analysis to provide insight into cellular structure and function.

FLUORESCENCE PHOTOBLEACHING


The ACAS puts photobleaching under your control. Use the Laser as a "photon gun" to photobleach or photoactivate for real-time measurements of molecular dynamics.

• Fluorescence Redistribution after Photobleaching (FRAP)

FRAP techniques involve photobleaching an area or volume of a fluorescently labeled cell with a high intensity laser pulse. The movement of unbleached fluorescent molecules from surrounding areas re-establishes fluorescence in the bleached zone.



Subcellular localization of protein kinase C in neuroblastoma-glioma (NG) cells using a fluorescent phorbol ester. (Courtesy of Dr. H. A. Pershad Singh, Univ. of California, San Francisco)

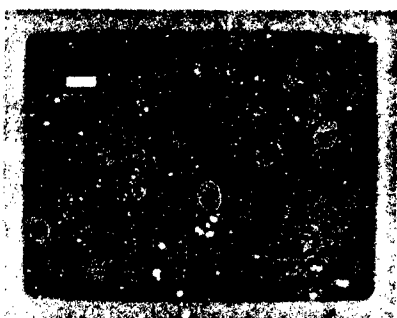

Analysis of the distribution of subcellular protein kinase C in NG cells labeled with a fluorescent phorbol ester. Polygon enclosed areas are monitored for distribution changes. (Courtesy of Dr. H. A. Pershad Singh, Univ. of California, San Francisco)

With the ACAS, you can directly measure molecular diffusion or flux, and degree of recovery, and use these rate determinations to explore cellular architecture and mechanics, including cytoskeleton organization, nuclear membrane structure, macromolecular assembly, and plant cell wall porosity.

• Cell-Cell Communication

Cell-cell communication techniques also rely on fluorescence photobleaching to measure intercellular communication between contacting plant or animal cells. Gap junction-mediated intercellular communication in animal cells is thought to play an important role in the regulation of cell proliferation and differentiation, providing a potential screening assay for environmental toxicants and drugs.

Use ACAS FRAP techniques to:

- Measure the lateral mobility of cell surface proteins
- Determine membrane potential
- Quantitate protein kinase C distribution and translocation
- Measure the modulation of gap junction-mediated cell-cell communication by tumor promoters, growth factors, and intracellular ions

Add the dimension of **Interaction** to your cell biology research. The ACAS 570 is the tool of choice for cell selection and fluorescence measurements in culture.

That's only one reason why the ACAS 570 is the **Cytometer For All Reasons.**

Find out more, call us at 800 267-5854, or Michigan call 517-540-7200.

MERIDIAN™


Instruments, Inc.

2310 Science Parkway
Okemos, Michigan 48864
TAX (517) 540-5967

© 1989 Meridian Instruments, Inc.




Calcium oscillations in freshly isolated rat aorta cells labeled with the fluorescent calcium probe Indo-1 (Courtesy of Dr. Oliver Flin, Bristol Myers Co.)

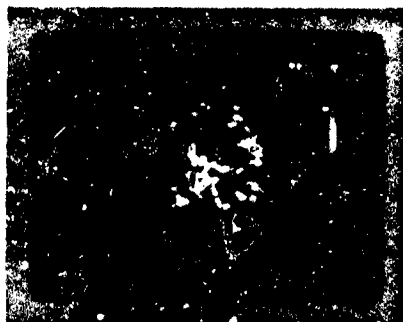

HIV-1RT infected CEM-SS cell stained by indirect immunofluorescence using anti- α p120 (Courtesy of Dr. Jonathan Warren, National Institutes of Health)



Cholesterol ester accumulation in human macrophages stained with Nile red from patients with coronary artery disease (Courtesy of Dr. Eugene Koren, Chukashima Medical Research Foundation)

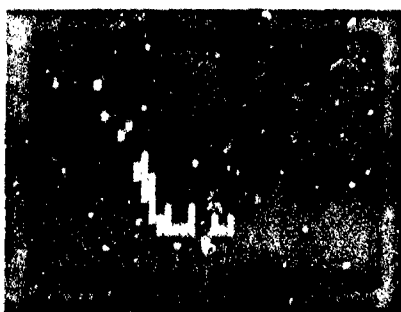

Histogram of cholesterol ester accumulation in human macrophages stained with Nile red (Courtesy of Dr. Eugene Koren, Chukashima Medical Research Foundation)

ACAS 570... The Cytometer for All Reasons

Reason #2 QUANTITATIVE FLUORESCENCE MEASUREMENTS

The ACAS 570 Interactive Laser Cytometer unleashes the analytical power of fluorescence, giving you reproducible **Quantitative** measurements of cellular structure, function, and response.

From single and multiple measurements, to simultaneous emission ratios, you get **real numbers** for fluorescence quantitation ranging from DNA content and molecular diffusion, to intracellular Ca^{2+} and pH.

QUANTITATIVE FLUORESCENCE APPLICATIONS

Physical and Biochemical Measurements

Using two dimensional fluorescence analysis you can easily identify and quantitate content, composition, and distribution of cellular components such as lysosomes, mitochondria, DNA/RNA, structural proteins, enzymes, and receptor molecules, in one cell or hundreds of cells. In addition, you can measure a variety of biochemical activities such as membrane potential, oxidation states, and ligand binding.

Current applications include:

- Quantitation of DNA, fluorescence and repair at the single cell level
- Dynamic *in situ* distribution of protein kinase C using fluorescent phorbol esters
- Monoclonal marker identification of human tumor cells in mixed cultures
- *In situ* hybridization

Calcium, pH and other Intracellular Ions.

With the ACAS you can rapidly measure single, multiple or simultaneous dual emissions in a point, line or two dimensional image scan using fluorescent dyes such as SNARE-1TM, Indo-1 and BCE-1. Powerful software helps you generate a

standard curve, and calculates ratio or dual label temporal and spatial concentration gradients.

Quantitative fluorescence measurements include:

- Cytosolic Ca^{2+} response to a variety of chemicals, including tumor promoters, growth factors, and hormones
- Repetitive transient Ca^{2+} oscillations in single heart muscle cells
- Simultaneous Ca^{2+} and pH measurements using the visible probes Fluo-3 and SNARE-1

Fluorescence Redistribution after Photobleaching (FRAP)

The ACAS puts photobleaching under your control. Use the unique tool of fluorescence photobleaching to create optical gradients at specific intracellular sites or in whole cells for real-time measurements of cellular dynamics and patterns.

FRAP techniques directly measure molecular diffusion or flux, and degree of recovery, important parameters in understanding cellular mechanics and architecture.

ACAS FRAP techniques measure:


- Lateral mobility of membrane proteins
- Plasmal cell wall pore diameter
- Plasmic membrane macromolecular mobility rates

Cell-Cell Communication


The ACAS is also used to measure intercellular communication between contiguous cells in animal and plant tissue.

This protocol, relying on fluorescence photobleaching, measures molecular transfer through cell gap junctions, gaining important insights into the regulation of cell proliferation and differentiation, and is a potential screening assay for environmental toxicants and drugs.

ACAS 570 is a registered trademark of

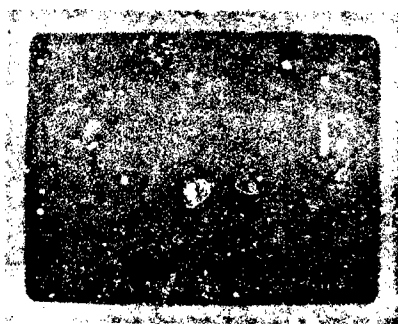

Detection of 16-41 photoproducts in UV-irradiated VP-1 cells by a monoclonal antibody (UV dose = 0.1 J/m²). Courtesy of Dr. Toshio Mori, Kanazawa University, Japan.

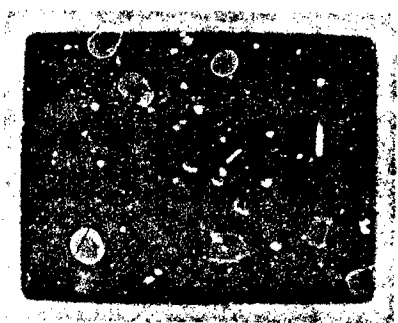

Same as above except UV dose equals 2.0 J/m².

Intercellular communication measurements include:

- Inhibition of gap-junction mediated intercellular communication by tumor promoters and growth factors
- Modulation of gap-junction function by intracellular levels of Ca²⁺, pH, and cAMP
- Control of transport between plant cells

• Quantitative Immunofluorescence

Immunofluorescence techniques can be used for sensitive, rapid screening assays on the ACAS. You can monitor antigen expression, cell binding and killing, as well as quantitative morphological information (i.e., mean fluorescence and cell size) for sensitive, non-subjective measurements.

Now you can use the ACAS for:

- Screening assays of anti-HIV compounds
- Cytotoxicity assays of NK cell binding to tumor cell targets
- Quantitative screening assays for antigen expression in tumor tissue sections

All the dimension of **quantitation** in your cell biology research. The ACAS 570 is the tool for rapid, accurate measurement of molecular concentration, distribution, and dynamics in cultured cells.

It's no wonder, therefore, the ACAS 570 is the **Cytometer For All Reasons.**

For a demonstration, call us at 800-247-8884. In Michigan, call 817-5549-7200.

MERIDIAN™

Instruments, Inc.

2310 Science Parkway
Okemos, Michigan 48864
FAX (517) 461-5065

ACAS 570...

The Cytometer for All Reasons



The ACAS Interactive Laser Cytometer is *the* instrument for quantitative fluorescence analysis and adherent cell sorting.

The power of the ACAS lies in the comprehensive, interactive software that turns laser cytometry with advanced computer technology and a fast, flexible operating system — giving you integrated data acquisition, analysis and management capabilities that are easy to learn and easy to use.

APPLICATIONS BASED PROGRAMS FOR THE CELL BIOLOGIST

You get a powerful assortment of flexible programs with applications ranging from cell and molecular biology to toxicology, immunology, and neurobiology.

- **Image Analyze** — Two dimensional fluorescence quantitation of single or dual emissions.
- **Auto Image** — Automatic scanning and fluorescence quantitation of a large number of the same field over time.
- **Ratio Analysis** — Quantitation of intracellular Ca²⁺, pH or other ratios.
- **Cell Sorting** — Automatic sorting of cells by laser addition.
- **FRAP** — Measurement of macromolecular diffusion.
- **Cell-Cell** — Measurement of inter-cellular communications via gap junctions, or general kinetic measurements.
- **Cell List** — Store locations of selected cells for automated measurements, screening, cell selection, and re-sorting.

SIMPLIFIED SYSTEM INTERACTION

Designed for simple operation right from the start, the intuitive, menu-driven ACAS software helps you learn the system easily, and become productive quickly. For example:

- Perform all operations from the logically organized main menu using easy pull-down menu functions.
- Use the convenient mouse or keyboard, working in parallel, for fast, simple interaction with instrument hardware *and* software.
- Follow the screen prompts and menu display for instant selection of your next step.
- Custom tailor directions, analysis parameters, palette colors, graphics and many other programs with the integrated editor functions.

TOTAL CONTROL OVER YOUR EXPERIMENT

The ACAS software makes it easy to design and conduct your experiment directly from the screen by integrating optical, mechanical and computer hardware functions.

With the touch of a finger you can:

- Move the mouse arrow to easily select experimental parameters from the program menu.
- Manually control instrument functions, such as, movement of the scanning stage, and laser intensity and duration.
- Perform instrument calibration and system diagnostics from the computer screen.
- Automatically store cell locations for repeated scans of the same cells over time.

COMPREHENSIVE ANALYSIS FUNCTIONS

Use the mouse driven, pull down menus to automatically perform a variety of image analysis and processing functions including:

- Threshold, threshold, and rescale data.
- Define cells with polygons for separate analysis.
- Eliminate cells for exclusion from final calculations.
- Query data from any area or at a specific point in an image.
- Create histograms based on cellular structure, for example, number of fluorescent cells.

average fluorescence, area, perimeter and shape factors) integrating hundreds of experiments.

- Choose data files and preset parameters for automatic analysis.

POWERFUL GRAPHICS CAPABILITIES

Create a variety of analytical displays with the powerful graphics functions for accurate, in-depth visual interpretation of your data. For example:

- Display fluorescence scans as 256 gray-levels or various pseudo-color palettes.
- Illustrate population differences using multiparameter scatterplots and histograms.
- Overlay dual detector images for visual comparison and discrimination of fluorescence values from each detector.
- Create composite displays of up to four screens for easy comparison of fluorescence images and corresponding graphics.
- Automatically create a graph of the fluorescence values along any specified cross section of a scan.

EASY, EFFICIENT DATA MANAGEMENT

The uniquely designed ACAS software can input, store and retrieve your data with speed and flexibility.

- Perform disk and file maintenance quickly and easily with the dedicated ACAS file utility system.
- Use the mouse or keyboard to tag files for group or individual operations.
- Automatically locate specific files on any disk.
- MS-DOS operating system ensures compatibility with a wide range of other available software.
- Easily integrated into networks via Ethernet*.

*Ethernet is a trademark of Digital Equipment Corporation.


LET'S DO AN EXPERIMENT

The ease and simplicity of the powerful ACAS software can be shown best by following an experiment

So let's examine the kinetics of calcium in living cells

The ACAS 570 can be used to measure intracellular calcium in single cells or groups of cells. The dual wavelength **Simultaneous Detection System** provides ratiometric measurements for a variety of experiments. Single point scans, integrated cross sectional line scans, and two dimensional spatial scans can be measured

The integrated software automatically calculates intracellular calcium by comparing the emission ratios of the free and bound forms of a fluorescent dye such as fura-1 to the ratios in a standard curve

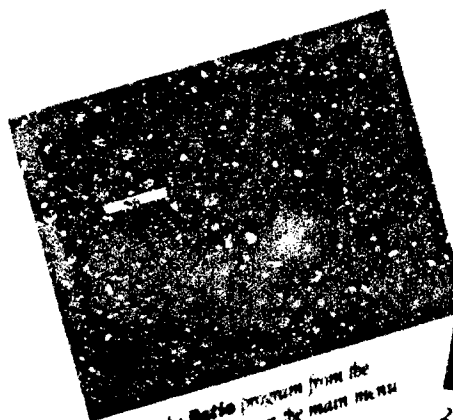


1 Access the **Ratio** program from the **Analysis** options on the main menu



2 Choose the **Calcium** option



3 You have the **Image Scan** will choose the

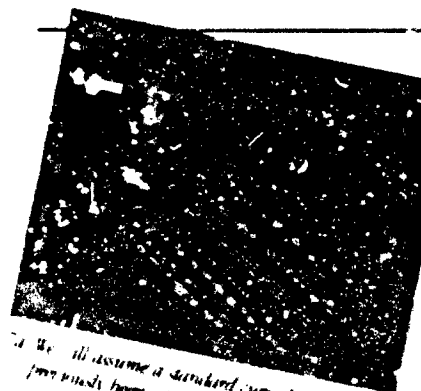


4 We will assume a standard curve has previously been established




5 Enter the curve name with the **Select Curve** option



6 The **Operation** contains a **Curve** menu option

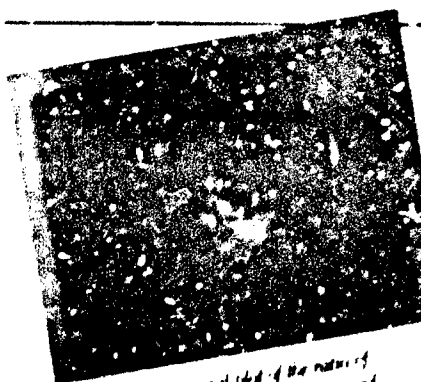

7 A two dimensional plot of the ratio of Ca^{2+} distribution can be displayed using the **Image** option. This can also be **Exported** from the map

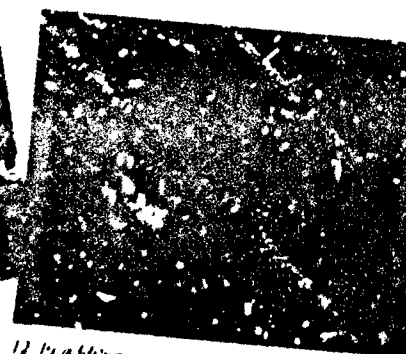


8 To achieve a pixel distribution of the Ca^{2+} distribution as seen on the map, can be **Exported** automatically




9 A graph of the change in Ca^{2+} distribution over time for the marked areas can also be automatically generated by accessing **Time Plot** from the **Operations** menu



10 The **Show** menu for **Time Plot** can be used to **Plot** the Ca^{2+} data

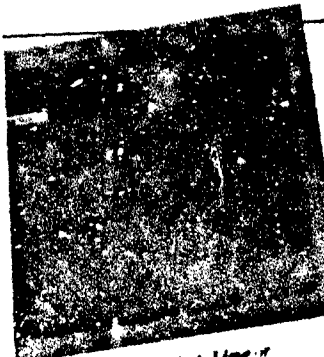

Save the choice of a **Point Line or Image Scan** for this experimental use. Choose the **Image Scan**.

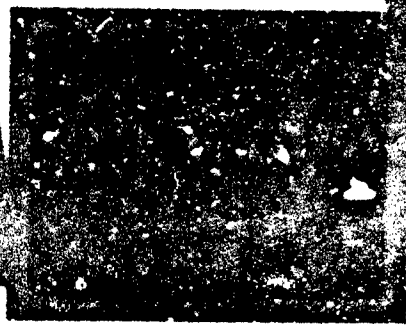

4. Begin the experiment by assigning a file name and setting instrument parameters for the scan from the options on the screen.



5. Cells loaded with a fluorescent calcium probe are placed on the scanning stage. Use the mouse to manually move the stage and select a field of cells. Begin scanning. Add the autofocus after the baseline Ca^{2+} levels have been established. Continue scanning. Data is saved to the hard disk or removable cartridge disk.




6. To analyze this data, access **Image Analyze** from the **Analysis** program on the main menu and enter the file name.



Operations and **Plot** menu's offer a variety of options for data reduction.



8. We will first **Smooth** **Rescale** and **Threshold** data.

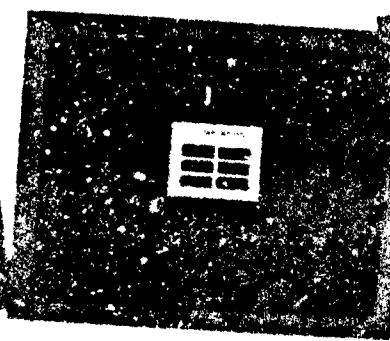


9. Use the mouse to mark specific cells with polygons for separate analysis.



10. Query raw data in several ways from the **Operations** menu. Use the mouse pointer to automatically report the calcium concentration within an entire scan area, within a polygon, in a cross-sectional line, or at a specific point.



12. Mouse provides a variety of options for data reduction.

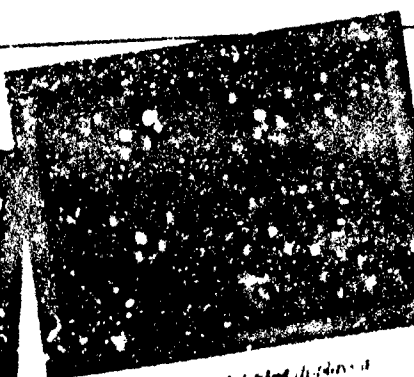

14. For example, **Multi plot** displays a summary of three scans and the time plot.




15. To end acquisition or analysis at anytime, use the **Control-Cancel** option always found on the menu.

As you can see, there's nothing h about our software. The unique design makes it easy to get the numbers you need for Calcium, and other cellular physical and biochemical measurements.




AUTOMATIC FILE INTEGRATION

Incorporating an innovative file utility system, the ACAS software automates data management for fast, simple file integration and statistical analysis. For example, you can easily integrate hundreds of individual experiments to create composite histograms and scatterplots illustrating a variety of analytical parameters.

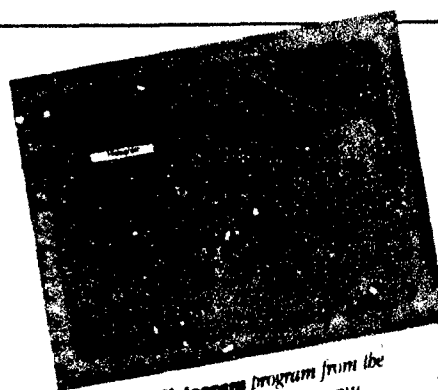


1. Enter the **Histogram** program from the **Analysis** option on the main menu.

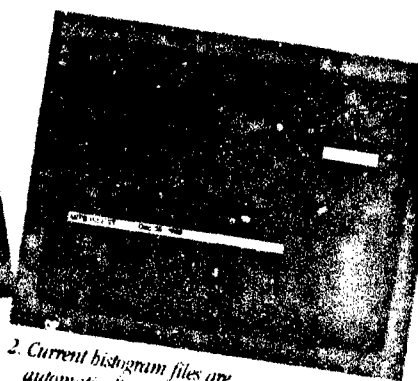


2. Current histogram files are automatically displayed. **Tag** the histogram files you want for integration and analysis.



3. Choose histogram display parameters from the options listed on the **Show** menu.

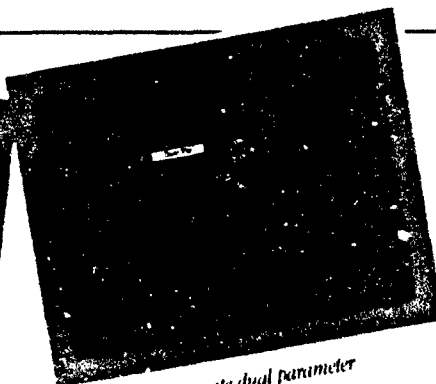

4a. You can also create dual parameter scatterplots.

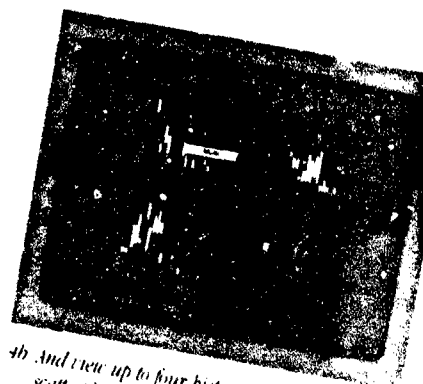


4b. And view up to four histograms and scatterplots simultaneously from the **Show** menu options.



5a. The **Operations** menu offers a variety of analysis and display options. For example, **Differentiate** assigns a specific color to individual files in the composite histogram for easy visual comparison.



5b. **Window** gives you statistics on up to four areas (windows) of a histogram.

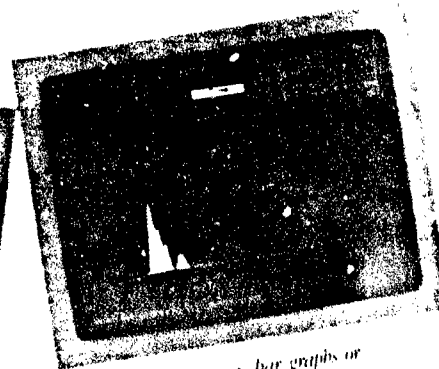

6. Display histograms as bar graphs or curved plots, using the **Plot** menu options. You can also export data to an ASCII format for further analysis.

These are just a few of the extensive data acquisition, analysis and management programs available with the **ACAS 570 Interactive Laser Cytometer**.

And, one more reason why the ACAS 570 is **The Cytometer For All Reasons**.


To see more, and find out how the ACAS can work for **you**, call us at 1-800-247-8084, or (517) 349-7200.

MERIDIAN
Instruments, Inc.

2410 Science Parkway • Okemos, MI 48864
517 349-7200 • 800 247-8084 • FAX 517 349-5967

DASY 9000—Data Analysis System

Extend and Enhance the Analysis Capabilities of Your ACAS

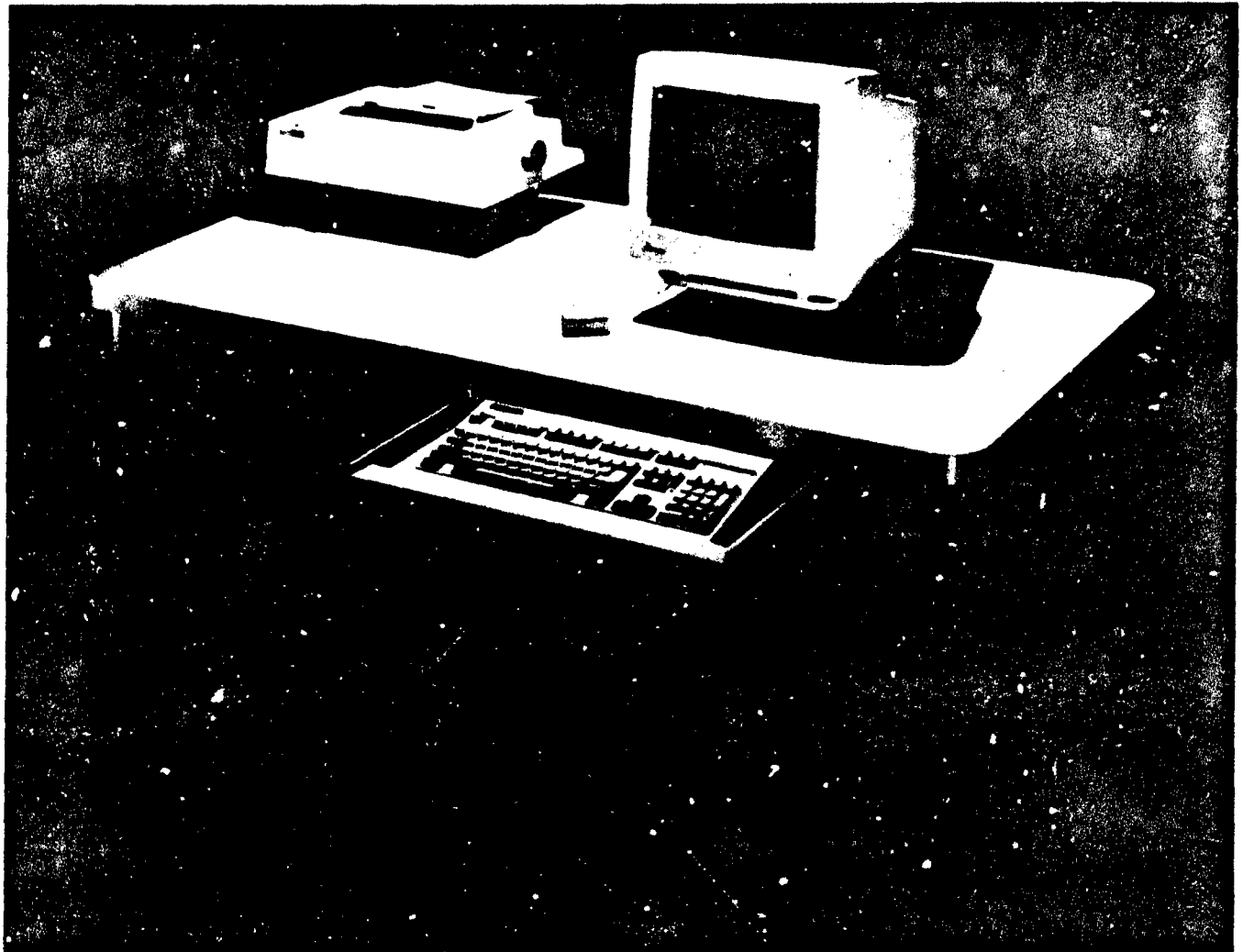

-
- Improves Laboratory Productivity
 - Duplicates ACAS 570 Analysis Capabilities
 - Provides Convenient Access to Data
 - Intuitive, Menu-Driven Software, Easy To Learn And Use
 - Convenient Mouse and Keyboard Control
 - Comprehensive Graphics Functions
 - High Resolution Gray-Scale or Dynamic Color Displays
 - High Quality Hard Copy Reproduction From Any Screen
 - MS-DOS Disk and File Maintenance
-

The **DASY 9000 — Data Analysis System** — is a complete data analysis workstation dedicated to the analysis and display of fluorescence data gathered on the ACAS 470 or 570 Interactive Laser Cytometer. As an independent computer workstation, the DASY provides convenient access to data, and rapid execution of analysis programs without interrupting data acquisition on the ACAS.

You also get all the power of the comprehensive versatile ACAS software, providing analysis, graphics and data management capabilities for a wide range of cell biology applications.

Now you can collect your data on removable, 20 Mbyte micro-bernoulli disks, or 1.2 Mbyte floppy disks for analysis on the DASY at any convenient time or location.

By relieving the ACAS of data analysis and reserving it for data acquisition, the DASY expands laboratory efficiency, saving you time and improving your productivity.

SPECIALIZED PROGRAMS FOR QUANTITATIVE FLUORESCENCE ANALYSIS

Specifically designed for cell and molecular biology applications, the DASY duplicates the powerful analytical capabilities of the ACAS, giving you accurate, in-depth analysis of *in vitro* adherent cell fluorescence measurements.

- Calculate fluorescence distribution within cells
- Quantitate cellular responses ratiometrically for Ca^{2+} , pH and other intracellular ions
- Calculate molecular diffusion across cell membranes (FRAP, Cell-Cell communication)
- Characterize cell populations based on physical and biochemical measurements



Two-dimensional pseudo-color fluorescent image of bladder tumor cell line MGH-U1 stained with rhodamine-123. Highly fluorescent areas represent regions containing mitochondria surrounding the less fluorescent nucleus. (Courtesy of Dr. Chi Wei Lin, Massachusetts General Hospital and Harvard Medical School)

VERSATILE SOFTWARE SIMPLIFIES INTERACTION

Designed for simple operation right from the start, the intuitive, menu-driven ACAS software helps you learn the system easily and become productive quickly.

All operations are performed with simple commands accessed by program selections from the main menu. Programs are logically organized and quickly chosen from pull-down menus using a keyboard or two-button mouse, working in parallel, to provide fast, easy system interaction.

For example: Use the mouse and pull-down menu operations to automatically perform a variety of image analysis functions including:

- Smooth and threshold data
- Mark selected areas for separate analysis
- Define cells with polygons for automated analysis
- Eliminate cells for exclusion from final calculations
- Calculate cellular statistics (i.e., number of cells, average cellular fluorescence, etc.)

POWERFUL ANALYTICAL GRAPHICS ILLUSTRATE FLUORESCENCE MEASUREMENTS

The comprehensive ACAS software gives you the capabilities to create detailed graphics for a variety of analytical displays.

- Automatically create histograms based on cellular statistics integrating hundreds of experiments
- Illustrate population differences (i.e., dual-detector scans) using multiparameter scatter plots, and bar or contour histograms
- Create composite displays of fluorescence images and corresponding graphics analysis for easy comparison
- Overlay dual detector images for visual comparison and discrimination of



Two-dimensional gray-scale fluorescent image of the bladder tumor cell line MGH-U1, stained with rhodamine-123.

- fluorescence values for each detector
- Query data from any area or at a specific point in an image

HIGH RESOLUTION COLOR IMAGES ENHANCE SIGNIFICANT DATA

Sharp, 640 × 480 pixel resolution gives you clear, vivid displays of contrasting cellular components and histogram populations.

Fluorescence scans are displayed as gray-scale or pseudo-color images composed from a variety of user defined color palettes.

High quality hard copies of any screen can be quickly obtained from the color printer, or through photographic reproduction via an optional in-line camera system for 35 mm slides and prints.

EASY, EFFICIENT DATA MANAGEMENT


Designed around a 32 bit, 80386 based microcomputer (IBM AT* compatible), the DASY can input, store, and retrieve your data with speed and flexibility.

Optimal program and data storage is achieved with a 60 Mbyte winchester hard disk and a 1.2 Mbyte floppy disk drive. (Additional storage capacity is available with an optional 20 Mbyte micro-bernoulli disk drive.)

• MS-DOS Operating System

By utilizing the widely accepted MS-DOS operating system, the ACAS software provides simple, easy to learn programs, automating data processing and analysis, and minimizing data management time. MS-DOS also permits the use of a variety of commercially available programs and easy integration into local area networks via Ethernet**.

To find out how the DASY 9000 can increase the analytical power of **your** ACAS, call us at 800-247-8084. In Michigan call 517-349-7200.

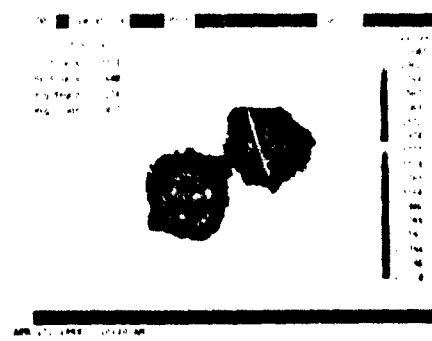

Color printer reproduction of two-dimensional pseudo-color fluorescent image of bladder tumor cell line MGH-U1 stained with rhodamine-123.

MERIDIAN™

Instruments, Inc.

2310 Science Parkway • Okemos, Michigan 48864

(517) 349-7200 • 800-247-8084

FAX (517) 349-5967

*IBM is a trademark of the International Business Machines Corp.

**Ethernet is a trademark of Digital Equipment Corp., Intel Corp., and Xerox Corp.

**ACAS INTERACTIVE LASER CYTOMETER™
BIBLIOGRAPHY**

AIDS

Abstracts

Warren, J.T. J.B. McMahon, O.S. Weislow, R. Kiser, and M.R. Boyd. National Cancer Institute, Frederick Cancer Research Facility (NCI-FCRF), Frederick, MD. (1988) **Temporal Expression of HIV-Associated Antigens Detected by Quantitative Immunofluorescence Using Interactive Laser Cytometry.**

CALCIUM

Burnatowska-Hledin, M.A. and W.S. Spielman. (1989) **Vasopressin V1 Receptors on the Principal Cells of the Rabbit Cortical Tubule: Stimulation of Cytosolic Free Calcium and Inositol Phosphate Production via Coupling to a Pertussis Toxin Substrate.** J. Clinical Investigation (in press).

Hosoi, Kazuo, and M. Edidin. (1989) **Exogenous ATP and Other Nucleoside Phosphates Modulate Epidermal Growth Factor Receptors of A-431 Epidermoid Carcinoma Cells.** Proc. Natl. Acad. Sci. 86: 4510-4514

Madhukar, B.V., S.Y. Oh, C.C. Chang, M.H. Wade and J.E. Trosko. (1989). **Altered Regulation of Intercellular Communication by Epidermal Growth Factor, Transforming Growth Factor-beta and Peptide Hormones in Normal Human Keratinocytes.** Carcinogenesis. 10(1):13-20.

Meridian Instruments, Inc. Application Notes

Wade, M.H., and S.A. McQuiston. (1988) **Calcium Measurements in Single Cells with the ACAS Interactive Laser Cytometer.**

Abstracts

Benjamins, J.A. and C.A. Dyer, Wayne State Univ., Dept. of Neurol., Detroit, MI (1989) **Glycolipids and Transmembrane Signaling in Oligodendroglia.** Nov. 16-18 meeting of New York Academy of Sciences on "Myelination & Demyelination".

Dyer, C.A., and J.A. Benjamins, Wayne State Univ., Dept. of Neurol., Detroit, MI. (1990) **Antibodies to Sulfatide and Galactocerebroside Increase Ca^{2+} in Oligodendroglia.** Am. Soc. Neurochem.

Dyer, C.A., J.A. Benjamins and S. Galoforo, Wayne State Univ., Dept. of Neurol., Detroit, MI. (1989) **Intracellular Ca^{2+} Levels in Oligodendroglia are Increased by Anti-galactocerebroside.** Nov. 16-18 meeting of New York Academy of Sciences on "Myelination & Demyelination".

Grillone, L.R., B.D. Jensen and S.T. Crooke. Smith Kline and French Laboratories, Philadelphia, PA. **Single Cell Calcium Measurements in A10 Smooth Muscle Cells.** Journal of Cell Biology, 105:186.

Jensen, B.D., L.R. Grillone and S.T. Crooke. Smith Kline and French Laboratories, Philadelphia, PA. **Low Concentrations of Vasopressin Induce Repetitive Transient Calcium Oscillations in Individual A10 Smooth Muscle Cells.** Society for Analytical Cytology. :74

Madhukar, B.V., S.Y. Oh, A.H. Jiwa* and M.H. Wade*. Michigan State University and *Meridian Instruments, Inc. (1989) **Inhibition of Intercellular Communication and Increase in Intracellular Calcium In Vitro by Xenobiotic Chemicals in a Liver Epithelial Cell Line.** In Vitro Toxicology.

Spiegel¹, S., S. McQuiston², and M.H. Wade³. ¹NIH, Bethesda, MD, ²University of California, San Francisco, CA and ³Meridian Instruments, Inc., Okemo, VT. (1988) **Stimulation of Quiescent Swiss 3T3 Cells with the B Subunit of Cholera Toxin Leads to Oscillations in Cytosolic Free Calcium.** FASEB

CELL-CELL COMMUNICATION & TUMOR PROMOTION ANALYSIS

Anders, J.J. Department of Anatomy, Uniformed Services University of the Health Sciences, Bethesda, Maryland. (1988) **Lactic Acid Inhibition of Gap Junctional Intercellular Communication in In Vitro Astrocytes as Measured by Fluorescence Recovery after Laser Photobleaching.** GLIA 1:371-379

Anders, J.J., M. Salopek, Department of Anatomy, Uniformed Services University of the Health Sciences, Bethesda, Maryland. (1989) **Meningeal Cells Increase in vitro Astrocytic Gap Junctional Communication as Measured by Fluorescence Recovery After Laser Photobleaching.** The Journal of Neurocytology. 18: 257-264

Anders, J.J., S. Neidermair, E. Ellis, M. Salopek, Dept. of Anatomy, Uniformed Services Univ. of the Health Sciences, Bethesda, Maryland. (1989) **The Response of Rat Cerebral Cortical astrocytes to Frex or Cobalt-induced Injury: An Immunocytochemical and Gap-FRAP Study.** GLISA

Baron-Epel, O., D. Hernandez, L.-W. Jiang, S. Meiners, and M. Schindler, Michigan State University, E. Lansing, MI. (1988) **Dynamic Continuity of Cytoplasmic and Membrane Compartments between Plant Cells.** Journal of Cell Biology. 106:715-721

Bombick¹, D.W., G. Zhang¹, and J.E. Trosko². ¹R.J. Reynolds Tobacco Co., Winston-Salem, NC and ²Michigan State University, E. Lansing, MI. (1987) **Use of Fluorescent Recovery After Photobleaching (FRAP) Analysis to Demonstrate Inhibition of Gap Junctional Intercellular Communication by Tumor Promoters.**

Bombick, D.W., R.J.R. Tobacco Co., Winston-Salem, NC. **Gap Junctional Communication in Various Cell Types after Chemical Exposure.** accepted for publication In Vitro Toxicology.

Evans, M.G., M.H. El-Fouly, J.E. Trosko, and S.D. Sleight. (1988) **Anchored Cell Analysis/Sorting Coupled with the Scrape-Loading/Dye Transfer Technique to Quantify Inhibition of Gap Junctional Intercellular Communication in WB-F344 Cells by 2,2', 4,4', 5,5'-Hexabromobiphenyl.** J. Toxicol. and Environ. Health. 24:261-271

Fitz, T.A., M.M. Marr, T.L. Walden, W.A. Schmidt and C.A. Winkel. (1989) **Gap Junctions in Transformed Rat Granulosa Cells.** Growth Factors and the Ovary, Plenum Publishing Corp. :245-250

Madhukar, B.V., S.Y. Oh, C.C. Chang, M.H. Wade and J.E. Trosko. (1989) **Altered Regulation of Intercellular Communication by Epidermal Growth Factor, Transforming Growth Factor-beta and Peptide Hormones in Normal Human Keratinocytes.** Carcinogenesis. 10(1):13-20.

Oh, S.Y., B.V. Madhukar and J.E. Trosko. (1988) **Inhibition of Gap Junction Blockage by Palmitoyl Carnitine and TMB-8 in a Rat Liver Epithelial Cell Line.** Carcinogenesis. 9(1):135-139.

Rezabek, M.S., J.E. Trosko, C. Jone, and S.D. Sleight. (1988) **Effects of Hepatic Tumor Promoters Phenobarbital and Polybrominated Biphenyls on Intercellular Communication Between Rat Liver Epithelial Cells.** In Vitro Toxicology. 2(1):45-58

Schindler, M., and L.-W. Jiang. (1987) **Fluorescence Redistribution After Photobleaching as a Tool for Dissecting the Control Elements of Nucleocytoplasmic Transport.** Methods in Enzymology B: Cellular Regulators (S.P. Kolowick, N.O. Kaplan, A.R. Mears, and T.M. Conn, Eds.), 141:447-459. Academic Press, New York.

Schindler, M., J.E. Trosko, and M.H. Wade. (1987) **Fluorescence Photobleaching Assay of Tumor Promotor 12-O-Tetradecanoylphorbol 13-Acetate Inhibition of Cell-Cell Communication.** Methods in Enzymology B: Cellular Regulators (S.P. Kolowick, N.O. Kaplan, A.R. Mears, and T.M. Conn, Eds.), 141:439-447. Academic Press, New York.

Suter, S., J.E. Trosko, M.H. El-Fouly, L.R. Lockwood and A. Koestner. (1987) **Dieldrin Inhibition of Gap Junctional Intercellular Communication in Rat Glial Cells as Measured by the Fluorescence Photobleaching and Scrape Loading/Dye Transfer Assays.** Fund. Appl. Toxicology 9:785-794.

Trosko, J.E., C.C. Chang, M.H. El-Fouly, R. Kulkarni, and R. Gera. Michigan State University, E. Lansing, MI. (1986) **Role of Gap Junctional Intercellular Communication in the Regulation of Growth and Differentiation of Normal and Malignant Human Epithelial Cells.** Fourth Decennial Review - Cell and Tissue Culture International Conference

Wade, M.H., J.E. Trosko, and M. Schindler. (1986). **A Fluorescence Photobleaching Assay of Gap Junction-Mediated Communication Between Human Cells.** Science 232:525-528.

Meridian Instruments, Inc. Application Notes

Trosko, J.E. (1986) **Rapid New Method to Study Gap Junction-Mediated Intercellular Communication in Mammalian Cells.**

Abstracts

- Anders, J.J. and S Woolery. Uniformed Services University of the Health Sciences, Bethesda, MD. **The Effects of Laser Injured Neurons on in vitro Astrocytic Gap Junctional Communication Measured by Fluorescence Recovery After Laser Photobleaching.** Anat. Rec., 223(4): 8A.
- Anders, J.J. Uniformed Services University of the Health Sciences, Bethesda, MD. (1988) **Effects of Lactic Acid on Astrocytic Gap Junctional Communication Measured by Fluorescence Recovery After Laser Photobleaching.** Lasers in Medicine and Surgery. 8(2): 147
- Anders, J.J. Uniformed Services University of the Health Sciences, Bethesda, MD. (1988) **Effects of Arachnoid-pial cells on In Vitro Astrocytic Gap Junctional Communication Measured by Fluorescence Recovery After Laser Photobleaching.** Journal of Neuroscience. 14(2): 1056
- Bombick, D.W. and J.E. Trosko. Michigan State University, East Lansing, MI. **Quercetin Reverses TPA Inhibition of Gap Junctional Intercellular Communication in Rat Glial Cells.**
- deFeyter¹, E., C.M. Weghorst², R. Rusch², J. Klaunig², C.C. Chang¹, B.V. Madhukar¹, M. El-Fouly¹, and J.E. Trosko¹. ¹Michigan State University, E. Lansing, MI and ²Medical College of Ohio, Toledo, OH. (1988) **Reduced Gap Junctional Intercellular Communication and Tumorigenicity of Rat Liver Epithelial ("Oval", F344-WB) Cells.** Toxicological Implications of Altered Gap Junctional Intercellular Communication.
- El-Fouly¹, M.H., J.E. Trosko¹, C. Nelson² and G.P. Dotto². ¹Michigan State University, E. Lansing, MI and ²Yale University, New Haven, CT. (1988). **Effects of Ca²⁺, TPA, and H-ras Oncogene on Gap Junction Intercellular Communication in Primary Epithelial Cells.** Toxicological Implications of Altered Gap Junctional Intercellular Communication.
- Evans, M.G. and J.E. Trosko. Michigan State University, E. Lansing, MI. (1987) **Dose-Dependent Inhibition of Metabolic Cooperation in Rat Epithelial (WB-F344) Cells by Hepatic Tumor Promoter 2,2', 4,4', 5,5'-Hexabromobiphenyl.** Society of Toxicology.
- Evans, M.G., and J.E. Trosko. (1988) **Concentration/Response Effect of 2,2', 4,4', 5,5'-Hexabromobiphenyl on Cell-Cell Communication in vitro: Assessment by Fluorescence Redistribution After Photobleaching ("FRAP").** J Cell Biology and Society of Toxicology. 4(2):
- Fitz, T.A., M.M. Marr, T.L. Walden, C.A. Winkel, W.A. Schmidt and R.C. Burghardt. Uniformed Services University of the Health Sciences, Bethesda, MD. (1989) **Intercellular Communication Persists Among Transformed Granulosa (DC3) Cells Following ⁶⁰Co Irradiation.** Eicosanoids & Bioactive Lipids in Cancer & Radiation Injury. 1st International Conference, Detroit, MI.
- Gera, R., R. Kulkarni, M.H. El-Fouly, C.C. Chang, and J.E. Trosko. Michigan State University, East Lansing, MI. (1987) **Inhibition of Intercellular Communication in Human Carcinoma Cell Lines of Various Origins as Determined by Three Newly Developed Methods.** AACR.

- Hasler, C.M., M.A. Frick, M.R. Bennink and J.E. Trosko. Michigan State University, East Lansing, MI. (1988). **TPA Inhibits Gap Junctional InterCellular Communication without Generating Free Radicals.** Toxicological Implications of Altered Gap Junctional InterCellular Communication.
- Madhukar, B.V., J.E. Trosko, C.C. Chang, M.H. El-Fouty and S.Y. Oh. Michigan State University, East Lansing, MI. (1988) **Inhibition of Gap Junctional Communication in Normal Human Keratinocytes by Various Growth Factors and Tumor Promoters.** Toxicological Implications of Altered Gap Junctional InterCellular Communication
- Madhukar, B.V., C.C. Chang, C. Hasler, E. deFeyter, S.Y. Oh, B. Lockwood, and J.E. Trosko. Michigan State University, E. Lansing, MI. (1988). **Development of an In Vitro Human Epithelial System to Study Gap Junctional InterCellular Communication.** Toxicological Implications of Altered Gap Junctional InterCellular Communication.
- Madhukar, B.V., S.Y. Oh, A.H. Jiwa* and M.H. Wade*. Michigan State University and *Meridian Instruments, Inc. (1989) **Inhibition of InterCellular Communication and Increase in Intracellular Calcium In Vitro by Xenobiotic Chemicals in a Liver Epithelial Cell Line.** In Vitro Toxicology.
- Madhukar¹, B.V., S.Y. Oh¹, J. Quensen¹, M.H. Wade², and J.E. Trosko¹. ¹Michigan State University, E. Lansing, MI and ²Meridian Instruments, Inc., Okemos, MI. (1988) **Inhibition of InterCellular Communication by Heptachlor and Heptachlor Epoxide In a Rat Liver Epithelial Cell Line.** 39th Annual Meeting of the Tissue Culture Association.
- Marr, M.M., T.L. Walden, W.A. Schmidt, C.A. Winkel and T.A. Fitz. (1987) **Gap Junctions In Transformed Rat Granulosa Cells.** Medical/Materials Application of Free Electron Lasers Second Annual Contractors' Meeting, Salt Lake City, UT.
- McKarns, S.C., D.W. Bombick, and D.J. Doolittle. R.J. Reynolds Tobacco Co., Winston-Salem, NC. (1988) **Kinetic Analysis of InterCellular Communication Inhibition by TPA.** Toxicological Implications of Altered Gap Junctional InterCellular Communication.
- Meda¹, P., R.E. Hensley², and J.J. Fischer². ¹University of Geneva, Switzerland and ²Michigan State University, E. Lansing, MI. (1988) **n-Heptanol Inhibition of Cell-Cell Communication Studied by Quantitative Fluorescence Recovery.** Toxicological Implications of Altered Gap Junctional InterCellular Communication.
- Meiners¹, S., M.H. Wade², and M. Schindler¹. ¹Michigan State University, East Lansing, MI, ²Meridian Instruments, Inc., Okemos, MI. **Measurement of InterCellular Communication in Plants Using Fluorescence Photobleaching.** Journal of Cell Biology. #8:2
- Oh¹, S.Y., B.V. Madhukar¹, C.C. Chang¹, J.E. Trosko¹, and E.C. Beyer². ¹Michigan State University, E. Lansing, MI and ²Harvard Medical School, Boston, MA. (1988) **Characterization of Some Mutants Lacking Gap Junction Communication.** Toxicological Implications of Altered Gap Junctional InterCellular Communication.

- Trosko, J.E. and C.C. Chang. Michigan State University, East Lansing, MI. (1986) **Role of Intercellular Communication in Modifying the Consequences of Mutations in Somatic Cells.** Antimutagenesis and Anticarcinogenesis Mechanisms.
- Trosko, J.E. and M.G. Evans. Michigan State University, East Lansing, MI. (1987). **Application of the "FRAP" Assay to Measure Inhibition of Intercellular Communication by 2,2', 4,4', 5,5'-Hexabromobiphenyl in WB cells.** Society of Toxicology.
- Trosko, J.E., C.C. Chang, M.H. El-Fouly, R. Kulkarni and R. Gera. Michigan State University, East Lansing, MI. (1986). **Role of Gap Junctional Intercellular Communication in the Regulation of Growth and Differentiation of Normal and Malignant Human Epithelial Cells.**
- Trosko, J.E., M.G. Evans and M.H. El-Fouly. Michigan State University, East Lansing, MI. (1987) **Fluorescence Quantitation of the Scrape-Loading/Dye Transfer Assay to Measure Inhibited Cell-Cell Communication in WB Cells by PBB and Dieldrin.** Society of Toxicology.
- Wade¹, M.H. and M. Schindler². ¹Meridian Instruments, Inc., Okemos, MI and ²Michigan State University, E. Lansing, MI. (1988) **Gap Junctional Mediated Intercellular Communication is a Function of Contacting Cell Number.** FASEB.

CELL SORTING

Meridian Instruments, Inc. Application Notes

Wade, M.H. (1987) **"COOKIE CUTTER" Method of Cell Selection.**

Jwa, A.H. and J.M. Wilson, Meridian Instruments, Inc., Okemos, MI., Howard Hughes Medical Institute, University of Michigan, Ann Arbor, MI. (1990) **Sorting and Selection of Mammalian Cells Expressing the E. Coli LacZ Gene Product Using the ACAS Interactive Laser Cytometer.** Society for Analytical Cytology.

CELLULAR GROWTH

Abstracts

- Jensen, B.D., P.K. Horan and G. Poste. Smith Kline and French Laboratories, Philadelphia, PA. (1988) **Fluorescent Analysis of the Cellular Growth Rate in Adherent Cell Systems.** Society for Analytical Cytology.
- Ronot, X. Institut Biomedical des Cordeliers, 75006, Paris. (1989) **Direct Quantification of In Vitro Mitochondrial Activity of Chondrocytes Through Anchorage Dependent Laser Cytometry.** Anal. Cell Path. 1(5/6):283.
- Ronot, X. Institut Biomedical des Cordeliers, 75006, Paris. **Fluorescence Imaging of Mitochondrial Compartment in Chondrocytes Cultured on Microporous Membranes.**

DNA

- Liu, P.K., B. Goudreau, and G.S. Hsu. (1989) **A Quantitative Assay of DNA Polymerase- α In Situ at Single cells using Fluorescence Pseudo-Color Image and ACAS 470.** Environmental and Molecular Mutagenesis (Alan R. Liss) **14**, Supp. 15
- Mori, T., A.A. Wani, S.M. D'Ambrosio, C.C. Chang, and J.E. Trosko. (1989). **In situ Pyrimidine Dimer Determination by Laser Cytometry.** Photochem. Photobiol. **49**(4): 523-526

FRAP

- Balint, E., A. Aszalos, P.M. Grimley, F. Edward Herbert Medical School USUHS; Bethesda, Md 20814. (1988) **Alpha Interferon Accelerates Lateral Diffusion of Daudi Cell Surface Differentiation Antigens: Measurement by Fluorescence Redistribution After Photobleaching** Biochemical and Biophysical Research Communications, **157**(2):808-815.
- Balint, E., A. Aszalos, P.S. Pine, P.M. Grimley. (1988) **Dynamic Analyses of Lymphoblast Membranes Exposed to Alpha Interferon Using Flow Cytometry and Fluorescence Recovery After Photobleaching.** Scanning Microscopy, **2**(4) 2153-2163.
- Liu, Z.-Y., M.E. Sanders, and V.W. Hu. (1989) **Effect of Complement on the Lateral Mobility of Erythrocyte Membrane Proteins. Evidence for Terminal Complex Interaction with Cytoskeletal Components.** J. Immunol. **142**(7):2370-2376.

Meridian Instruments, Inc. Application Notes

- Schindler, M. (1987) **A Rapid Method to Analyze Diffusion of Membrane Components and Cytoplasmic Structures.**

Abstracts

- Balint, E., A. Aszalos and P.M. Grimley, F. Edward Hebert School of Medicine, USUHS, Bethesda, MD, and FDA, Washington, D.C. **Increased Lateral Mobility of Leu and IGM Surface Antigens in B Lymphoblastoid Cells (Daudi) After Binding of Alpha Interferon.** Membrane Signal Transduction.
- Bombick, D., G. Zhang, and J.E. Trosko. Michigan State University, E. Lansing, MI. (1988) **Cell Type Response to Different Chemical Modulators of Gap Junctional Function.** Society of Toxicology.
- Schindler, M., Gharyal, P., Jiang, L., Michigan State University, E. Lansing, MI. (1989) **The Dynamic Parameter Fluorescence Photobleaching as a Tool to Dissect Space in Biological Systems.** Fluorescence in Biochemistry and Cell Biology.

IMMUNOLOGY

- Falk, L.A., L.M. Wahl, and S.N. Vogel. (1988) **Analysis of Ia Antigen Expression in Macrophages Derived from Bone Marrow Cells Cultured in Granulocyte-Macrophage Colony Stimulating Factor (GM-CSF) or Macrophage Colony Stimulating Factor (CSF-1).** Journal Immunology. April 15, 140:2652-2660
- Gilbert, C.W., M.H. Zaroukian, and W.J. Esselman. (1988) **Poly-N-Acetyllactosamine Structures on Murine Cell Surface T200 Glycoprotein Participate in NK Cell Binding to yac-1 Targets.** Journal Immunology. April 15, 140:2821-2828
- Hogan, M. Michele, Pin Yu Perera and Stefanie N. Vogel. (1989) **Examination of macrophage cell surface antigen regulation by rIFN-a/B utilizing digital imaging by a novel laser detection system.** Journal Immunology. (in press)
- Melnicoff, M.J., P.S. Morahan, B.D. Jensen, E.W. Breslin, and P.K. Horan. (1987). **In vivo Labeling of Resident Peritoneal Macrophages.** Journal Leukocyte Biology 43:387-397.

Abstracts

- Frey, C.L., B.D. Jensen and D.J. Drutz. Smith Kline and French Laboratories, Philadelphia, PA. (1988) **Localization and Distribution of iC3b Binding Sites on *Candida albicans*.**
- Simonsen¹, M.G., K. Crickard¹, M.H. Wade², and R.J. Bernacki^{3,1} ¹Baxter Health Care Corp., Irvine, CA, ²Meridian Instruments, Inc., Okemos, MI, ³Roswell Park Memorial Institute, Buffalo, NY. (1987) **Identification of Human Tumor Cells in Mixed Cultures with Tumor Specific Monoclonal Markers.**
- Takahama, Y., S. Ono, K. Ishihara and T. Hamaoka, Division of Oncogenesis, Biomedical Research Center, Osaka University Medical School, Osaka, Japan. (1989) **Cluster Formation Among Small Residing B Lymphocytes Leading to B Cell Activation.** Immunology 1(1):36-42

In Situ HYBRIDIZATION

- Bresser, J. and M.J. Evinger-Hodges (1987) **Comparison and Optimization of In Situ Hybridization Procedures Yielding Rapid, Sensitive mRNA Detections.** Gene Anal Technology 4:89-104

INSTRUMENT DESCRIPTION

- Schindler, M., M.L. Allen, M.R. Olinger and J.F. Holland. (1985) **Automated Analysis and Survival Selection of Anchorage-Dependent Cells Under Normal Growth Conditions.** Cytometry 6:368-374
- Schindler, M., L.-W. Jiang, M. Swaisgood and M.H. Wade. (1989) **Analysis, Selection, and Sorting of Anchorage Dependent Cells Under Growth Conditions.** Methods in Cell Biology 32

Yentsch, C.M. and S.A. Pomponi. (1986) **Automated Individual Cell Analysis in Aquatic Research.** International Review of Cytology. 105:183-243

Abstracts

Holland¹, J.F., J.E. Trosko¹, and M.H. Wade^{2, 1} Michigan State University, East Lansing, MI, ²Meridian Instruments, Inc., Okemos, MI. **An Instrument System for Fluorescence Measurements within the Living Cell: A Possible Alternative to Animal Studies for Tumorigenic Screening.**

KINETIC

Abstracts

Frey, C.L., B.D. Jensen and D.J. Drutz. Smith Kline and French Laboratories, Philadelphia, PA. (1988) **Time Resolved Analysis of Expression of iC3b Binding Sites of *Candida albicans* by Fluorescence Mapping.** Society of Analytical Cytology.

Grillone, L.R., B.D. Jensen and S.T. Crooke. Smith Kline and French Laboratories, Philadelphia, PA. **Vasopressin-Induced Changes in F-actin Architecture: Image Analysis of Fluorescently Labeled Filaments.** Cytometry. (1) 78

Grillone, L.R., B.D. Jensen, and S.T. Crooke. Smith Kline and French Laboratories, Philadelphia, PA. **Imaging of Vasopressin-Induced f-Actin Changes in Smooth Muscle Cells.** Journal of Cell Biology. 103(5):265

LASER SURGERY ON CELLS

Joshi, H.C., D. Chu, R.E. Buxbaum and S.R. Heidemann. (1985) **Tension and Compression in the Cytoskeleton of PC 12 neurites.** Journal of Cell Biology. 101:697-705

Abstracts

Grillone, L.R., B.D. Jensen and S.T. Crooke. Smith Kline and French Laboratories, Philadelphia, PA. (1988) **Regulation of pH in Single Smooth Muscle Cells after Vasopressin Stimulation.** Journal for Analytical Cytology

Grillone, L.R., B.D. Jensen and S.T. Crooke. Smith Kline and French Laboratories, Philadelphia, PA. (1988) **Vasopressin Induces a Cytosolic pH Response in a Smooth Muscle Cell Line, A-10, that Occurs Subsequent to the Calcium Spike.** Fourth Int'l. Congress of cell Biology.

Stevens R.J., Kenton Sanders. University of Nevada, Reno, NV. **Sodium-Hydrogen Exchange in Canine Colonic Smooth Muscle Cells.**

PHORBOL ESTERS

Abstracts

Bernstein, R.L., H.A. Pershadsingh. University of California, San Francisco, CA. (1988) **Fluorescent Phorbol Esters as Probes of Brain Protein Kinase C Activity and Translocation in Neural Cells.** ASCB.

Bernstein, R.L., M. Balazs¹, J. Szollosi¹, S. Damjanovich², R.P. Haugland³, A.P. Guzikowski³, S.A. McQuiston⁴, M.J. Fulwyler, and H.A. Pershadsingh. University of California, San Francisco, CA, ¹University Medical School, Debrecen H-4012, Hungary, ³Molecular Probes, Inc. Eugene, OR, ⁴Meridian Instruments, Inc., Okemos, MI. (1988) **Monitoring Novel Fluorescent Phorbol Esters in NG 108-15 Cells.**

Pershadsingh¹, H.A., M. Balazs¹, J. Szollosi¹, A.P. Guzikowski², R.P. Haugland², S.A. McQuiston³, S.C. Griffin³, and M.J. Fulwyler. UCSF, San Francisco, CA, ¹Univ. Medical School, Debrecen H-4012, Hungary, ²Molecular Probes, Inc. Eugene, OR, ³Meridian Instruments, Inc., Okemos, MI. (1988) **Microscopic Imaging of Novel Fluorescent Phorbol Esters in Diverse Cell Lines.** Fourth Int'l Congress of Cell Biology.

Pershadsingh¹, H.A., G. Gold³, R.P. Haugland², S.A. McQuiston¹ and M.J. Fulwyler¹. University of California¹, San Francisco, CA, Molecular Probes, Inc.², Eugene, OR, Lilly Research Laboratories³, Indianapolis, IN. (1989) **Probing the Subcellular Distribution of Protein Kinase C in Insulin-Secreting Cells using a Novel Fluorescent Phorbol Ester.** ASCB

Pershadsingh¹, H.A., Balazs, Margit, Szollosi, Janos, Haugland, P. Richard, Guzikowski, P. Anthony, Fulwyler, J. Mack, Damjanovich, Sandor, Feuerstein, Burt, G. University of California, San Francisco, CA, Department of Laboratory Medicine. (1989) **Fluorescent Tetradecanoylphorbol Acetate As A Probe Of Protein Kinase C in the P₃HR-1 Human B Lymphocyte Tumor Line.** Submitted Journal of Cell Biology.

PHOTOACTIVATION

Schneider, M. (1989) **Dynamic Imaging-creating subcellular views by Photoactivation and Photobleaching.** Digitized Video Microscopy Alan R. Liss (in press).

SUBCELLULAR LOCALIZATION

Liu, Z.-Y., R. Solow and V.W. Hu. (1988) **Fluorescence Analysis of Size Distribution and Mode of Dye Release from Carboxyfluorescein-Loaded Vesicles: Application to the Study of Complement-Membrane Interaction.** Biochem Biophys. Acta 945:253-262.

Matsuura, T., S. Nagamori, K. Fujise, S. Hasumura, S. Homma, H. Sugino, K. Shimizu, M. Niya, H. Kameda, and K. Hirose. (1989) **Retinol Transport in Cultured Fat-storing Cells of Rat Liver: Quantitative Analysis by Anchored Cell Analysis and Sorting System.** Laboratory Investigation 61(1):107-115

Mori, T.¹, A.A. Wani¹, S.M. D'Ambrosio², C.-C. Chang¹ and J.E. Trecko¹, ¹Department of Pediatrics and Human Development, Michigan State University and ²Department of Radiology, Ohio State University (1989) **In situ Pyrimidine Dimer Determination by Laser Cytometry.** Photochemistry and Photobiology 49(4): 523-526.

Moutsatsos, I.K., M.H. Wade, M. Schindler, and J.L. Wang. (1987) **Endogenous Lectins from Cultured Cells: Nuclear Localization of Carbohydrate-Binding Protein 35 in Proliferating 3T3 Fibroblasts.** Proc. Natl. Acad. Sci. USA, 84:6452-6456.

Ronot, X. (1989) **Fluorescence Imaging of Mitochondrial Compartment in Chondrocytes Cultured on Microporous Membranes.** Presented at the First Conference of the European Society for Analytical Cellular Pathology

Shuck, J.R., M.H. Wade and C.W. Lin. (1988) **Subcellular Localization of Tetraporphyrin Derivative in Bladder Tumor Cells in Culture.** (Submitted for publication)

Yang, H.M., and R.A. Reistad. (1987) **Doxorubicin Conjugated with a Monoclonal Antibody Directed to a Human Melanoma-Associated Proteoglycan Suppresses the Growth of Established Tumor Xenografts in Nude Mice.** Proc. Natl. Acad. Sci. USA 85:1189-1193.

Moran Instruments, Inc. Application Notes

Bimock, D.W. (1988) **Detection and quantitation of Peroxides and Hydroperoxides in Attached Cells Using the Dye 2,7-Dichlorofluorescein.**

Wang, J.L. (1987) **Digital Image and Statistical Analysis of Immunofluorescence at the Level of Single Cells.**

Abstracts

DeGiandro, D. and M. Schindler, Michigan State University, East Lansing, MI. **Nuclear Glycosylation Appears Unaffected in Lec 1.3c CHO Cells.**

Ko, V. Uniformed Services University of the Health Sciences, Bethesda, MD. **Fluorescence Analysis of Size Distribution and State of Aggregation of Carboxyfluorescein-Loaded Vesicles: Application to the Study of Complement-Membrane Interactions.**

Jensen, E.D. Zynaxis Cell Science, Inc., Malvern, PA. (1988) **Measurement of Cycling among Plasma Membrane and Interacellular Membrane Compartments in Smooth Muscle Cells.** ASCB

Koren, E., L. Hunt, C. Corder and G. Wen, Oklahoma Medical Research Foundation, Oklahoma City, OK. (1988) **Cholesterol Ester Accumulation in Macrophages from Patients with Coronary Artery Disease.** 8th International Symposium on Atherosclerosis

McConathy, W.J., E. Koren, V.N. Troy and D.L. Stiers, Oklahoma Medical Research Foundation, Oklahoma City, OK. **Plasma Protein Binding and Macrophage Uptake of Cholesterol Crystals.** FASEB, 1988, ASCB/ASBMB, 1989.

Wang¹, J.L., M.H. Wade², and M. Schindler³. ¹The Rockefeller University, New York, NY, ²Meridian Instruments, Inc., Okemos, MI, ³Michigan State University, East Lansing, MI. **Fluorescence Image Analysis of Statin and Cyclin-like Protein Expression In Replicating and Nonreplicating Cells.**

TRANSPORT PROCESSES

Baron-Epel, O., D. Hernandez, L.-W. Jiang, S. Meiners, and M. Schindler. (1988) **Dynamic Continuity of Cytoplasmic and Membrane Compartments Between Plant Cells.** Journal of Cell Biology. 106:715-721.

Jiang, L.-W., and M. Schindler. (1987) **Nuclear Transport in 3T3 Fibroblasts: Effects of Growth Factors and Transformation and Cell Shape.** Journal of Cell Biology. 106:13-19.

Abstracts

Jiang, L.-W., and M. Schindler, Michigan State University, E. Lansing, MI. **Poipeptide Growth Factors Enhance Macromolecular Nuclear Transport.**

CONFOCAL

Kanistanaux, D., and P.H. Burrill. (1989) **Confocal Fluorescence Microscopy: A New Dimension.** Biomed. Prod. Sept:54,72.

Rev. 12/06/89

APPENDIX B

Application Notes on new techniques that the Upgrade of ACAS can perform:

1. Calcium measurements.
2. Free Radicals.
3. Reprint "In situ (6-4) photoproduct determination by laser cytometry and autoradiography".
4. Reprint "Characterization of an in vitro human kidney epithelial system to study gap junctional intercellular communication".

ACAS 570 INTERACTIVE LASER CYTOMETER

MERIDIAN
Instruments, Inc.

Application Note Number E-2
Calcium Analysis

DETECTING CALCIUM RESPONSE IN CULTURED CELLS USING THE VISIBLE WAVELENGTH CALCIUM PROBE, FLUO-3

Margaret H. Wade, Meridian Instruments, Inc.,
Okemos, MI.

The measurement and quantitation of intracellular calcium responses within cells is of great importance to cell biologists (see Application Note E-1). Small changes in free intracellular calcium can effect gap junction function, cell division, muscle contraction, neutrophil stimulation and other cell functions. Recently several fluorescent probes to measure calcium have been developed by Roger Tsien. Indo-1 and Fura-2 were developed as dual emission or dual excitation calcium sensitive probes, and more recently the single excitation and emission probes, Fluo-3 and Rhod-1 have been made available (1,2).

Fluo-3 is a new fluorescent probe which can be excited by the 488 nm line of the argon laser, and emits in the visible spectrum similar to fluorescein (520 nm). The advantage of this probe is that UV excitation is not required; however a major disadvantage is that it is difficult to quantitate the absolute calcium levels in cells because it is not possible to use a ratio technique. A calcium response, however, is relatively easy to detect in cells, and further quantitation could be done with the calcium ratio probe, Indo-1. Fluo-3 also offers the possibility of dual labeling experiments where one may wish to monitor another property of the cell simultaneously with calcium levels, such as pH or membrane potentials.

A solution of the Fluo-3 free acid in high potassium buffer containing 1 mM EGTA, pH 7.45 was examined using the ACAS 570 Interactive Laser Cytometer. The sample was excited with 488 nm light from the argon ion laser, and the emission above 515 nm was detected with the photomultiplier tube and digitized. The absolute fluorescent units are plotted in Figure 1 as a function of free calcium (Orion standard) calculated from total calcium using the K_d of EGTA (105 nM) in this buffer. There is clearly a relationship between added calcium and increased fluorescence. However, as the fluorescent value in cells may depend upon dye leakage, photobleaching, degree of labelling, etc., this type of calibration curve cannot be used for quantitation.

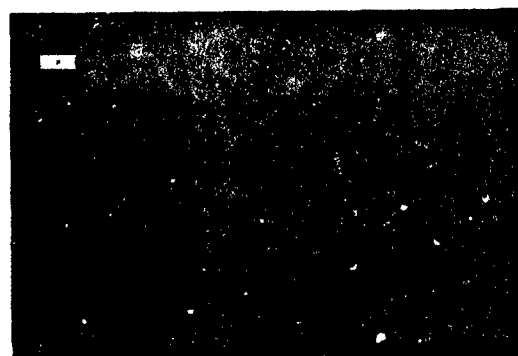


Figure 1. Standard curve showing Fluo-3 free acid fluorescence as a function of calcium concentration.

Human teratocarcinoma cells (HT) or normal diploid fibroblasts were grown in 35 mm tissue culture dishes and labelled with the cell permeable form of the dye, Fluo-3 AM, in PBS containing magnesium and calcium. Although the manufacturer suggests the addition of the solubilizing agent, pluronic acid (3), to the incubation media, we found this addition unnecessary for these cells (see Table 1). After rinsing, the dish was placed on the X-Y scanning stage and either a fluorescent image of a field of cells or a line scan across one cell was generated as a function of time. In some experiments, ionomycin was added and the calcium response was measured.

The concentration of Fluo-3 necessary to label cells was examined (Table 1). Cells incubated with 1 μ M Fluo-3 were very poorly labeled and not included in the figure; those incubated with 3 or 5 μ M for 1 or 2 hours showed good labelling without .05% pluronic acid, and less labelling with pluronic acid. We chose 5 μ M Fluo-3 for 1 hour with no pluronic acid to perform the remaining experiments.

Average Fluorescence of Fluo-3

	<u>3 μM</u>	<u>5 μM</u>
1 hr	285	1190
1 hr +Pluronic	212	598
2 hr	254	821
2 hr +Pluronic	127	703

Table 1. Average fluorescence per cell after various labelling conditions.

The photobleaching characteristics of cells labelled with Fluo-3 were determined by repetitively scanning the same area of the cells and plotting the integrated value as a function of scan number. Less than .1% fluorescence per scan is lost when using the ACAS (data not shown). This is a result of the use of a focussed laser beam (~1 micron in diameter) and Acousto-Optic Modulator (AOM). The AOM pulses the laser such that the sample is only illuminated for a brief period of time over a discrete area.

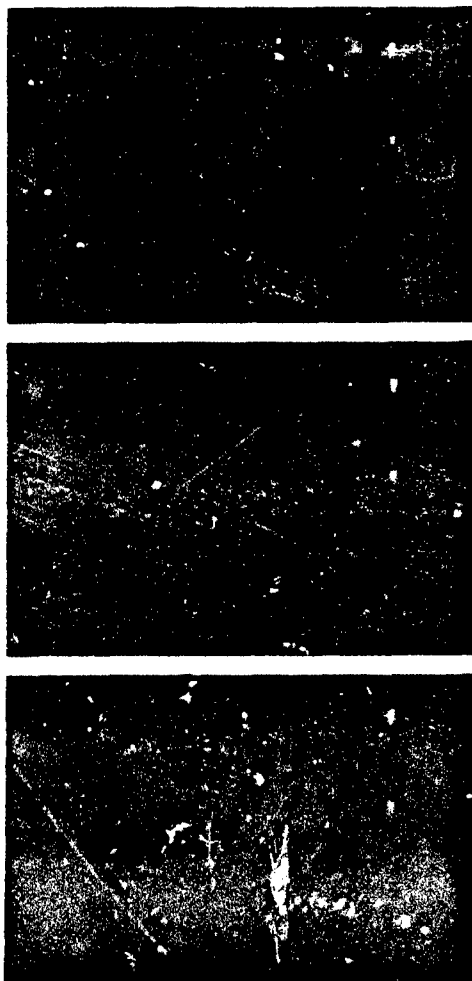


Figure 2. Two dimensional images of Fluo-3 incorporation in normal diploid fibroblasts A) before ionomycin B) immediately after ionomycin and C) 4 minutes after ionomycin.

The response of cells to the ionophore, ionomycin, is shown in Figures 2 and 3. HT cells were labelled with Fluo-3, and after the second image scan, 3 μ M ionomycin in PBS containing calcium and magnesium, was added while the scanning continued. Figure 2a is a two dimensional pseudo-color image of normal diploid fibroblasts before ionomycin is added, Figure 2b is the same field immediately after ionomycin, and Figure 2c is the same field about 4 minutes later when the cells

have returned to baseline values. Figure 3 is the time plot of a similar experiment with HT cells showing the time response to ionomycin. It is clear that an increase in fluorescence was detected, with a gradual return to baseline values. No attempt was made to correlate these fluorescent values with absolute calcium values.

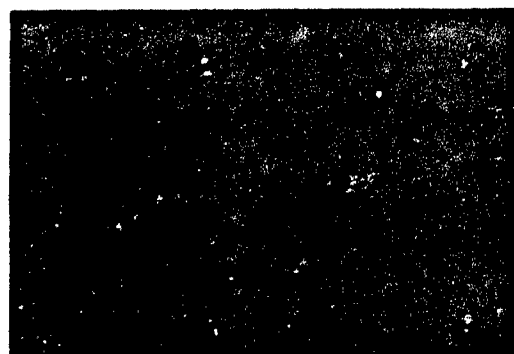


Figure 3. Time plot of ionomycin response in HT cells.

The new calcium sensitive probe, Fluo-3, will be very useful for determination of qualitative calcium changes in cells as the result of agonist treatment, such as growth factors, antibodies, etc. Although various methods have been suggested for quantitating calcium based on a fluorescent signal (3), one will always be concerned about the effect of dye leakage, photobleaching of the signal, unequal labelling among cells, etc. Fluo-3 is particularly useful as a first approximation of a calcium response because it uses visible, fluorescein-like, excitation and emission, and is suitable for simultaneous measurements of calcium and another fluorescent probe, i.e. pH.

References:

1. Minta, A., J. P. Y. Kao and R.Y. Tsien (1989) *J. of Biol. Chem.* 264:8171-8178.
2. Kao, J. P. Y., A. T. Harootunian and R.Y. Tsien (1989) *J. Biol. Chem.* 264:8179-8184.
3. Haugland, R. (1989) *Molecular Probes Catalog*

MERIDIANTM
Instruments, Inc. 2310 Science Parkway
Okemos, Michigan 48864
517/349-7200

Detection and Quantitation of Peroxides and Hydroperoxides in Attached Cells Using the Dye 2,7-Dichlorofluorescein

David W. Bombick, R. J. Reynolds Tobacco Company, Winston-Salem, NC 27102

INTRODUCTION

Several active oxygen species are important in affecting biological systems. Peroxides and hydroperoxides (e.g. hydrogen peroxide) may form free radicals (such as the hydroxyl radical) which can react with many bio-organic compounds including proteins, nucleic acids and lipids. Since some of these reactions are harmful, measurement of active oxygen species can be important in toxicological research in areas including genetic toxicology, *in vitro* toxicology, immunotoxicology, reproductive toxicology and neurotoxicology.

Methods to quantitate hydroperoxides and peroxides in living cells are limited at the present time. Most techniques require cell-free systems and involve chromophores generated by reactions of thiobarbituric acid with malondialdehyde (1,2) or I^- with the peroxide (3). An assay to detect hydroperoxides based on the conversion of 2,7-dichlorofluorescein to the fluorescent 2,7-dichlorofluorescein by hydroperoxides has recently been developed (4,5,6). This reaction was sensitive to picomole levels and linear with respect to varying concentrations of a number of hydroperoxides and peroxides. The chemical 2,7-dichlorofluorescein diacetate, which freely crosses cellular membranes was used to facilitate loading into cells. Cellular esterases cleave the acetate groups and the ionic 2,7-dichlorofluorescein, which cannot cross cellular membranes, is trapped in cells. Once the 2,7-dichlorofluorescein is trapped within a cell it is available for oxidation to 2,7-dichlorofluorescein by the peroxides or hydroperoxides present. This method has been adapted for use with the Meridian Instruments ACAS 470 Interactive Laser Cytometer™ to allow quantitation of peroxides and hydroperoxides on a single cell basis.

METHODS

The dye 2,7-dichlorofluorescein diacetate

(Molecular Probes, Eugene, OR) was dissolved in ethanol at a concentration of 1 mM. The dye was loaded into cells (5 ul/ml) in culture media for 30 minutes at 37° C. Chemical treatment of the cells may precede, follow, or occur simultaneously with the loading of the dye as long as this is kept constant in the experimental design and the treatment does not inhibit esterase activity. After 30 minutes of dye exposure, the plate of cells is rinsed with PBS containing 5 mM glucose (PBSg) approximately five times. The cells are maintained in a small amount of the PBSg (1-2mls) and examined with the ACAS 470 Interactive Laser Cytometer. The instrument is equipped with an argon ion laser tuned to 488 nm, and emission above 515 nm is quantitated from two dimensional image scans generated by a 1 micron laser beam and an X-Y scanning stage. An average fluorescence value from a confluent cell area or the fluorescence from single cells can be obtained. A standard curve (Figure 1) can be generated by adding varying concentrations of hydrogen peroxide (e.g., 10 ng/ml to 1 mg/ml) in PBSg with or without a confluent cell monolayer, depending upon the experimental design. This standard curve can be used to express fluorescence values as a function of hydrogen peroxide quantity.

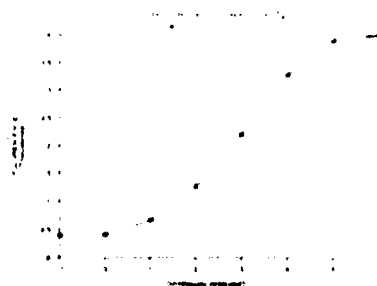


Figure 1. Standard curve showing hydrogen peroxide concentration (log ug/ml) as a function of fluorescence intensity.

RESULTS AND DISCUSSION

Quantitative analyses of peroxides and hydroperoxides in several cell types have been done using 2,7-dichlorofluorescein diacetate and the ACAS 470. Figure 2 illustrates the respiratory burst of neutrophils after addition of 10 nm/ml of PMA (phorbol-12-myristate-13-acetate). The cells were plated in multiple

dishes at equal densities and the average fluorescence values for neutrophils was calculated. At 0, 1, 5 and 10 minutes after PMA treatment those values were 305, 430, 931 and 1551, respectively. These values can be related to hydrogen peroxide concentrations with a standard curve. The rapid increase in peroxide formation parallels results obtained with flow cytometry (5). The grey areas within some of the images indicate the presence of cells which have not responded and are therefore not fluorescent. The images from the ACAS 470 demonstrate that the response is heterogeneous among the neutrophil population and that not all cells are activated to the same degree. This also demonstrates the importance of analyzing single cells, rather than a whole population. Figure 3 illustrates the fluorescent images generated when WB rat liver oval cells are exposed to a 12.45 mM concentration of hydrogen peroxide for 5 minutes. The conversion of 2,7-dichlorofluorescein to the fluorescent 2,7-dichlorofluorescein in these cells is rapidly increased when hydrogen peroxide is exogenously added. Again, the response in the WB cells is heterogeneous which may indicate different subpopulations with varying protective mechanisms against hydrogen peroxide.

This method appears useful in determining whether particular toxicants generate peroxide or hydroperoxide species. Formation of these reactive oxygen species may be correlated to cytotoxicity parameters or other toxic manifestations. A further use of this method may include the examination of lipid peroxidation and subsequent membrane damage. Finally, an *in vitro* assay could be developed to test a chemical's anti-oxidant capability.

In summary this report discusses a method to quantitate the formation of reactive oxygen species *in vitro*. The advantages include the capability of carrying out the analyses in attached monolayers of cells rather than cell suspensions. In addition, collection of these peroxides and hydroperoxides can be done quickly and on a single cell basis.

REFERENCES

1. Dahle, I.K., Hill, E.G. and Holman, R.T. (1962) *Anal. Biochem. Biophys.* 98, 253.
2. Pryor, W.A., Stanley, J.P. and Blair, E. (1976) *Lipids* 11, 370.
3. Hicks, M., and Babicki, J.M. (1979) *Anal. Biochem.* 99, 249.

4. Cathcart, R., Schwiart, E., and Ames, B.N. (1983) *Anal. Biochem.* 134, 111.
5. Bass, D.A., Parce, J.W., Dechatelet, L.R., Szejda, P., Seeds, M.C. and Thomas, M. (1983) *J. Immunol.* 130, 1910.
6. Szedo, P., Parce, J.W., Sees, M.S. and Bass, D.A. (1984) *J. Immunol.* 133, 3303.

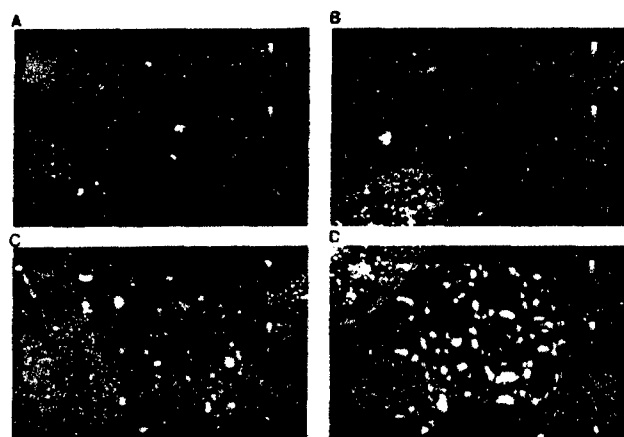


Figure 2. Two dimensional pseudo-color images depicting fluorescence in cells treated with 10 ng/ml TPA for 0, 1, 5 and 10 minutes (panels A, B, C and D respectively).

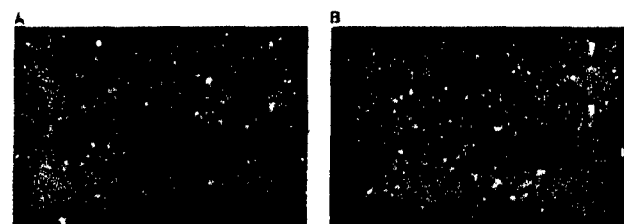


Figure 3. Fluorescent images of rat liver WB cells exposed to no (left panel) or 12.45 mM hydrogen peroxide (right panel).

MERIDIAN™
Instruments, Inc. 2310 Science Parkway
Okemos, Michigan 48864
517/349-7200

In situ (6-4) photoproduct determination by laser cytometry
and autoradiography

Toshio Mori^{1,3}, Tsukasa Matsunaga¹, Chia-Cheng Chang²,
James E. Trosko², and Osamu Nikaido¹

¹Division of Radiation Biology, Faculty of Pharmaceutical
Sciences, Kanazawa University, Kanazawa 920, Japan

²Department of Pediatrics and Human Development, Michigan State
University, East Lansing, MI 48824, U.S.A.

Running title: In situ (6-4) photoproduct determination

Key words: Monoclonal antibody, (6-4)photoproduct, UV damage,
Laser cytometry, In situ determination

³To whom all correspondence should be sent:

Dr. Toshio Mori
Division of Radiation Biology
Faculty of Pharmaceutical Sciences
Kanazawa University
Kanazawa 920, Japan
Telefax. 81-762-64-1088

Summary

The UV-induced (6-4) photoproducts and repair in individual human cells were quantitatively determined by using argon-laser imaging microspectrofluorometry or autoradiography with a well-characterized monoclonal antibody against (6-4) photoproducts. (6-4) Photoproduct induction curves were linear as a function of UV dose, using both methods. The formation of (6-4) photoproducts was detected in the cells irradiated with as low as 10 and 25 J/m² of UV by autoradiography and laser cytometry, respectively. Normal cells repaired more than 80% of the initial damage within 4 hr post-irradiation. In contrast, almost no repair was observed in xeroderma pigmentosum cells (complementation group A) within 8 hr.

Introduction

Available evidence has not ruled out (6-4) photoproduct as a premutagenic lesion (Franklin and Haseltine, 1986). On the contrary, there is evidence that both cyclobutane pyrimidine dimers and (6-4) photoproducts are responsible for UV-induced mutations in *E. coli* and in human cells (Brash, 1988). Additional evidence indicating a major biological role for (6-4) photoproducts came from studies of a chemical mutagen-induced revertant of xeroderma pigmentosum cells (x.p., complementation group A) which is normal in repair of (6-4) photoproducts but remains defective in repair of pyrimidine dimers. The revertant was found to be normal in repair replication and UV-induced sister chromatid exchanges and mutations (Cleaver et al 1987, 1988).

Although there is increased interest in (6-4) photoproducts, a sensitive method for detecting (6-4) photoproducts in DNA induced by UV at physiological dose is not readily available except the immunological method (Mitchell and Clarkson, 1984; Eggset et al., 1987; Mori et al., 1988). Besides the advantage of sensitivity, the immunological method also allows researchers to study the distribution of (6-4) photoproducts in individual cells within a tissue or in a population of cultured cells in vitro. For this purpose the antibody used should be very specific for (6-4) photoproducts. The monoclonal antibody specific for (6-4) photoproducts appears to be the best choice for this kind of study.

In this report, we developed a sensitive method to measure the (6-4) photoproducts and their repair in individual human cells in culture by using an argon-laser imaging microspectrofluorometry (Meridian ACAS 470) or autoradiography with a monoclonal antibody against (6-4) photoproducts.

Materials and Methods

Cells and media

Normal human skin fibroblasts (AH) and xeroderma pigmentosum (XP) cells, XP1EH (complementation group A), were cultured in Dulbecco's modified Eagle's medium (DMEM, Nissui Seiyaku, Tokyo), supplemented with 10% fetal bovine serum (FBS, M.A. Bioproducts, Walkersville, MD). Another XP cell strain, XP12BE (CRL1223, complementation group A), was obtained from the American Type Culture Collection (Rockville, MD). XP12BE cells were grown in a

modified Eagle's medium (Eagle, 1959) (Earle's balanced salt solution with a 50% increase in all vitamins and essential amino acids except glutamine) and supplemented with nonessential amino acids (100% increase), 1 mM sodium pyruvate, and 10% FBS (Hazleton Research Products, Denver PA). Cells were grown in humidified air with 5% CO₂ at 37°C.

Preparation of a monoclonal antibody specific for (6-4) photoproduct

The details of the establishment of a monoclonal antibody (64M-1) specific for (6-4) photoproducts have been published previously (Mori *et al.*, 1988). Briefly, UV-irradiated single-stranded DNA (20kJ/m², UVssDNA) conjugated with methylated bovine serum albumin (mBSA) was injected 3 times into BALB/c mice. Three days after a booster injection of UV-poly(dT)-mBSA, spleen cells were fused with mouse myeloma cells (P3-NSI/1-Ag4-1). Hybridomas secreting antibody recognizing UV-irradiated DNA were screened and the cells in the promising cells were cloned 3 times by limiting dilution. Subsequently, 64M-1, a monoclonal antibody recognizing TT(6-4) photoproducts and TC(6-4) photoproducts was established (Matsunaga *et al.*, In press).

In situ (6-4) photoproduct determination by laser cytometry.

XP12BE cells (2×10^5) were directly plated in 35 mm plastic dishes and cultured overnight. After washing with phosphate buffered saline (PBS), cells were irradiated with UV and then fixed with cold methanol : acetic acid (3:1) for 20 min on ice, washed by ethanol and dried. Cells were then treated with 0.07 M NaOH/70%

ethanol for 2 min to denature DNA, followed by washing with PBS (5 times). Cells were incubated with 1% BSA in PBS for 30 min at room temperature to prevent non-specific antibody binding. After washing 5 times with PBS, cells were then incubated with the monoclonal antibody (64M-1) for 30 min at room temperature and washed again with PBS. Cells were then incubated with sheep anti-mouse IgG conjugated with biotin, F(ab')₂ fragment (Sigma, St. Louis, MO, final concentration: 25 µg/ml) for 30 min at room temperature and washed five times with PBS. The last step involves the incubation of cells with streptavidin conjugated with FITC (Sigma, St. Louis, MO, final concentration: 10 µg/ml) for 15 min at room temperature, and washing five times with PBS and twice with distilled water. All antibodies and fluorescent chemicals were diluted with PBS containing 0.1% BSA and 0.1% NaN₃. After drying, the cells were covered by a drop of PBS/glycerol (containing 5% wt/vol n-propylgallate), and a clean coverglass with the edges sealed with nail polish. Samples were observed using the Meridian ACAS 470 workstation (Wade *et al.*, 1986, Mori *et al.*, 1989, Meridian Instruments, Okemos, MI). For each treatment, the average fluorescent intensity of about 200 cells in 10 randomly selected areas is presented.

Preparation of ³H-monoclonal antibody

We followed the method described by Galfre and Milstein (1981). Briefly, 1.2×10^7 hybridoma cells secreting 64M-1 were washed twice with lysine-free DMEM (Nissui Seiyaku, Tokyo) and incubated for 20 hr at 37°C in the incorporation medium: lysine-free DMEM,

3.0 ml; L-[4,5-³H] lysine monohydrochloride (37 MBq/ml, 2.7 TBq/mmol, Amersham, England), 2.16 ml; 10 times Hank's balanced salt solution, 0.24 ml; dialyzed FBS, 0.6 ml. Supernatant was collected and applied on a column containing protein A-Sepharose CL-4B (Pharmacia, Sweden). IgG_{2b} fractions were collected and dialyzed overnight against PBS.

In situ (6-4) photoproduct determination by autoradiography

The method has been described previously (Mori *et al.*, 1988). Briefly, normal skin fibroblasts or XP1EH cells (2×10^5) were incubated for 24 hr in a Lab-Tek chamber (No. 4802, Naperville, IL) and washed twice with PBS. Cells were irradiated with various doses of UV at room temperature at a dose rate of 1.13 J/m²/s. Immediately after irradiation or after post-irradiation incubation, cells were fixed with ice-cold ethanol:acetic acid (3:1) for 20 min at 4°C. Slides with cells were passed through 0.07M NaOH in 70% ethanol for 2 min, rinsed with 70% ethanol and dried. Cells were incubated with PBS containing 10% calf serum (CS) at 37°C for 30 min and subsequently with the ³H-monoclonal antibody (13.9 µg/ml, 3.8×10^5 dpm/µg protein) for 30 min. Cells were then serially washed with PBS containing 10% CS, PBS and water. Dried slides were processed for autoradiography using NR-M2 emulsion (Konishiroku Photo Inc., Tokyo). The slides were kept in a freezer at -20°C for exposure for 30 days, developed, fixed and stained with Giemsa solution.

Results

In situ (6-4) photoproduct determination by laser cytometry

Fig. 1 shows the typical fluorescent images of XP cells irradiated with various doses of UV obtained from the ACAS 470 workstation. In the irradiated cells, there was a weak fluorescence on the whole cells. In the cells irradiated with UV, there was a strong fluorescence over the nuclei, in addition to a weak fluorescence over the whole cells. The fluorescent intensities on the nuclei clearly increased with increasing UV doses.

The software associated with the Meridian ACAS 470 workstation allows one to eliminate the background level fluorescence and quantitates the fluorescent intensity per nucleus from the original fluorescent image shown in Fig. 1. Using this method, we measured the fluorescent intensity per cell as a function of UV dose in XP cells (Fig. 2). There was a little background level fluorescence in unirradiated cells and the fluorescent intensity per cell increased linearly with increasing UV doses. We could detect the formation of (6-4) photoproducts in individual cells irradiated with a UV dose as low as 25 J/m^2 . The method, however is not sensitive enough to detect (6-4) photoproducts in individual cells induced at lower but biologically significant dose level. To increase the sensitivity, we tried the autoradiography with the ^3H -labelled monoclonal antibody as described in the following.

In situ (6-4) photoproduct determination by autoradiography

Fig. 3 shows the autoradiographs of UV-irradiated and

unirradiated normal human cells treated with ^3H -monoclonal antibody. In the unirradiated cells, some background level grains were observed on the whole cells. In the cells irradiated with 50 J/m^2 , many grains appeared only on the nuclei.

To quantitate the formation of (6-4) photoproducts in individual cells irradiated with various UV doses, we counted the grain number per nucleus under a microscope. Fig. 4 shows the induction of (6-4) photoproducts in the DNA of UV-irradiated normal and XP cells. The number of grains per nucleus increased linearly with increasing UV doses. As expected, there was no difference in the induction of DNA damage between the two types of cells. Even at as low as 10 J/m^2 of UV dose, there was a significant increase of average grain number per nucleus compared with control cells.

Fig. 5 and 6 show the time course of repair in normal and XP cells after irradiation with 20 and 40 J/m^2 of UV. Normal cells repaired more than 80% of the initial (6-4) photoproducts within 4 hr. On the other hand, XP cells were incapable of repairing (6-4) photoproduct within 8 hr.

Discussion

The results of this study indicate that the two methods, laser cytometry or autoradiography combined with a monoclonal antibody against (6-4) photoproducts, are suitable for measuring (6-4) photoproducts induced by physiological UV doses in individual human cells. The linear (6-4) photoproducts induction curves as a function of UV dose were obtained using both methods. The

formation of (6-4) photoproducts in the cells irradiated with as low as 10 and 25 J/m² of UV were detected by autoradiography and laser cytometry, respectively. One third of thymine dimer induction (Mitchell, 1988), 3.2×10^{-19} mole (1.9×10^5 /cell) and 7.9×10^{-19} mole (4.8×10^5 /cell) of (6-4) photoproducts could be detected by autoradiography and laser cytometry respectively, assuming 10 pg DNA/cell, and 0.0025% of total thymines (unpublished data) and (6-4) photoproduct induction.

Although XP cells could repair almost no (6-4) photoproducts within 8 hr, normal cells repaired more than 80% of the initial damage within 4 hr. These results were quite different from the repair patterns of pyrimidine dimers detected by laser cytometry (Mori *et al.*, 1989), in that normal cells repaired 50% of the initial pyrimidine dimers within 8 hr and 60% at 24 hr after UV irradiation. The fast repair pattern of (6-4) photoproducts in normal cells resembles the results obtained by Mitchell *et al.* (1985). In their study, the removal of 75% of (6-4) photoproducts was observed within 4 hr in human cells. These results, in combination with the results obtained by Cleaver *et al.* (1987, 1988), suggest that there might be different short repair processes for (6-4) photoproducts and pyrimidine dimers, in addition to the long common process.

We believe this is the first report of *in situ* (6-4) photoproduct determination with a well-characterized monoclonal antibody. The monoclonal antibody we used is very specific for (6-4) photoproducts formed in TT or TC sequence (Matsunaga *et al.*,

In press). Polyclonal antibodies against UV-irradiated DNA do not seem to be suitable for in situ determination of a specific damage, because they might contain mixed populations of antibodies recognizing various kinds of damages; i.e. (6-4) photoproducts, pyrimidine dimers and thymine glycol (Mitchell and Clarkson; 1984, Eggset et al.; 1987).

Recently, we have established four new monoclonal antibodies specific for (6-4) photoproducts (unpublished data). Preliminary data show that three of them seem to be more sensitive for detection of (6-4) photoproducts than the 64M-1 antibody. We should be able to increase sensitivity for in situ (6-4) photoproducts determination using these new antibodies.

Acknowledgements

We wish to thank Mrs. Beth Lockwood for her technical assistance. This work was carried out under the NIBB Cooperative Research Program for the Okazaki Large Spectrograph (88-505). This work was supported by a grant from the NCI (CA21104) and NIH (1S10RR03842-01) to J.E.T. The work was also supported by a Grant-in Aid for Scientific Research (No.) from the Ministry of Education, Science and Culture of Japan to O.N.

References

- Brash, D.E. (1988) UV mutagenic photoproducts in Escherichia coli and human cells: A molecular genetics perspective on human skin cancer, Photochem. Photobiol., 48, 59-66.




- Cleaver, J.E., F. Cortes, L.H. Lutze, W.F. Morgan, A.N. Player and D.L. Mitchell (1987) Unique DNA repair properties of a xeroderma pigmentosum revertant, *Mol. Cell. Biol.*, 7, 3353-3357.
- Cleaver, J.E., F. Cortes, D. Karentz, L.H. Lutze, W.F. Morgan, A.N. Player, L. Vuksanovic and D.L. Mitchell (1988) The relative biological importance of cyclobutane and (6-4) pyrimidine-pyrimidone photoproducts in human cells: Evidence from a xeroderma pigmentosum revertant, *Photochem. Photobiol.*, 48, 41-49.
- Eagle, H. (1959) Amino acid metabolism in mammalian cell culture, *Science*, 130, 432-437.
- Eggset, G., G. Volden and H. Krokan (1987) Characterization of antibodies specific for UV-damaged DNA by ELISA, *Photochem. Photobiol.*, 45, 485-491.
- Franklin, W.A. and W.A. Haseltine (1986) The role of the (6-4) photoproduct in ultraviolet light-induced transition mutations in *E. coli*, *Mutation Res.*, 165, 1-7.
- Galfre, G. and C. Milstein (1981) Preparation of monoclonal antibodies: Strategies and procedures, in: J.J. Langone and H.V. Vunakis (Eds.), *Methods in Enzymology*, 73, Part B, Academic Press, New York and London, pp. 3-46.
- Matsunaga, T., T. Mori and O. Nikaido (1989) Base sequence specificity of a monoclonal antibody binding to (6-4) photoproducts, *Mutation Res.*, In press.
- Mitchell, D.L. and J.M. Clarkson (1984) Use of synthetic polynucleotides to characterize an antiserum made against UV-

- irradiated DNA, Photochem. Photobiol., 40, 743-748.
- Mitchell, D.L., C.A. Haipek and J.M. Clarkson (1985) (6-4) photoproducts are removed from the DNA of UV-irradiated mammalian cells more efficiently than cyclobutane pyrimidine dimers, Mutation Res., 143, 109-112.
- Mitchell, D.L. (1988) The relative cytotoxicity of (6-4) photoproducts and cyclobutane dimers in mammalian cells, Photochem. Photobiol., 48, 51-57.
- Mori, T., T. Matsunaga, T. Hirose and O. Nikaido (1988) Establishment of a monoclonal antibody recognizing ultraviolet light-induced (6-4) photoproducts, Mutation Res., 194, 263-270.
- Mori, T., A.A. Wani, S.M. D'Ambrosio, C.-C. Chang and J.E. Trosko (1989) In situ pyrimidine dimer determination by laser cytometry, Photochem. Photobiol., 49, 523-526.
- Wade, M.H., J.E. Trosko and M. Schindler (1986) A fluorescent photobleaching assay of gap junction mediated communication between human cells, Science, 232, 525-528.

Figure Legends

- Fig. 1. The typical fluorescent images of XP12BE cells irradiated with various doses of UV obtained from the ACAS 470 workstation. The values shown in the pictures are UV doses (J/m^2).
- Fig. 2. The induction of (6-4) photoproducts as measured by the in situ immunofluorescent laser cytometry in XP12BE cells irradiated with various UV doses.
- Fig. 3. The autoradiographs of UV-irradiated and unirradiated normal human cells treated with ^3H -monoclonal antibody. The values shown in the pictures are UV doses (J/m^2).
- Fig. 4. The induction of (6-4) photoproducts as measured by autoradiography in normal and XP1EH cells irradiated with various UV doses. 100 cells were counted at each point.
- Fig. 5. The repair of (6-4) photoproducts in normal and XP1EH cells irradiated with $20 \text{ J}/\text{m}^2$ as measured by autoradiography. 100 cells were counted at each point.
- Fig. 6. The repair of (6-4) photoproducts in normal and XP1EH cells irradiated with $40 \text{ J}/\text{m}^2$ as measured by autoradiography. 100 cells were counted at each point.

Insitu.ms

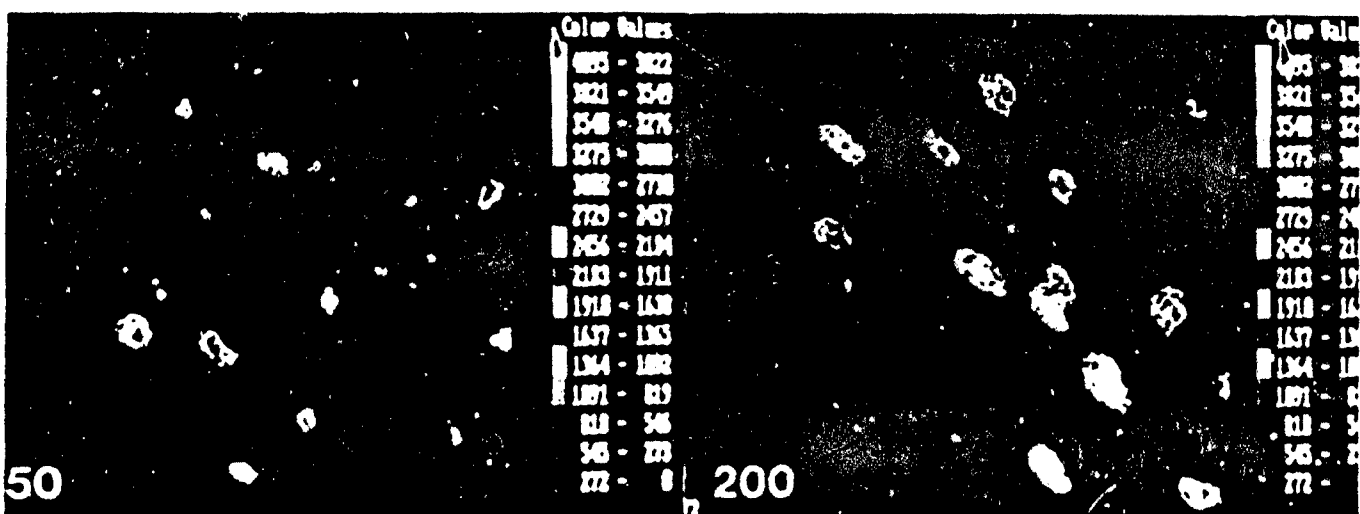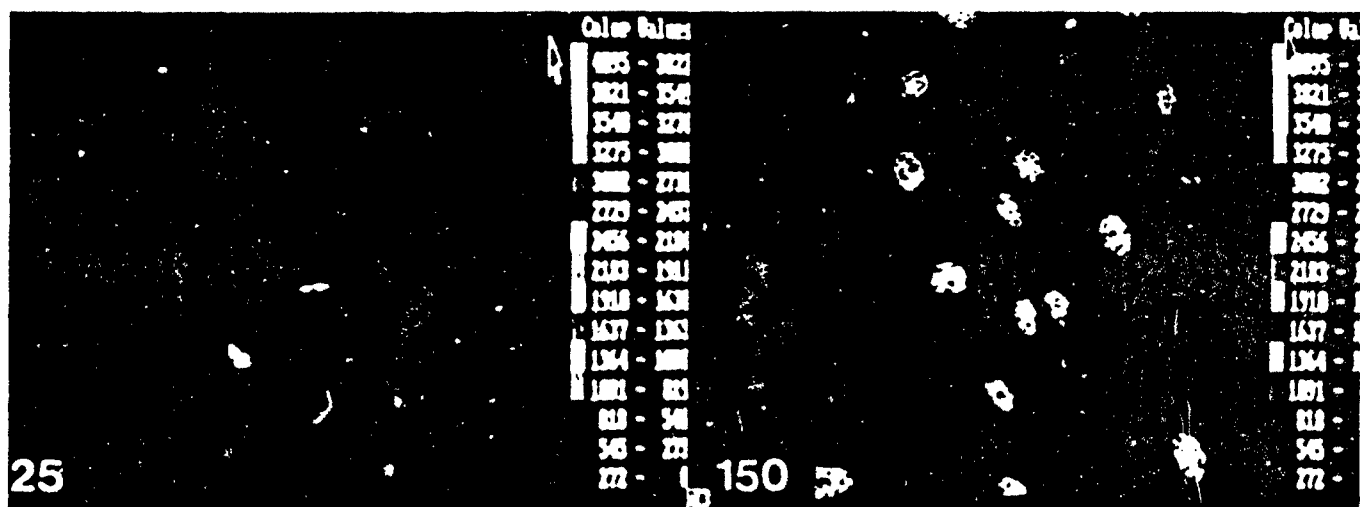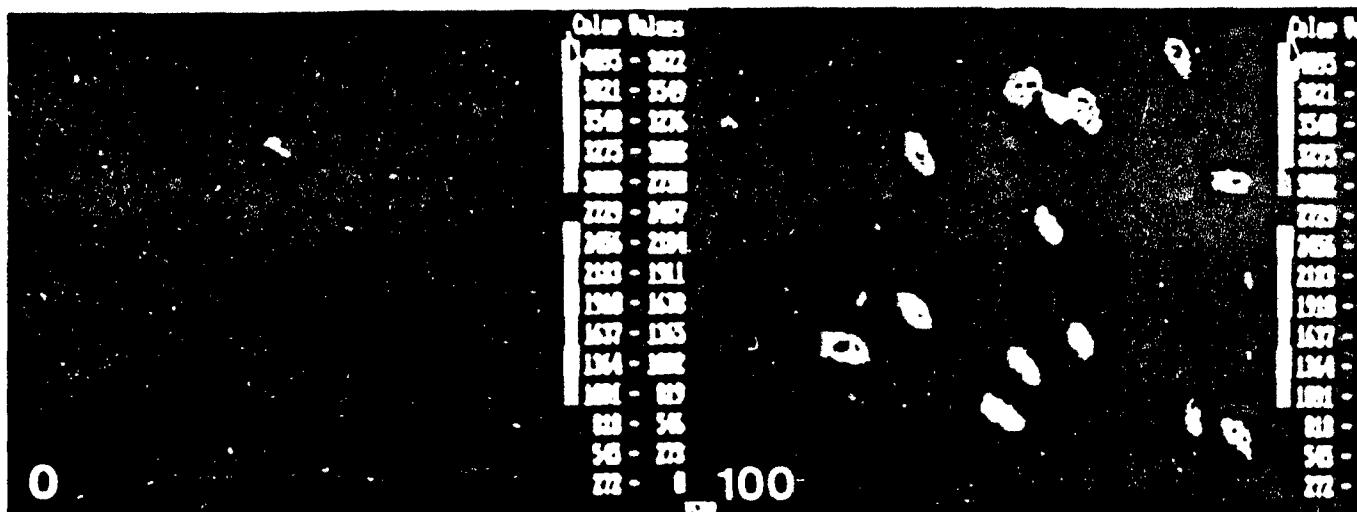

Characterization of an In Vitro Human Kidney
Epithelial System to Study Gap Junctional
Intercellular Communication

B.V. Madhukar, Ph.D.
A.W. de Feijter
C.M. Hasler, Ph.D.
B. Lockwood
S.Y. Oh, Ph.D.*
C.C. Chang, Ph.D.
E. Stanbridge, Ph.D.†
J.E. Trosko, Ph.D.

Department of Pediatrics and Human Development
College of Human Medicine
Michigan State University
East Lansing, Michigan 48824

†Department of Microbiology and Molecular Genetics
College of Medicine
University of California
Irvine, California 92717

*Present Address: Department of Clinical Pharmacology
Flinders Medical Center
Bedford Park, S. Australia 5042

ABSTRACT

In the present study we have investigated if a human cell culture model can be used to study the interaction of xenobiotic chemicals to inhibit gap junctional intercellular communication (GJIC). The G401.2/6TG.1 epithelial cell line derived from the kidney of a Wilms' tumor patient exhibited extensive GJIC as measured using FRAP analysis and scrape-loading dye transfer techniques. Several known hepatic tumor promoting chemicals, skin tumor promoter, TPA and polyhalogenated biphenyls were tested for their effect on IC. TPA completely abolished GJIC of these cells at concentrations as low as 0.3 ng/ml and continued to down modulate cell-cell communication for up to 3 days following a single treatment. Among the other chemicals studied, dieldrin, heptachlor and heptachlorepoxyde were most effective as inhibitors GJIC; DDT, lindane mirex and Fire Master BP-6 were ineffective. Two fatty acids, linoleic acid and decasohexaenoic acid, were also effective in down regulating cell-cell communication. The data also showed that while TMB-8 antagonized TPA's action on GJIC, dieldrin synergistically interacted to abolish GJIC. The results of the study indicate that the response of these cells may be different from other animal derived cell culture models and suggest the need for using human cell culture models to study toxic interactions of xenobiotic chemicals.

INTRODUCTION

In all metazoans, coordination of cell growth, differentiation and

Key Words: Gap junctional intercellular communication, cell-cell communication, heptachlor, heptachlor epoxide, TPA, dieldrin, lindane, mirex, PBB's, PCB's, fatty acids, human kidney epithelial cells, free radicals, FRAP analysis, scrape-loading/dye transfer assay.

development, as well as activities of cells within tissues, are mediated by two major forms of intercellular communication of ionic and molecular messages

(Loewenstein, 1981; Potter, 1983; Snyder, 1985; Pitts and Finbow, 1986; Green, 1988; Hertzberg and Johnson, 1988). One form of intercellular communication depends on the protein structure, the gap junction, found on the plasma membrane between contiguous cells which functions as a channel for ions and molecules of approximately 1000 daltons (Loewenstein, 1979; Evans, 1988). Gap junctions have been shown to be modulated by both endogenous (i.e., hormones, neurotransmitters, growth factors) and exogenous (i.e., biological toxins, drugs, pollutants, nutrients, food additives, pesticides, etc.) agents (Loewenstein and Risinger, 1985; Trosko et al, 1988; Elmore et al, 1987; Neyton and Trautman, 1986). Many cancer cells have reduced or selective gap junctional intercellular communication (Kanno, 1985; Yamasaki and Fitzgerald, 1988). In addition, many chemicals having been shown to be teratogens (Trosko et al, 1982; Welsch and Stedman, 1984), tumor promoters (Yotti et al, 1979; Murray and Fitzgerald, 1979), and neurotoxicants (Trosko et al, 1987) have been linked to their ability to inhibit gap junction function.

To date many techniques, such as electrocoupling (Enomoto et al, 1981), microinjection/dye transfer (Enomoto and Yamasaki, 1984), metabolic cooperation (Yotci et al, 1979; Williams et al, 1981; Kavanagh et al, 1986; Ruch et al, 1987; Jone et al, 1987; Gupta et al, 1985; Hooper, 1981; Davidson et al, 1985; Mosser and Bols, 1982), fluorescence recovery after photobleaching (Wade et al, 1986), and scrape-loading/dye transfer (El-Fouly et al, 1986), have been utilized to measure gap junction function. In addition, many different cell types, strains, and lines have been studied. An emerging picture indicates that not all gap junctions in different cell types respond to the same chemicals in the same manner, possibly due to species, or tissue type receptor differences, as well as differences in gap junction proteins (Saez et al, 1989). While teratogens, tumor promoters and neurotoxicants show species, tissue and developmental stage differences, it would be surprising that a single cell type in vitro system could be predictive of all in vivo conditions (Trosko et al, 1988).

Clearly, if in vitro assays are to be used for risk assessment purposes, as well as for studies of mechanisms of action by which various toxicants work, a series of human epithelial and mesenchymal assays from various organs will eventually be required to minimize the errors of extrapolating data from animal cells to human beings. Consequently, this study was designed to develop a human kidney epithelial cell system to detect chemicals which could modulate gap junction function. The results indicate that this cell line will be an excellent system for further study.

MATERIALS AND METHODS

Cell Culture

The human kidney epithelial cell line C401.2/6TG.1 used in these experiments was originally derived from a Wilms' tumor patient and established as described previously (Weissman et al., 1987). Cells were maintained in D-medium, a modified Eagle's medium containing Earle's balanced salt solution with a 50% increase of vitamins and essential amino acids except glutamine, a 100% increase of non-essential amino acids and 1 mM sodium pyruvate, 5.5 mM glucose, 14.3 mM NaCl and 11.9 mM NaHCO₃ (pH 7.3). The medium was supplemented with 5% fetal bovine serum (Gibco Laboratories) and 50 µg/ml gentamicin (Quality Biological, Inc.). The cells were grown at 37°C and 5% CO₂ in a humidified incubator. Under these conditions C401.2/6TG.1 cells had a doubling time of 27 hours.

Chemicals

TPA (Sigma Chemical Co.) was dissolved in ethanol. Dieldrin, p,p'-DDT, lindane, and mirex were generously provided by Dr. Richard E. Leavitt, Pesticide Research Center, Michigan State University. Heptachlor and heptachlor epoxide, 99+% pure were a generous gift of Velsicol Chemical Corporation, Chicago, IL. Aroclor-1254, Firemaster BP-6, were kindly provided by Dr. Matthew Zabik, Pesticide Research Center, Michigan State University. 2,4,5,2',4',5'- and 3,4,5,3',4',5' hexachlorobiphenyl (HCB) and 2,4,5,2',4',5' hexabromobiphenyl (HBB) were purchased from Analabs, North Haven, CN and were 99+% pure. Stock solutions and subsequent dilutions of the test chemicals were made in ethanol or dimethyl sulfoxide (DMSO) and were added directly to the culture dishes to obtain the required concentrations. Hydrogen peroxide (30% H₂O₂), palmitic acid (PAL, 16:0), linoleic acid (LA, 18:2), linolenic acid (LN, 18:3) and docosahexaenoic acid (DHA, 22:6) were obtained from Sigma Chemical Co. (St. Louis, MO). Stock solutions of fatty acids were dissolved in ethanol and kept frozen under an atmosphere of nitrogen at -20°C. To achieve the desired working concentration of fatty acid, aliquots of stock fatty acids were added directly to complete media supplemented with 5% delipidated serum (Pelfreeze Biologicals) and briefly sonicated immediately prior to treatment of the cells. The final concentration of solvent in the medium was 0.1%. For the scrape-loading/dye transfer (SL/DT) experiments Lucifer Yellow CH (Sigma Chemical Co., MW 457.2) and tetramethyl rhodamine dextran (Molecular Probes, Inc., MW 10,000) were dissolved in phosphate-buffered saline (PBS) at a concentration of 0.5 mg/ml (El-Fouly *et al.*, 1987). For the Fluorescence Redistribution After Photobleaching (FRAP) assay 1 mg of 5,6-carboxyfluorescein diacetate (Molecular Probes, Inc.) was dissolved per ml of ethanol (Wade *et al.*, 1986).

Measurement of Gap Junctional Intercellular Communication (GJIC)

Fluorescence Redistribution After Photobleaching (FRAP). For the FRAP assay, originally described by Wade *et al.* (1986), G401.2/6TG.1 cells were plated at low density in 35 mm plates. After exposure to the test chemicals the cells were rinsed several times with phosphate buffered saline with Ca and Mg (Ca,Mg-PBS) and maintained in 2 ml of Ca,Mg-PBS during FRAP analysis. Fourteen μ l of the dye, 5,6-carboxyfluorescein diacetate were added to each plate. After 15 minutes the plates were rinsed several times with Ca,Mg-PBS. 2 ml of Ca,Mg-PBS were kept on the plates and the cells were subjected to analysis using the Meridian ACAS 470 Workstation (Meridian Instruments, Okemos, MI). Selected cells (5 to 8 per scan) were photobleached and monitored for the return of fluorescence at 5 minute intervals for a total period of 15 minutes (4 scans). The experiments were repeated on 4-5 plates at each concentration or treatment. In each plate a single cell was selected but not photobleached to determine the background decline of fluorescence. Occasionally, an isolated cell was photobleached as a negative control for fluorescence redistribution. Results of the FRAP-experiments were usually expressed as the average (percent) recovery of fluorescence $\pm$ standard error of the mean (SEM). While the short-term treatments were performed at room temperature, for long-term exposure, cells were incubated at 37°C and the medium was changed every 2 days followed by addition of fresh test chemical(s).

Scrape-loading/dye transfer (SL/DT). The SL/DT technique was performed as described previously. After treatment with the test chemicals, confluent

cultures in 35 mm plates were rinsed several times with phosphate buffered saline without Ca or Mg (PBS) and two ml of the dye mixture (0.5 mg/ml of each Lucifer Yellow and tetramethyl rhodamine dextran in PBS) was added. Three or four scrape lines were made on the monolayer with a surgical blade. After 3 minutes, the cells were rinsed several times with PBS and examined for dye transfer under a Nikon epifluorescence microscope and photographs were taken. Because of its low molecular weight Lucifer Yellow can be transmitted between adjacent cells via gap junctions, but it does not diffuse through intact plasma membranes. The high molecular weight rhodamine dextran (M.W. 10,000 daltons) however, can neither diffuse through intact plasma membranes nor cross the junctional channels, and therefore serves to identify the primary loaded cells. For the long-term experiments, the medium was changed every 2 days and test chemicals were added. All experiments were repeated at least once.

RESULTS

The human kidney epithelial cell line G401 2/6TG.1 expresses a very high level of gap junctional intercellular communication as shown by both FRAP (Fig. 1) and scrape-loading/dye transfer (Fig. 2) techniques. Fifteen minutes after photobleaching of selected cells, the average recovery of fluorescence in G401.2/6TG.1 cells was 53.3%. In the SL/DT assay the dye, Lucifer Yellow distribution extended 8 to 10 rows of cells beyond the primary recipient cells during the 3 minutes after scrape-loading.

Quantitation of the SL/DT assay was done following two methods: The semiquantitation method of the effect of the chemicals on GJIC was done by assigning an arbitrary scale for communication in control cells and comparing the level of communication after chemical treatment with that of the controls. In the other method a more precise quantitative analysis was done by counting the number of cells with Lucifer yellow in a 5 cm^2 area of the scrape-loading dye transfer picture.

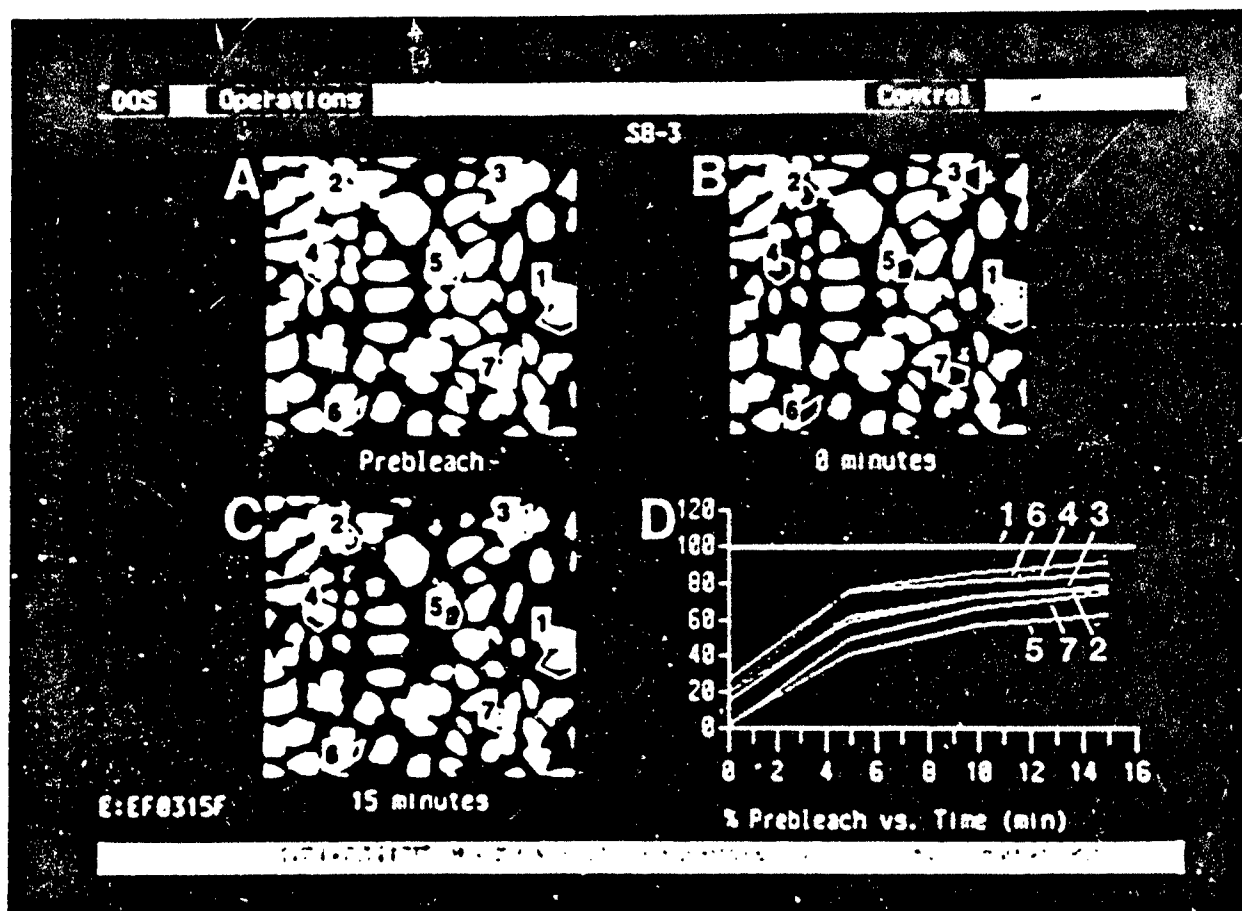


FIGURE 1. Meridian ACAS images of human kidney epithelial C401.2/6TG.1 cells, fluorescently labelled with 5,6-carboxyfluorescein. Digitized pseudo-images (A) prior to, (B) immediately after and (C) 15 minutes after photobleaching of selected cells are displayed. (D) shows a plot in which the percent recovery of fluorescence is expressed in time for each cell selected. Four scans were made at 5 minute intervals after photobleaching. Data were corrected for the background loss of fluorescence in one area (area 1) with two isolated cells. FRAP was performed as described in Materials and Methods.

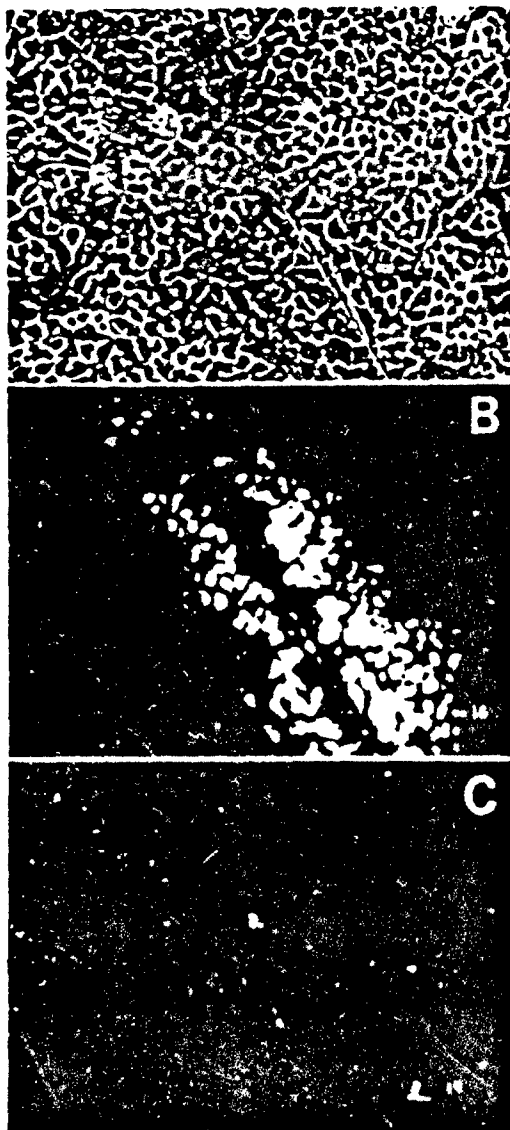


FIGURE 2. CJIC in C401.2/6TG.1 cells as measured by the Lucifer yellow dye transfer technique. Scrape loading of the cells with Lucifer yellow/Rhodamine dextran was done as described. (A) Phase contrast picture of the cells showing normal morphology. Arrow indicates the scrape line for introducing the dye mixture. (B) Lucifer yellow distribution after scrape loading. (C) High mol. wt. Rhodamine dextran is limited to the primary loaded cells only. 300 X

Characterization of TPA-effect on CJIC.

TPA-induced inhibition of intercellular communication in C401.2/6TG.1 cells was characterized using the FRAP assay. Both effects of concentration and exposure time were established. Fig. 3 displays the dose-response curve of TPA

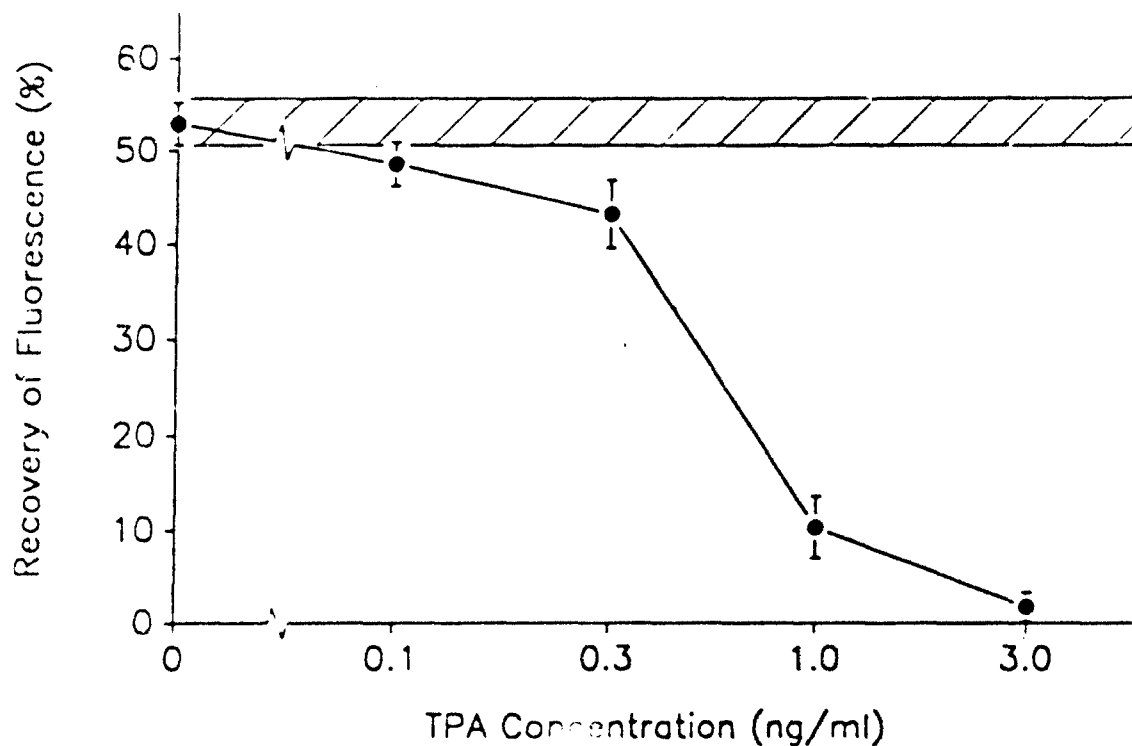


FIGURE 3. Dose-response relationship of TPA-effect on gap junctional communication by FRAP. C401.2/6TG.1 cells were exposed to different concentrations of TPA for 1 hour and fluorescence recovery was determined over a 15-minute period. The shaded area shows the recovery of fluorescence in untreated control cells. Each point in the graph represents the mean recovery $\pm$ SEM of 5 experiments.

In C401.2/6TG.1 cells after a 1 hr. exposure. Up to 0.3 ng/ml the compound had no significant effect on the recovery of fluorescence. From 0.3 up to 3.0 ng/ml communication was inhibited in a dose-dependent manner with complete blockage at the highest dose tested. The recovery of fluorescence after different exposure times of C401.2/6TG.1 cells to 3 ng/ml TPA is shown in Fig. 4. When C401.2/6TG.1 cells were continuously treated with TPA, communication was completely blocked up to 1 day. At continued exposure, a time dependent return was observed up to 40% of the control level after 4 days. In ethanol, the solvent for TPA, showed no effect on GJIC up to 4 days of exposure at the concentration applied (Fig. 4).

Characterization of Dieldrin-effect by scrape-loading/dye transfer.

The neurotoxic insecticide, dieldrin was studied for its inhibitory effect on GJIC in G401.2/6TG.1 cells using the dye-transfer assay. Data presented (Fig. 5) show

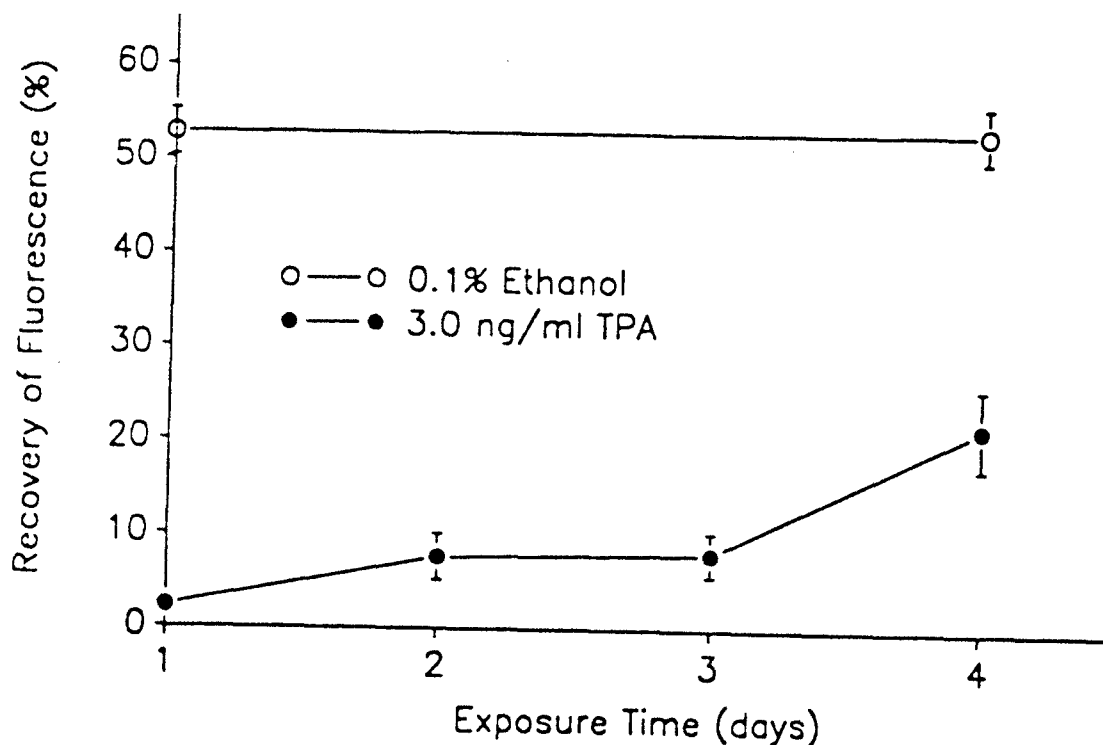


FIGURE 4. Time-response relationship of TPA-effect on gap junctional communication by FRAP. G401.2/6TG.1 cells were treated with 3 ng/ml TPA or the solvent (0.1% ethanol) for 1 to 4 days and fluorescence recovery was determined over a 15-minute period. Each point in the graph represents the mean recovery $\pm$ SEM of 5 experiments.

that intercellular communication was inhibited in a time dependent manner. Loss of communication was elicited above 2 μ g/ml at one hour after treatment and progressively increased at higher doses up to 10 μ g/ml. A progressive loss of GJIC was seen beginning 10 min after treatment with 10 μ g/ml dose and was complete at 1 hour. Inhibition sustained for up to 2 hours, following which, the cells gradually regained communication, although the extent of communication was somewhat lower than in controls at the end of 24 hours.

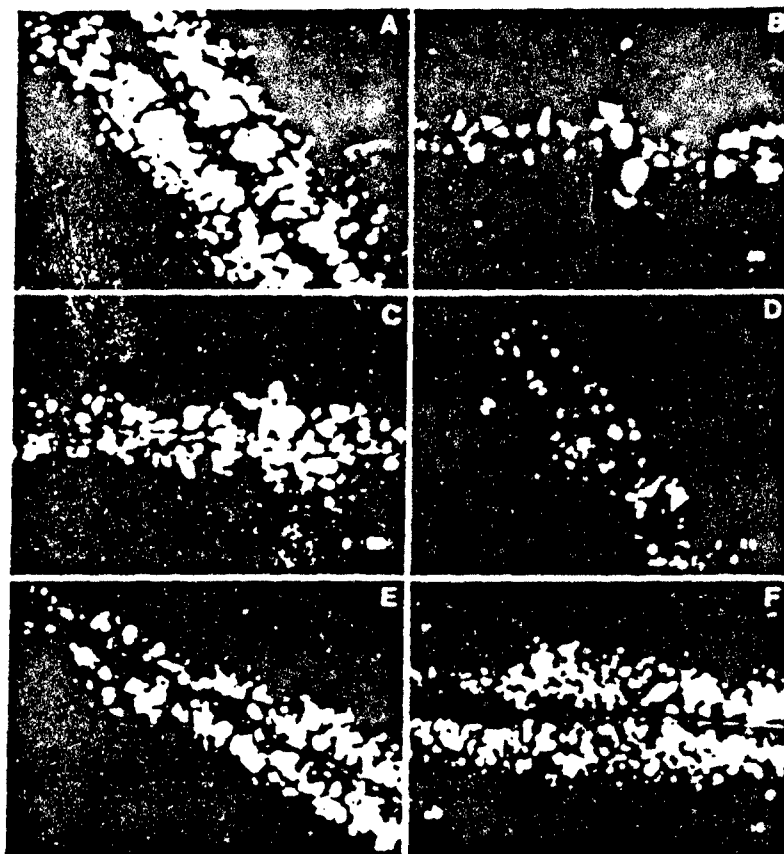


FIGURE 5. Temporal changes in GJIC of C401.2/6TG.1 cells treated with dieldrin. GJIC was determined by the dye transfer method at the required time point after treatment. A. Control. B-F dieldrin (10 $\mu\text{g/ml}$) at 1, 4, 6, 8, 12 and 24 hrs respectively. 300 X

Effects of other xenobiotics

We have also tested a number of xenobiotic chemicals for their ability to inhibit GJIC using the dye transfer assay. Data presented in Table 1 clearly suggest that there are differences in the action of these agents to down regulate GJIC in C401.2/6TG.1 cells. While the chlorinated insecticides, dieldrin, heptachlor and heptachlor epoxide almost completely abolished communication, p,p' DDT, lindane and mirex were ineffective. The polychlorinated biphenyl (PCB) mixture, Aroclor-1254 effectively blocked junctional communication at 10 $\mu\text{g/ml}$, while the polybrominated biphenyl (PBB) mixture, Firemaster BP-6, did not block GJIC. Interestingly, neither of the hexachlore biphenyl isomers nor the

hexabromobiphenyl tested were able to down regulate GJIC.

Effects of Fatty Acids

Confluent monolayers of G401.2/6TG.1 cells were treated with DHA, LA, Ln, and PAL at concentrations of 25, 100 and 250 μ M for various time periods ranging from 5 min up to 24 hr. GJIC was assessed initially utilizing the SL/DT technique. Quantitation of GJIC was done with FRAP analysis. The unsaturated fatty acids of GJIC were done with FRAP analysis. The unsaturated fatty acids DHA and LN significantly inhibited (81%) GJIC at a concentration of 250 μ M after 1 hr as assessed by FRAP analysis (Table 2). LA decreased GJIC slightly (25.5%), but fluorescence recovery values in cells treated with LA did not differ statistically from

TABLE 1

Effect of Several Xenobiotics on GJIC in G401.2/6TG.1 Cells¹.

¹GJIC was measured by scrape-loading/dye transfer technique as described.

²Communication in control cells is given an arbitrary scale of +++, which is equal to the spread of Lucifer yellow into 8-10 rows of cells beyond the primary loaded cells. + indicates significant loss of communication, indicated by spread of Lucifer yellow to 1 layer of cells at the most. The extent of communication was confirmed in 3 different plates in each treatment.

Treatment (ng or μ g/ml)	Relative Communication ²
Control - solvent	++++
TPA - 10 ng/ml	<+
dieldrin - 10 μ g/ml	+
Heptachlor - 10 μ g/ml	+
p,p'-DDT - 10 μ g/ml	+++
Lindane - 10 μ g/ml	++++
Mirex 10 μ g/ml	++++
Aroclor - 1254 10 μ g/ml	+
2,4,5,2',4',5' HCB 10 μ g/ml	++++
3,4,5,3',4',5' HCB 10 μ g/ml	++++
Firemaster BP-6 10 μ g/ml	++++
2,4,5,2',4',5' HBB 10 μ g/ml	++++

TABLE 2

Effect of Various Fatty Acids on Dye Recovery in G401.2/6TG.1 Cells^a^aCells were treated with fatty acids (250 μ M) for 1 hr.^b% dye recovery (normalizing using a transformation procedure) is expressed as means $\pm$ SE for n = 38 (control), n = 7 (LA), n = 9 (DHA), n = 10 (PAL), n = 12 (LN) plates.^cSignificantly different from untreated controls, $p < 0.05$, Tukey's test.

Treatment	Dye Recovery ^b
Control	50.3 $\pm$ 2.5
Linoleic Acid (LA)	37.5 $\pm$ 4.5
Linolenic Acid (LN)	9.6 $\pm$ 1.4 ^c
Palmitic Acid (PAL)	49.6 $\pm$ 4.1
Docosahexaenoic Acid (DHA)	9.3 $\pm$ 2.5 ^c

untreated controls. GJIC in G401.2/6TG.1 cells treated with PAL was similar to that of control cells. The fatty acid-induced inhibition of GJIC by LN and DHA was rapid, occurring within 5 min, and reversible, with GJIC returning to control levels within 1 hr upon removal of fatty acid treatment from cells (data not shown).

Effect of H₂O₂

H₂O₂ concentrations were determined by freshly diluting the 30% stock solution with distilled water and measuring absorbance on a uv/visible spectrophotometer at 240 nm. Exposure of confluent monolayers of SB-3 cells to 100, 250, 500, 1000 and 2000 μ M H₂O₂ in serum-free D media for 1 hr had no effect on GJIC. These findings differ with results obtained utilizing liver epithelial WB-F344 cells, (Hasler, unpublished observations) in which H₂O₂ at concentrations of 250 and 500 μ M had inhibited GJIC.

Antagonists to chemical inhibition of gap junctional communication.

In a previous study we have reported that in the rat liver epithelial cell line, WB F-344, TMB-8 negated the inhibitory action of TPA on GJIC (Oh et al, 1988). Since the G401.2/6TG.1 cell line appeared to respond differently from the rat liver cell line to the action of xenobiotics, we investigated if TMB-8 would exert the same protective action in G401.2/6TG.1 cells by preventing the down modulatory effect of TPA. As shown in Fig. 6, when the G401.2/6TG.1 cells were pretreated with TMB-8 (50 μ M) 10 min prior to treatment with TPA (10 ng/ml), GJIC was only minimally reduced as opposed to complete inhibition by TPA alone. In contrast, TMB-8 appeared to have no effect on dieldrin inhibited GJIC under the experimental conditions used. TMB-8 itself, under these conditions had no effect on GJIC (data not shown). However, the protective action of TMB-8 appeared to be transient under these conditions, as we have noted complete loss

appeared to be transient under these conditions, as we have noted complete loss of communication after 2 hours in cells treated with TMB-8 and TPA.

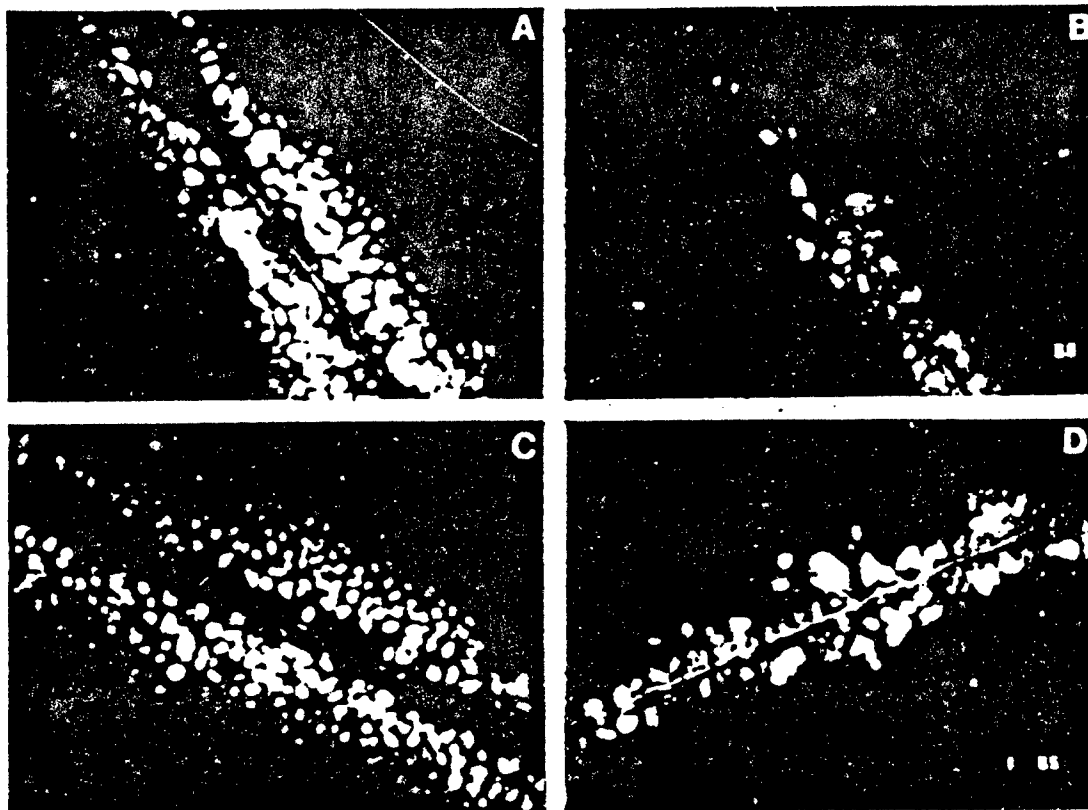


Fig. 6 Antagonistic action of TMB-8 on the inhibition of CJIC by TPA or dieldrin in C401.2/6TC.1 cells. Communication was assayed by Lucifer yellow distribution after scrape loading. Treatment with TMB-8 was done 10 min before adding TPA or dieldrin. A. TMB-8 (50 μ M). B. TPA 100 ng/ml. C. TMB-8 + TPA 100. D. TMB-8 + dieldrin (10 μ g/ml).

Synergism between TPA and dieldrin.

Synergism between chemicals is a phenomenon by which two mechanistically different chemicals potentiate each others action when added together at levels at which either agent alone has only minimal effect. We have examined if TPA and dieldrin exhibit a synergistic interactions to down-modulate CJIC. As presented in Fig. 7 (A-C) TPA at 1 ng/ml or dieldrin at 1 μ g/ml had only a marginal effect in abolishing CJIC under the experimental conditions used. When added together at the same concentrations, under similar conditions, communication was significantly inhibited (Fig. 7D) by these agents. The data thus suggest that these agents exert their inhibitory action of CJIC through divergent mechanisms.

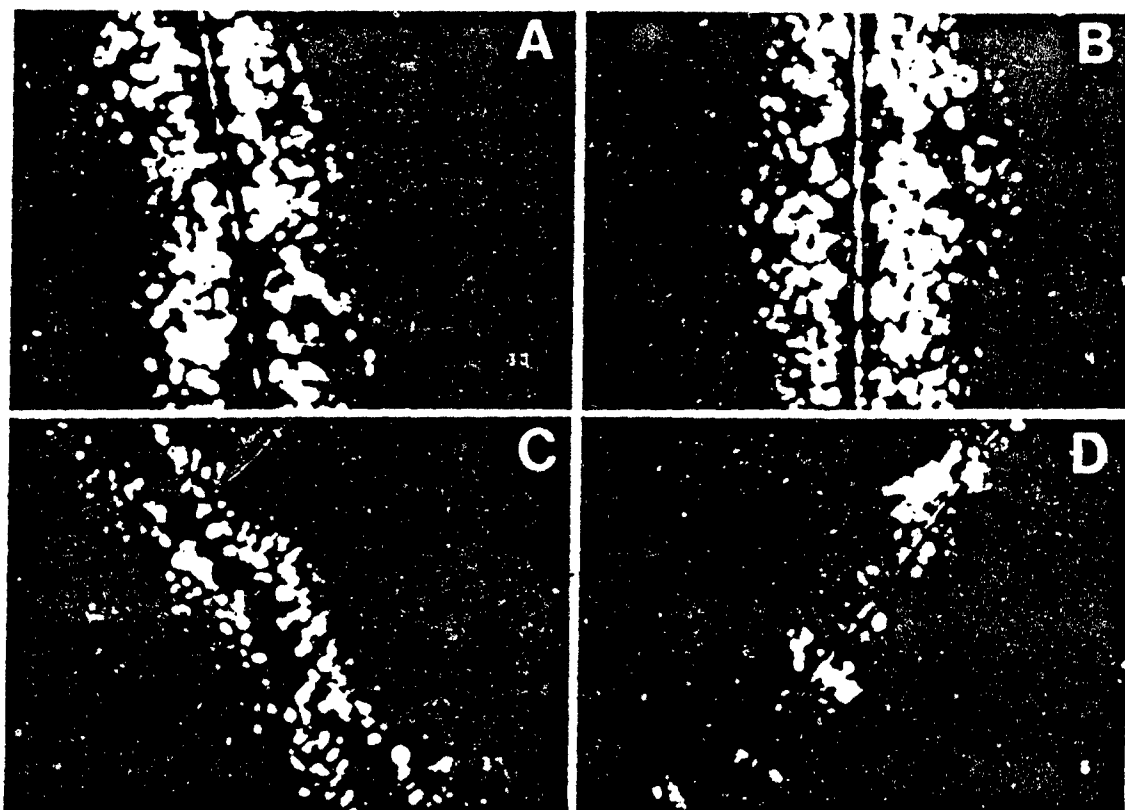


Fig. 7 Synergistic interaction between TPA (1.0 ng/ml) and dieldrin (1 μ g/ml) to inhibit GJIC in G401.2/6TG.1 cells. Communication was determined 1 hr after treatment by the dye transfer technique as described. Lucifer yellow distribution in A. Control. B. Dieldrin (1 μ g/ml). C. TPA (1.0 ng/ml). D. TPA + dieldrin. 300 X

DISCUSSION

In the present study we have investigated the nature of modulation of GJIC in a human tissue derived cell culture system, with a view to develop a *in vitro* model to study the toxic interactions of xenobiotic chemicals. Clearly, *in vitro*

systems derived from human tissue are advantageous to draw meaningful conclusions on the possible effects of xenobiotic exposure of humans. The G401.2/6TG.1 cell line is derived from a Wilm's nephroblastoma tumor and is of epithelial origin. Since gap junctional intercellular communication has been implicated to have a major role in embryonic development, normal cellular growth and proliferation, perturbations in the normal pattern of GJIC by xenobiotics can be of significant importance in determining the potential toxicity of chemicals to humans. In this context, the G401.2/6TG.1 cell line is an excellent model, since these cells have extensive GJIC and disturbances in the extent of intercellular communication by chemical toxicants can be evaluated with relative ease although the results are not a true replication of the situation with primary cells of human origin. The extrapolation of these results should be considered with caution.

The results of the present study clearly suggest that while GJIC in this system is down regulated by xenobiotics, not all the chemicals studied inhibited intercellular communication in these cells. While TPA, a well known tumor promoting agent and a potent inhibitor of gap junctional function, blocked dye transfer of G401.2/6TG.1 cells, there were clear differences in the response of several halogenated hydrocarbon pesticides. Dieldrin and heptachlorepoxyde, two neurotoxic insecticides, inhibited GJIC but three other neurotoxic chemicals, DDT, lindane and mirex, did not have any effect. Previous studies from this laboratory, using V79 Chinese hamster lung fibroblasts and WB rat liver epithelial cells, have reported a number of these chemicals which were ineffective in G401.2/6TG.1 cells that have abolished metabolic cooperation and/or dye transfer in these cell lines (Trosko et al, 1987; Suter et al, 1987; Rezebek et al, 1987; Tsushimoto et al, 1987).

An interesting observation in the present study was the difference between the action of TPA and dieldrin in inhibiting GJIC. In a previous study using rat liver epithelial cell system, WB F-344, we noted that TPA's inhibitory action on GJIC was transient while that of dieldrin was sustained for more than 24 hrs (Oh et al, 1988). In contrast, in G401.2/6TG.1 cells, TPA down modulated gap junction function for more than 2 days after a single treatment whereas the effect of dieldrin was reversed as early as 4 hours. It is well known that many of the biochemical effects of TPA are mediated by the activation of a calcium- and phospholipid-dependent protein kinase, protein kinase C (PKC) (Jong and Blumberg, 1989). This kinase is activated upon translocation into the plasma membrane from the cytosol. It is likely that in the case of G401.2/6TG.1 cells, TPA translocated PKC remains plasma membrane bound for a longer period of time than it is in other cell lines including the WB cells (Oh et al, 1988). Indeed, our preliminary observation on the translocation of PKC indicates this to be the case (unpublished observations). While the inhibitory action of TPA on GJIC can be explained by activation and translocation of PKC in G401.2/6TG.1 cells, the mechanism by which the neurotoxic agent, dieldrin, inhibited intercellular communication is less clear. Since increased intracellular free calcium, (Ca^{++})_i has been implicated in junctional uncoupling, it is possible that dieldrin treatment might lead to increased [Ca^{2+}]_i through inhibition of the Ca^{2+} pump.

In the present study, we have also shown that the unsaturated fatty acids LN and LA differ in their action on GJIC in these cells. In a previous study Aylsworth et al have shown that in Chinese hamster lung fibroblast cell line V79, many of these unsaturated fatty acids have abrogated metabolic cooperation again emphasizing the differences in the response of various cell types to the same chemical(s). Since the unsaturated fatty acids are of considerable health concern their action on GJIC may be significant as a biochemical marker.

In summary, whatever mechanisms might be involved in the uncoupling of gap junctions by xenobiotics, it is clear from our present study that there are

intrinsic differences in the response of various tissue cell types to the inhibitory action of chemicals, thus warranting the need to study more than one cell culture system in order to ascertain the toxic potential of xenobiotics. Moreover, with the demonstration that there is a family of gap junction proteins (see Saez et al, 1989), there is clearly a need to study how specific gap junction proteins might be regulated by cellular responses to different xenobiotics.

ACKNOWLEDGEMENTS

Supported by a grant from the Air Force Office of Scientific Research (AFOSR-89-0325), a NIEHS Superfund grant ES-04911, a gift from R.J. Reynolds and a grant from the National Cancer Institute (CA19104). The authors wish to thank Mrs. Darla Conley for her excellent secretarial skills in preparing the manuscript.

REFERENCES

- AYLSWORTH, C.F., WELSCH, C.W., KABARA, J.J., and TROSKO, J.E. (1987). Effects of fatty acids on gap junctional communication: Possible role in tumor promotion by dietary fat. *Lipids* 22, 445-454.
- DAVIDSON, J.S., BAUMGARTEN, I.M., and Harley, E.H. (1985). Effects of TPA and retinoids on intercellular junctional communication measured with a citrulline incorporation assay. *Carcinogenesis* 6, 645-650.
- EL-FOULY, M.H., TROSKO, J.E., and CHANG, C.C. (1987). Scrape-loading and dye transfer: A rapid and simple technique to study gap junctional intercellular communication. *Exp. Cell Res.* 168, 422-430.
- ENOMOTO, T., SASAKI, Y., SHIBA, Y., KANNO, Y., and YAMASAKI, H. (1981). Tumor promoters cause a rapid and reversible inhibition of the formation and maintenance of electrical cell coupling in culture. *Proc. Natl. Acad. Sci.* 78, 5628-5632.
- ENOMOTO, T., and YAMASAKI, H. (1984). Lack of intercellular communication between chemically transformed and surrounding nontransformed BALB/c 3T3 cells. *Cancer Res.* 44, 5200-5203.
- EVANS, W.H. (1988). Gap junctions: Towards a molecular structure. *BioEssays* 8, 3-6.
- GREEN, C.K. (1988). Evidence mounts for the role of gap junctions during development. *BioEssays* 8, 7-10.
- GUPTA, R.S., SINGH, B., and STETSKO, D.K. (1985). Inhibition of metabolic cooperation by phorbol esters in a cell culture system based on adenosine kinase deficient mutants of V79 cells. *Carcinogenesis* 6, 1359-1366.
- HART, R.W., and HOERGER, F.D., eds. (1988). Banbury Report 31: Carcinogen Risk Assessment: New Directions in the Qualitative and Quantitative Aspects. Cold

Spring Harbor Laboratory, Cold Spring Harbor, New York, pp. 139-168.

HERTZBERG, E.L., and JOHNSON, R.G., eds. (1988). Modern Cell Biology, Vol. 7, Gap Junctions. Alan R. Liss, Inc., New York.

HOOPER, M.L., and SUBAK-SHARPE, J.H. (1981). Metabolic cooperation between cells. Intern. Rev. Cytol. 69, 45-105.

JENG, A.Y., and BLUMBERG, P.M. (1989). Biochemical mechanisms of action of the phorbol ester class of tumor promoters. In: The Pathobiology of Neoplasia. Sirica, A.M., ed., Plenum Press, New York, pp. 371-383.

JONE, C., TROSKO, J.E., and CHANG, C.C. (1987). Characterization of a rat liver epithelial cell line to detect inhibitors of metabolic cooperation. In Vitro Cell Devel. Biol. 23, 214-220.

KANNO, Y. (1985). Modulation of cell communication and carcinogenesis. Japan. J. Physiol. 35, 693-707.

KAVANAGH, T.J., CHANG, C.C., and TROSKO, J.E. (1986). Characterization of a human teratocarcinoma cell assay for inhibitors of metabolic cooperation. Cancer Res. 46, 1359-1366.

LANGENBACH, R., and ELMORE, E., eds. (1988). Tumor Promoters: Biological Approaches for Mechanistic Studies and Assay Systems. Raven Press, New York, pp. 131-147.

LARSEN, W.J., and RISINGER, M.A. (1985). The dynamic life histories of intercellular membrane junctions. Mod. Cell Biol. 4, 151-216.

LOEWENSTEIN, W.R. (1979). Junctional intercellular communication and the control of growth. Biochim. Biophys. Acta 560, 1-65.

LOEWENSTEIN, W.R. (1981). Junctional intercellular communication. The cell to cell membrane channel. Physiol. Rev. 61, 829-913.

MATSUMURA, F., and NARAHASHI, T. (1971). ATPase inhibition and electrophysiological changes caused by DDT and related neuroactive agents in lobster nerve. Biochem. Pharmacol. 20, 825-837.

MILMAN, H.A., and ELMORE, E., eds. (1987). Biochemical Mechanisms and Regulation of Intercellular Communication. Princeton Scientific Publ., Inc., Princeton, NJ, pp. 265-292.

MOSSER, D.D., and BOLS, N.C. (1982). The effect of phorbols on metabolic cooperation between human fibroblasts. Carcinogenesis 3, 1207-1212.

MURRAY, A.W., and FITZGERALD, D.J. (1979). Tumor promoters inhibit metabolic cooperation in cocultures of epidermal and 3T3 cells. Biochem. Biophys. Res. Commun. 91, 395-401.

NEYTON, J., and TRAUTMANN, A. (1986). Acetylcholine modulation of the conductance of intercellular junctions between rat lacrimal cells. J. Physiol. 377, 283-295.

PITTS, J.D., and FINBOW, M.E. (1986). The gap junction. *J. Cell Sci.* 4, 239-266.

POTTER, V.R. (1983). Cancer as a problem in intercellular communication: Regulation by growth-inhibiting factors. *Prog. Nucleic Acid Res.* 29, 161-173.

REZABEK, M.S., TROSKO, J.E., JONE, C., and SLEIGHT, S.D. (1988). Effects of hepatic tumor promoters, phenobarbital and polybrominated biphenyls on intercellular communication between rat liver epithelial cells. *In Vitro Toxicol.* 2, 45-58.

RUCH, R.J., KLAUNIG, J.E., and PEREIRA, M.A. (1987). Inhibition of intercellular communication between mouse hepatocytes by tumor promoters. *Toxicol. Appl. Pharm.* 87, 111-120.

SNYDER, S.H. (1985). The molecular basis of communication between cells. *Sci. Amer.* 253, 132-141.

SUTER, S., TROSKO, J.E., EL-FOULY, M.H., LOCKWOOD, L.R., and KOESTNER, A. (1987). Dieldrin inhibition of gap junctional intercellular communication in rat glial cells as measured by the fluorescence photobleaching and scrape-loading/dye transfer assays. *Fund. Appl. Toxicol.* 9, 785-794.

TROSKO, J.E., Chang, C.C., and NETZLOFF, M.L. (1982). The role of inhibited cell-cell communication in teratogenesis. *Teratogen. Carcinogen. Mutagen.* 2, 31-45.

TROSKO, J.E., JONE, C., and CHANG, C.C. (1987). In vitro detection of chemicals which inhibit gap junctional-mediated intercellular communication: Aldrin, dieldrin and toxaphene as examples. *Molec. Toxicol.* 1, 83-93.

TSUSHIMOTO, G., TROSKO, J.E., CHANG, C.C., and AUST, D. (1982). Inhibition of metabolic cooperation in Chinese hamster cells in culture by various polybrominated biphenyl (PBB) congeners. *Carcinogenesis* 3, 181-185.

WADE, M.H., TROSKO, J.E., and SCHINDLER, M. (1986). A fluorescent photobleaching assay of gap junction-mediated communication between human cells. *Science* 232, 525-528.

WELSCH, F., and STEDMAN, D.B. (1984). Inhibition of metabolic cooperation between Chinese hamster V79 cells by structurally diverse teratogens. *Teratogen. Carcinogen. Mutagen.* 4, 283-301.

WILLIAMS, G.M., TELANG, S., and TONG, C. (1981). Inhibition of intercellular communication between liver cells by the liver promoter, DDT. *Cancer Lett.* 11, 339-344.

YOTTI, L.P., CHANG, C.C., and TROSKO, J.E. (1979). Elimination of metabolic cooperation in Chinese hamster cells by a tumor promoter. *Science* 206, 1089-1091.


Corresponding Author

B.V. Madhukar, Ph.D.
Department of Pediatrics/
Human Development
Michigan State University
E. Lansing, Michigan
48824
Tel: 517/353-5431

APPENDIX C

New technique to measure heterologous and homologous cell-cell communication using the ACAS:

1. Actual ACAS-Image data.
2. Abstract.

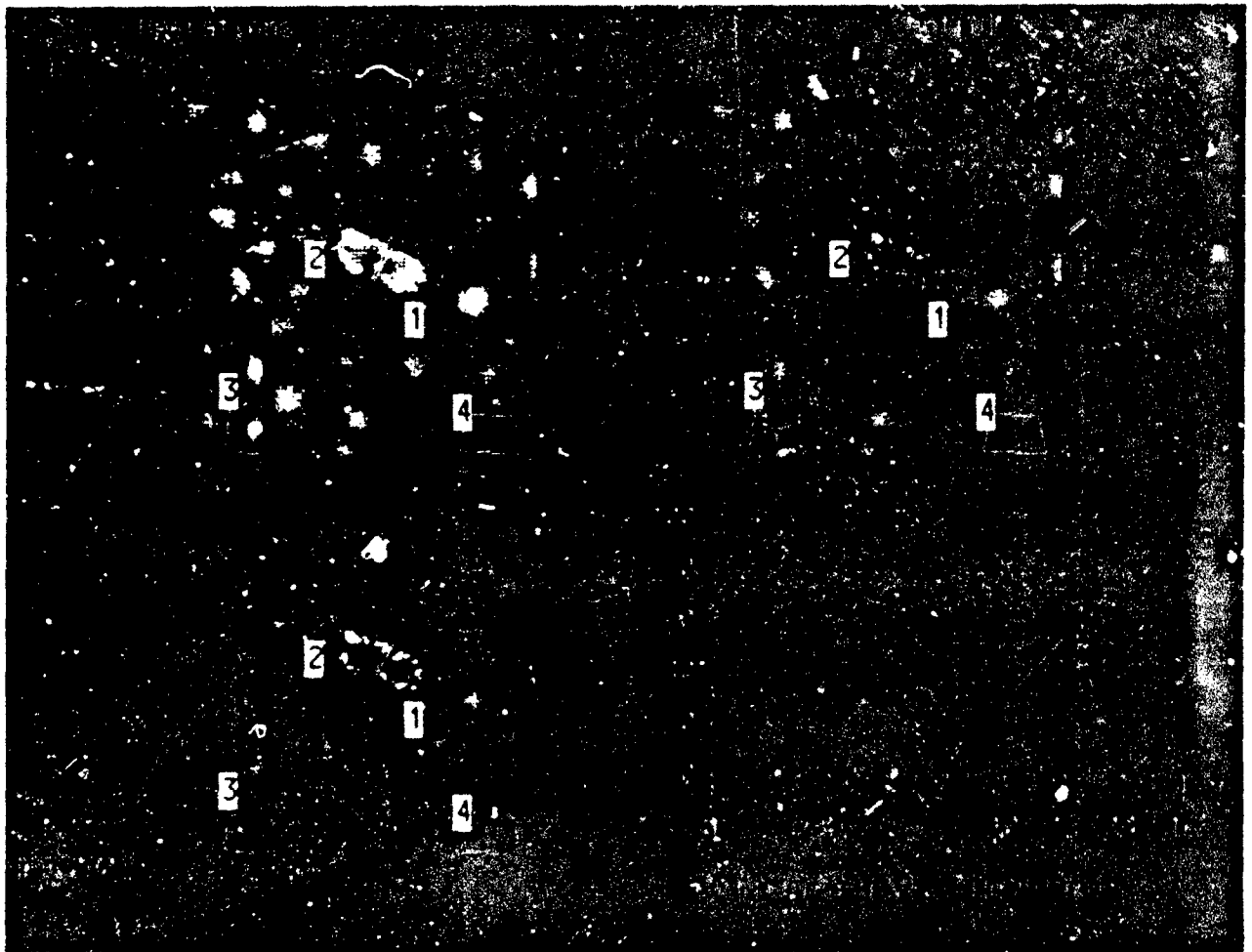

JAN 15, 1990 10:17 AM

Homologous and Heterologous Communication in Co-cultures of Normal and raf-transfected Rat Liver Epithelial Cells. Panels show FRAP analysis, using the Meridian ACAS-570 laser cytometer, of homologous communication between normal RLE cells as indicated by cell #3 and heterologous communication between normal and raf transformed RLE cells as indicated by cells #2 and #3. V-raf transformed cells are loaded with 0.7 um fluorescent beads and then mixed with an excess of normal rat liver epithelial cells. Cells #1 and #2 are two such bead loaded cells which are surrounded by normal RLE cells. Cell #3 is a normal RLE cell. Cell #4 is an unbleached control.

DEADLINE FOR RECEIPT:
DECEMBER 8, 1989

ENCLOSE:

- Check or Purchase Order for US\$25.00, payable to American Association for Cancer Research, Inc.
- Original Form plus 20 photocopies (Items 1-8 below must be answered completely and signed where requested on all copies.)
- Acknowledgment Card (completed and stamped)
- A self-addressed, stamped envelope

DO NOT FOLD

FOR OFFICE
USE ONLY

Perm. Abst. No.
PO CK Other NP

1990 AACR ABSTRACT FORM

TEMPORARY ABST. NO. 4559

Type abstract within blue lines below. See sample abstract.

Gap junctional intercellular communication in raf and raf/myc transformed rat liver epithelial cell lines. G.H. Kalimi¹, J.E. Trosko¹, L.L. Hampton¹, S.S. Thorgeirsson², and A.C. Huggett². Michigan State Univ., East Lansing, MI 48824 and ²National Cancer Inst., Bethesda, MD 20892.

Rat liver epithelial cell lines transformed with v-raf and v-raf/v-myc represent a good in vitro model of tumor progression (Worland et al., 1990, Mol. Carcinogenesis, in press). We studied homologous gap junctional intercellular communication (GJIC) in these cell lines using dye transfer and fluorescence recovery after photobleaching (FRAP) assays, and heterologous GJIC with normal cells using a new bead labelling-FRAP assay. GJIC was found to be reduced in both series of cell lines and showed a good correlation with increasing tumorigenicity. Northern analysis showed reduced levels of connexin 43 in cell lines exhibiting lower GJIC. While raf transformed cell lines showed good heterologous GJIC with normal RLE cells, the raf/myc transformed cells had complete absence of heterologous communication. These results imply the cooperation of the two oncogenes in abolishing heterologous communication and further suggest the involvement of both homologous and heterologous communication in tumorigenicity of RLE cell lines.

Please verify that you have adhered to all AACR abstract submission rules. Infractions of authorship regulations will result in rejection and return of the abstract to the SPONSOR without consideration by the Program Committee. Infractions of format, including abstracts submitted without appropriate signatures, will result in a fee of \$35 which will be billed to the SPONSOR.

1. Category and Subclassification:

(Type the two-digit code for one sub-classification in the space provided.)

Biochemistry BC: — —

Biology BL: — —

Carcinogenesis CG: 0 8

Clinical Investigations CI: — —

Endocrinology EN: — —

Epidemiology EP: — —

Immunology IM: — —

Molecular Biology/Virology MV: — —

Preclinical Pharmacology
Experimental Therapeutics PT: — —

2. Optional Minisymposium Preference:

(If presentation in a minisymposium is preferred, enter the appropriate code number in the space below. Final decisions must be made at the discretion of the Program Committee.)
Minisymposium Code 24

3. Abstract is SPONSORED by

James E. Trosko, Ph.D.

(typed name of SPONSOR)

B240 Life Sciences Bldg., M.S.U. E. Lansing, MI 48824

(517) 353-6346

(address)

(city, state, zip code, country)

(area code, telephone number)

4. Abstract is to be PRESENTED by

Ghulam Kalimi, Ph.D.

(typed name of PRESENTER)

B240 Life Sciences Bldg., M.S.U. E. Lansing, MI 48824

(517) 353-5431

(address)

(city, state, zip code, country)

(area code, telephone number)

5. ☒ The PRESENTER of this abstract is a predoctoral student, postdoctoral fellow, or physician in training

6. As the SPONSOR of this abstract, on behalf of all of the authors I hereby indicate my support for the data contained herein and transfer its copyright to the American Association for Cancer Research, Inc.

Signature of SPONSOR

James E. Trosko

7. As the PRESENTER of this abstract, on behalf of all of the authors I hereby give permission to the American Association for Cancer Research, Inc., to record my presentation at its 41st Annual Meeting for audiocassette sales during and after the meeting.

Signature of PRESENTER

Ghulam Kalimi

8. ☐ The PRESENTER and SPONSOR of this abstract is an Associate Member. The content of the abstract is endorsed by