



Flexural Fatigue Behavior of New All-Ceramic CAD/CAM Materials

Daniel A. Reid, DMD
Capt, USAF, DC
23 Mar 2022

Disclaimer

- The views expressed are those of the authors and do not reflect the official views or policy of the Uniformed Services University, Department of Defense, or its Components. The authors do not have any financial interest in the companies whose materials are discussed in this abstract. The views of Dentsply/Sirona, Ivoclar Vivadent, Kuraray, or VITA are not necessarily the official views of, or endorsed by, the U.S. Government, the Department of Defense, or the Department of the Air Force. No Federal endorsement of Dentsply/Sirona, Ivoclar Vivadent, Kuraray, or VITA is intended.

Co-Contributors

- Wen Lien, DDS, MS,
 - Col, USAF, DC
- Christopher Raimondi, DDS, MS
 - Col, USAF, DC
- Stephen Arnason, DDS, MS
 - Maj, USAF, DC
- Kraig Vandewalle, DDS, MS
 - Col (Ret), USAF, DC



Outline

- Background
- Objectives
- Materials and Methods
- Results
- Conclusion

Background

- Increased demand for efficiency, productivity, and naturalistic esthetics
 - Ceramics
 - Lithium disilicate
 - Zirconia oxide

Lithium Disilicate

- Glass ceramic (Kelly, 2008)
 - Randomly organized crystals
 - Deflect crack propagation (Shenoy and Shenoy, 2010)
 - Provide adequate flexural strength
 - As high as 500MPa (IPS e.max, 2019)

Zirconia Oxide

- Polycrystalline oxide material
 - Exists in three phases, each with distinct properties
 - Monoclinic, tetragonal, and cubic (Sanal and Killinc, 2020)
 - Tetragonal (best properties) stabilized with yttrium oxide (Jansen et al, 2019)
 - Excellent mechanical properties
 - 1200 MPa for tetragonal (Jansen et al, 2019)
 - Questionable esthetic properties
 - Esthetics can be improved at the expense of mechanical properties (Jansen et al, 2019)

New materials

- Tessera
 - New novel lithium disilicate and lithium aluminum silicate (Virgilite) material
 - Reported material has a crystallization time of only 4.5 minutes in the CEREC SpeedFire
 - Reported biaxial flexural strength greater than 700 MPa in a 1.0mm preparation design



New materials

- IPS e.max ZirCAD Prime and MT Multi
 - Gradient technology is a unique manufacturing process that uses special powder conditioning to combine 3Y-TZP and 5Y-TZP or 4Y-TZP and 5Y-TZP
 - Seamless progression of strength with 3Y-TZP or 4Y-TZP in the dentin zone and 5Y-TZP in the incisal or occlusal translucent zone



Objective

- The purpose of this study was to compare the flexural fatigue behavior of new chairside and lab CAD/CAM materials to established materials on the market.

Null Hypotheses

- There will be no difference in in flexural strength, flexural fatigue behavior, and endurance between the different ceramic materials.

Materials and Methods

- Materials

- 2 lithium disilicates

- IPS e.max CAD (control) (Ivoclar Vivadent)
 - Tessera (Dentsply)

- 4 zirconia oxides

- IPS e.max ZirCAD MT (control) (4Y-TZP) (Ivoclar Vivadent)
 - Katana STML (4Y-TZP) (Kuraray Noritake)
 - IPS e.max ZirCAD MT Multi (4Y-TZP/5Y-TZP) (Ivoclar Vivadent)
 - IPS e.max ZirCAD Prime (3Y-TZP/5Y-TZP) (Ivoclar Vivadent)
 - VITAYZ ST (4Y-TZP) (VITA Zahnfabrik)




Materials and Methods

- Specimen preparation
 - Specimen dimensions:
 - Length and width: 18.0 mm x 4.0 mm
 - Thickness: 1.3 mm
 - 30 specimens were fabricated for each group(7)
 - 90 specimens for Tessera
 - 30 no glaze regular fire
 - 30 glazed regular fire
 - 30 glazed speed fire




Materials and Methods

- Specimen preparation
 - Lithium disilicate materials
 - Milled from blocks in a IsoMet 5000, (linear precision saw, Buehler)
 - Crystallized according to manufacturers specifications
 - Zirconia oxide materials
 - Specimens were milled from zirconia discs
 - The two gradient materials was centered within the gradient
 - Sintered according to manufacturers specifications

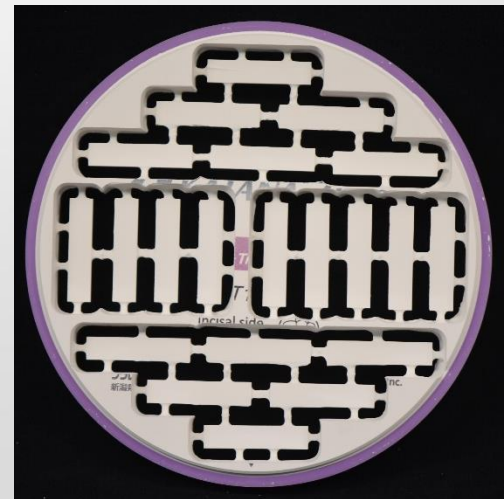

Materials and Methods

- Flexural Strength
 - Fifteen beam specimens per group
 - Fractured in universal machine
 - 3-point bending-test device
 - Central load were applied with a head diameter of 2.0mm
 - Flexural strength calculation
 - $FS = 3Fl/2bd^2$

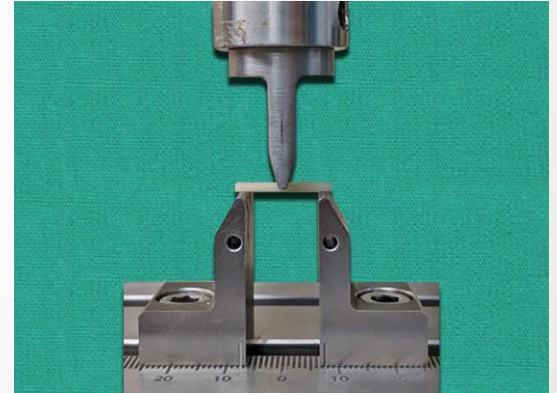

Flexural Fatigue Strength

- Staircase method
 - Fifteen beams per group
 - 3-point test device
 - Load profile was determined based on the single load to fracture data
 - 100,000 cycles until specimen failure or survival at 20 Hz
 - 20% increase/decrease of initial starting force value
 - Based on survival or failure



$$\sigma_{ff} = X_0 + d \left(\frac{\sum in_i}{\sum n_i} \pm 0.5 \right)$$


$$SD = 1.62d \left(\frac{\sum n_i \sum i^2 n_i - (\sum in_i)^2}{(\sum n_i)^2} + 0.029 \right)$$

Data Analysis

- Means and standard deviation for flexural strength was calculated for each ceramic material
 - Kaplan-Meier survival test
 - Analyzed using one-way ANOVA and Tukey's post hoc tests to ($\alpha=0.05$)
- Weibull Modulus


Results

- Significant differences between groups
 - Silica based
 - IPS e.max CAD vs Tessera
 - Flexural strength of e.max 384.14 compared to 206.07-243.61 for Tessera
 - No significant differences in flexural fatigue limits

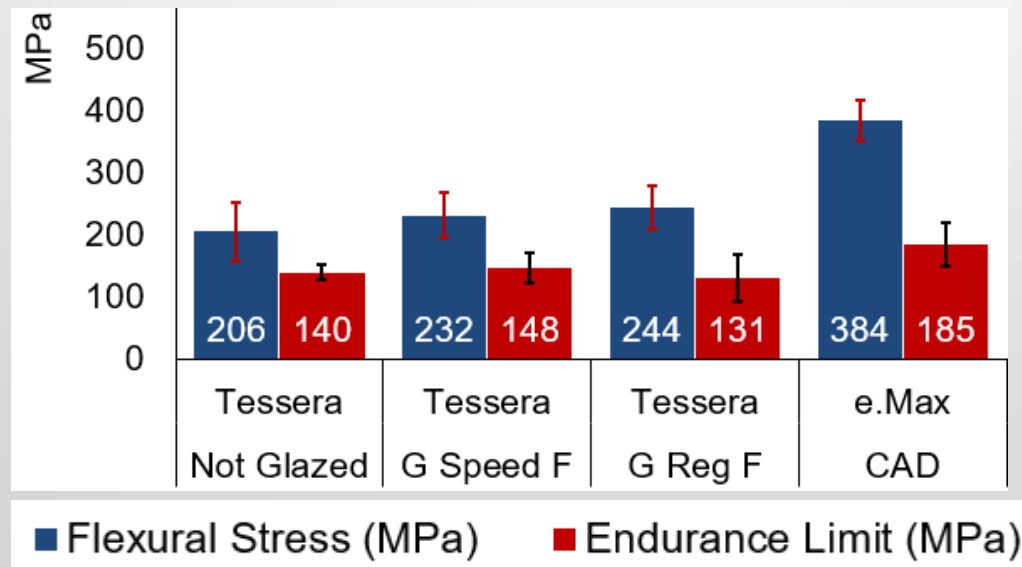

Results

- Significant differences between groups
 - Oxide-based ceramics
 - VitaYZ had the greatest flexural strength and fatigue

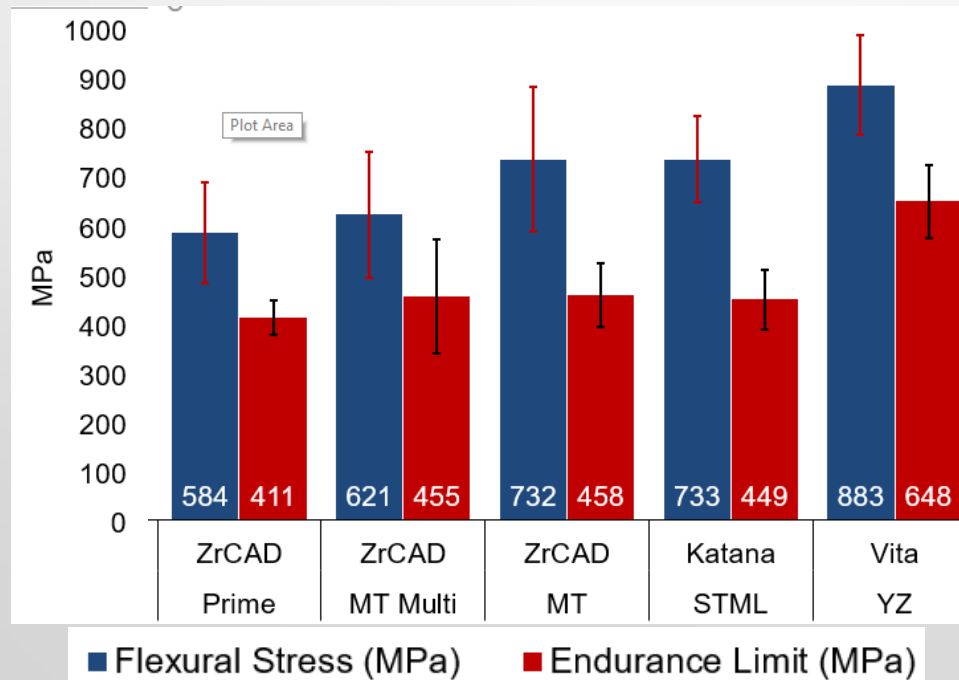

Conclusion

- Oxide-based, zirconia exhibited higher values for all properties evaluated when compared to silica-based ceramics
- Tessera vs IPS e.max CAD
 - Tessera has significantly flexural stress value
 - Tessera is comparable to IPS e.max CAD is fatigue behavior and endurance limits

The background is split diagonally. The top-left portion is dark grey and filled with a dense, overlapping pattern of 3D dollar signs (\$). The bottom-right portion is a light grey gradient. A series of parallel diagonal lines, consisting of a blue line, a dark grey line, and a thin white line, separate the two background sections.

Questions?